

-200 to 670°C

- Useable Range -200°C to 670 °C
- 25.5 Ohm SPRT
- Outstanding Performance

The 670SQ is our latest thermometer, to be specifically designed to give optimum performance up to the aluminium point. Its construction permits the four internal platinum lead wires to expand and contract in the same manner as those of silver-point thermometers. The all-quartz construction of the support members gives the most stable performance with minimal drift, and a unique platinum radiation shield prevents heat radiating up the inside of the sheath.

The temperature range and design of this new unit means that we can now offer 25.5 Ohm (Ro) and 100 Ohm (Ro). The construction, including the coiled sensing element, heat-shunt baffles and light scattering barriers, creates a thermometer of unsurpassed stability.

Because the 670SQ goes beyond the temperature range of oxide growth to the level at which the oxide dissociates, the 670 is filled with a unique argon/oxygen mixture. A 2.5 metre length of low thermal EMF, high temperature, screened cable is connected in the handle, via a strain-relieving transition, to the all pure platinum construction of the thermometer.

Gold-plated U-shaped terminals complete the cable construction, and the 670SQ is delivered in an elegant soft lined carry-case of our own design. A 670SQ is supplied only after a stabilising process which is complete when the reproducibility of R_{TPW} is within 0.0005°C after excursions to the extremes of its temperature range. Values of R_{TPW} and W_{ga} are routinely provided with the 670SQ.

The 670SQ can be supplied with R_{TPW} and W_{ga} only or with full UKAS calibration. "With calibration" means that you will get an Internationally accepted Fixed Point calibration. For best accuracy, recommended maximum measuring currents for the 670SQ are 1mA for the 25.5Ω (Ro) and 0.5mA for the 100Ω (Ro).

A comprehensive handbook and tutorial will help you get the very best performance and stability from your 670SQ.

The 670SQ 650mm long is our recommended SPRT offering ultra stability, vibration, shock, immersion and self heating characteristics. From the success of the original Model 670 SPRTs we have introduced new models into the 670 range offering metal sheathed and low temperature models.

Primary SPRT 670 Family



*Super Stable Standard Platinum
Resistance Thermometer*

| | |
|--------------------|---|
| Models | 670SQ, 670SH, 670SL |
| Measuring Range | -200°C to 670°C |
| Nominal Resistance | 25.5Ω Ro or 100Ω Ro |
| Resistance Ratio | $W_{ga} > 1.11807$ as required by ITS-90 |
| Sensitivity | 0.1Ω / °C (25.5Ω) 0.4Ω / °C (100Ω) |
| Long term drift | from 0.001°C / year depending on use |

How to Order

Model 670SQ, 670SH, 670SL / 25.5 or
Model 670SQ, 670SH / 100

State "with UKAS Calibration" or "without UKAS Calibration".

Refer to introduction for Calibration Uncertainties.


Whilst metal sheathed thermometers appear more robust than the quartz glass models it should be noted that ALL SPRTs are fragile devices and must be handled with care. Quartz glass thermometers have the advantage that the internal components are visible and can be inspected and continue to be our recommended models.

The low temperature models have excellent immersion characteristics and a significant cost saving when compared to the higher temperature models.

Having selected the finest cells and apparatus select Isotech SPRTs, with over thirty years experience of manufacturing platinum resistance thermometers it is not surprising we have developed the finest standard thermometers.

Our preferred standards in our UKAS lab are the 670SQ for up to 660°C and the 96178/0.25 for use up to the silver point, 962°C.

Benefits to your laboratory, ultra stability, best SPRT vibration and shock resistance, best immersion characteristics, low self heating, longer life and less contamination.

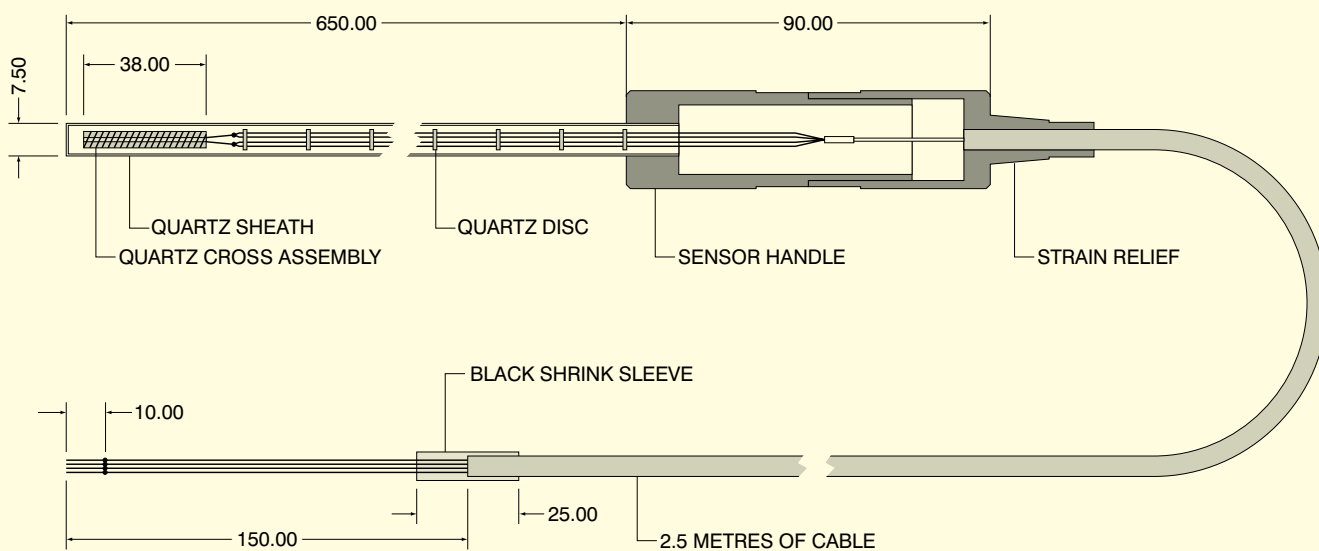

| Model | Range (°C) | Ro (Ohms) | Sheath | Diameter | Length | Sensing Length | Comments |
|-------|-------------|-----------|--------|----------|----------------------|----------------|--|
| 670SQ | -200 to 670 | 25.5 | Quartz | 7.5mm | 650mm or 480mm | 35mm | Recommended for wide range use in the Primary Laboratory |
| 670SQ | -200 to 550 | 100 | Quartz | 7.5mm | 650mm or 480mm | 35mm | 100 Ohm version |
| 670SH | -80 to 670 | 25.5 | Metal | 6mm | 650mm or 480mm | 35mm | Metal sheathed high temperature model |
| 670SH | -80 to 550 | 100 | Metal | 6mm | 650mm or 480mm | 35mm | 100 Ohm version |
| 670SL | -200 to 165 | 25.5 | Metal | 6mm | 480mm | 35mm | Metal sheathed low temperature model |