


Calibration Solutions for

- Temperature Sensors - Infrared Thermometers - Thermocouple Referencing Systems


Temperature Calibration Equipment & Services

Thermometer & Thermocouple Reference

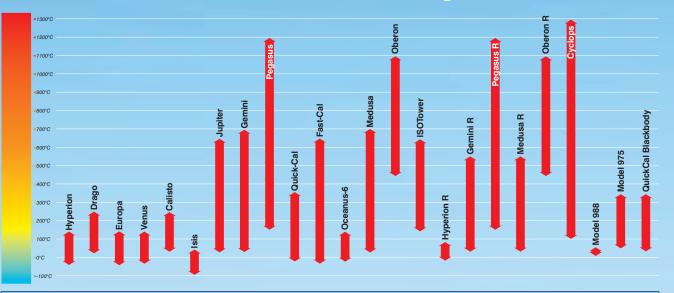

Minimum Calibration Points to Enable Coefficients to be calculated for a Platinum Resistance Thermometer

IPRTs	Required P	oints				Number of Points
Temp > 0°C	0°C and two	positive values	Rt = Ro(1 + A)	t + Bt²)		3
Temp < 0°C to > 0°C	0°C, two pos one negativ	sitive values and e value	Rt=Ro[1 + A	$t + Bt^2 + C(t-100^{\circ}C)$	C) t ³]	4
ITS-90	Required P	oints				
0.01 – 29.76°C	0.01°C	29.76°C				2
0.01 – 156°C	0.01°C	156°C				2
0.01 – 232°C	0.01°C	156°C	232°C			3
0.01 – 419°C	0.01°C	232°C	419°C			3
0.01 – 660°C	0.01°C	232°C	419°C	660°C		4
0.01 – 962°C	0.01°C	232°C	419°C	660°C	962°C	5
-38 – 29.76°C	-38°C	0.01°C	29.76°C			3
-189 – 0.01°C	-189°C	-38°C	0.01°C			3

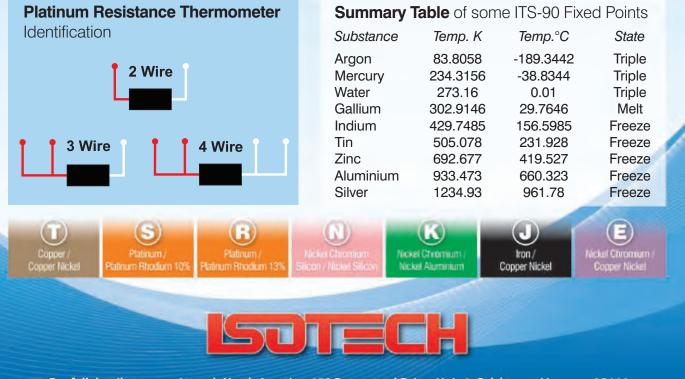

For full details, contact Isotech North America, 158 Brentwood Drive, Unit 4, Colchester, Vermont 05446 Telephone +1 (802) 863.8050 Fax +1 (802) 863.8125 Email sales@isotechna.com Web www.isotechna.com



Contents

Thermometer & Thermocouple ReferenceIns Introduction	
A Guide to Industrial Temperature Calibration	6
Alternative methods of using Portable Calibrators	7 - 9
Introducing the 4000 Series. Hyperion / Drago Europa / Venus / Calisto Jupiter Gemini Pegasus	. 16 - 19 . 20 - 23 . 24 - 27 . 28 - 31
Fast Calibrators Quick-Cal Fast-Cal	. 38 - 39
Cryogenic Dry Block Isis	. 42 - 43
ITS-90 Fixed Point Cells ISOTower ITS-90 Slim Cells Oceanus-6 510 Medusa & 511 Medusa 3 Oberon	46 - 48 49 - 50 51 - 53 54 - 56
Precision Thermometers TTI-10 millik TTI-22 Terminal Adaptor Minature Fixed Resistor - Model 836 Selector Switch - 8 Way Surface Temperature Measurement - Model 944 . Small Hot Plate	. 60 - 61 . 62 - 65 . 66 - 67 . 68 - 69 70 . 71 . 72 - 73 74

	Platinum Resistance Thermometers		
	Thermocouples	79 -	81
Ter	nperature Calibration Software		82
	Cal NotePad		
	I-Cal Easy	. 84 -	85
	S.V.		
Bla	ickbody Sources		
	Blackbody Fixed Point Cells		
	Gas Flow System - Model 984		
	Hyperion R		
	Gemini R		
	Pegasus R		
	Medusa R		
	Oberon R		
	Cyclops		
	Blackbody Source - Model 988		
	Greybody Source - Model 975		
	QuickCal Blackbody		99
Th	ermocouple Referencing Techniques	1	00
	Isotech Custom Engineering		
	Cold Junction Thermocouple Probes		
	TRU - Model 938		
	TRU - Model 937		
	TRUrac - Model 847		
	Hotbox - Model 830	1	06
	Isobox - Model 842		
	Isorac - Model 844		
Th	ermocouple Reference Chart	1	09
Lal	boratory Services 1	10 - 1	111
lso	tech Training Course 1	12 - 1	13
Ind	lustrial Platinum Resistance		
	Thermometer Tables 1	14 - 1	15



The company is always willing to give technical advice and assistance where appropriate.

Equally because of the program of continual development and improvement, we reserve the right to amend or alter characteristics and design without prior notice.

This publication is for information only.



Introduction Industrial Temperature Calibration

...our knowledge, product quality and integrity are the reasons for buying from Isotech

This catalogue describes Isotech products designed for Industrial Temperature Calibration. It has equipment for calibrating industrial sensors of widely varying size, temperature range, accuracies and type. In this book calibration solutions exist for accuracies from 1°C to less than 0.001°C, temperature ranges from below -100°C to 1200°C, from handheld to bench top devices. Products include Dry Block Calibrators, Liquid Baths, High Temperature Furnaces, Precision Temperature Indicators, Temperature Sensors with supporting services, UKAS Calibration and software.

Benefit from Isotech's premier UKAS Calibration Laboratory

Many manufacturers of temperature calibration equipment do not operate accredited laboratories. At lsotech we remain unique in that in addition to manufacturing equipment we also operate a full scale UKAS accredited laboratory calibrating both our own equipment and a full range of temperature equipment for our clients, from industrial thermocouples to primary standards.

Benefit from Experience

Isotech fully understand the calibration needs of our customers. We design and specify products accordingly. Isotech make available full evaluation reports with uncertainty calculations and publish uncertainty graphs. Dry Blocks from other manufacturers are sometimes specified in interesting ways, such as Minimum Operating Temperature -30°C but with an ambient of 5°C, Isotech labs and workshops are kept close to 20°C! We don't specify performances of thermometers only at 0°C we also tell you what it is like at high temperatures. Isotech understand the difference between resolution and accuracy, between accuracy and uncertainty, between "specmanship" and the firm evidence an auditor would expect to see.

As well as solutions for industrial sensors we have a full range of laboratory equipment and provide Primary Standards and Equipment to National Metrology Institutes all around the world

Solutions for Temperature Laboratories

A separate publication "Volume 1: Solutions for Primary & Secondary Laboratories" contains details of our laboratory equipment including Fixed Point Cells, Thermometry Bridges, SPRTs, Comparison Calibration Equipment and Software.

Isotech Innovation

Isotech is renowned for innovation, with products like the ISOCAL-6 - the multi function calibrator providing Dry Block, Liquid Bath, Surface Sensors, Blackbody and Fixed Point calibration in one piece of apparatus.

UKAS Accredited Laboratory (ISO/IEC 17025:2005)

Benefit from Isotech's experience and expertise, our laboratory was first accredited in 1985. The Isotech team of calibration engineers are available for pre and post sales support.

http://www.isotechna.com



The world's leading National Metrology Institutes choose Isotech - shouldn't you?

Industrial Calibrators -100°C to 1200°C

Isotech's range of calibrators for industrial probes includes a range of multi function calibrators that can calibrate all types of thermometer, resistance thermometers, thermocouples, surface sensors and infrared thermometers, the ISOCAL-6 Family.

A range of Dry Block Calibrators offer a variety of volumes from 25 x 115mm in a handheld device to 50 x 300mm in bench top apparatus.

■ ITS-90 Fixed Points -38.8344°C to 961.78°C

A range of Slim Fixed Point cells provide calibration uncertainties to less than 1mK, 0.001°C. These cells are lower cost, more robust and easier to use than the larger cells used by National Metrology Institutes. The affordable prices justify the use of ITS-90 Primary Standards to a wider range of laboratories.

The new ITS-90 Isothermal Towers are integrated devices combining cell and apparatus in an innovative design, the ultimate in simple to use high performance ITS-90 fixed point calibration.

■ Thermometer Read Outs -250°C to 2315°C

For Reference Platinum Resistance Thermometers and Thermocouples there is a range of precision thermometers and thermometry bridges offering performance to 1mk, 0.001°C. There is also a Surface Temperature Thermometer that uses a thermally compensated probe to permit high accuracy surface measurement.

■ Reference Probes -200°C 1300°C

A range of probes covering from -200°C to 1300°C exist for use alongside the calibration products and for general purpose applications.

Software

Isotech have software to totally automate the calibration of temperature sensors, from recording the sensor output to printing certificates and calculating coefficients.

Blackbody Sources

For calibrating infrared thermometers a dedicated series of blackbody sources are available covering up to 1300°C. A range Blackbody Cells, from Gallium 29.7646°C to Silver at 961.78°C allow calibration to ITS-90 Fixed Points.

Thermocouple Referencing

Models include those for use with Laboratory Standard Thermocouples to industrial systems and multi channel reference systems suitable for the largest of installations applications including power stations, boiler rooms, environmental monitoring and research.

For temperature calibration equipment for laboratories refer to our separate publication, "Solutions for Primary & Secondary Laboratories".



ISOTECH

A Guide to Industrial Temperature Calibration: **Traceable Calibration**

For best practice, the thermometer (or thermometers) under test are placed into the calibration volume alongside a calibrated standard. This is so that the test thermometers "can be related to appropriate standards, generally international or national standards, through an unbroken chain of comparisons". This "traceability" meets the requirements of quality systems including that of ISO 9000.


Using the Calibrator itself as the Reference (or standard) raises a number of issues, such as how is temperature difference between the test thermometer and the calibrator display determined - how can this 'uncertainty value' be known?

International Guidelines have been published from EURAMET that give guidance, and requirements, for the calibration of Dry Blocks EURAMET/ cg-13/v.01 (formerly EA10-13) and Isotech Calibrators meet the calibration capacity requirements of EURAMET/cg-13/v.01.

We recommend that a reference probe is used, the same method as used in secondary temperature laboratories. For less demanding calibration, and the quick testing of sensors, the Calibrator can be used without a reference probe, refer to the units Evaluation Report for typical performance data.



Dry Block Calibrator of Poor Design

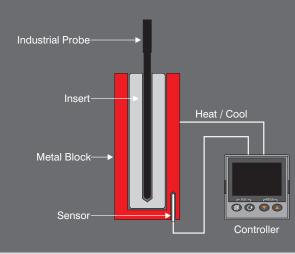

Pre-purchase check list



1 Does the supplier have an accredited laboratory? UKAS accreditation, "the means by which, in the public interest, the integrity and competence of independent evaluators is confirmed and declared". Isotech can issue a UKAS certificate with the performance expressed in the manner that you will need, not to some confusingly expressed specification that is made with no confirmation of integrity and competence.

2 Experience

Does the producer have experience? Do they understand the difference between accuracy and uncertainty? Can they tell you how to calculate the uncertainty of a probe being calibrated in the block? Isotech can.

3 Expandable

Can the Dry Block be used with other sensors? Are there accessories available for future expansion? With Isotech products they are.

4 PC Support

Can it be connected to a computer? Is there software available, can it be automated? Isotech Dry Block Calibrators have a range of software options.

5 Documented

Is the bath fully documented? Can you download a full evaluation report from the Web Site? Does it come with a comprehensive handbook and tutorial? Is training available? Isotech provide all of these free of charge.


6 Practical

Isotech Dry Blocks are practically designed with a strong metal case, and are a compact portable size. If you are going to carry it around don't forget to check the size and weights. It is surprising how large some other blocks are, even though they take the same number of probes. Beware if the specification does not include the weight.

7 Value

Check the prices, all the above come at an amazingly competitive price when you choose **Isotech**.

Dry Block Calibrator Meeting ISO9000 Requirements

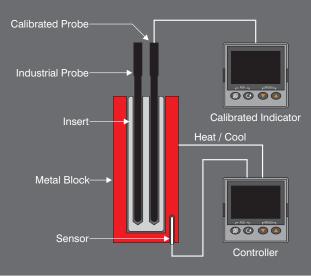

Alternative methods of using **Isotech Portable Calibrators**

A Basic Dry Block Calibrator


The thermometer under test is compared to the dry block controller value. Useful for moderate temperature ranges and quick testing.



A control system with digital display shows the set and measured value from the calibration volume.

An ISO 9000 Calibration System

A thermometer under test is compared to a calibrated standard, for true traceability and clearly meets the requirements of ISO9000

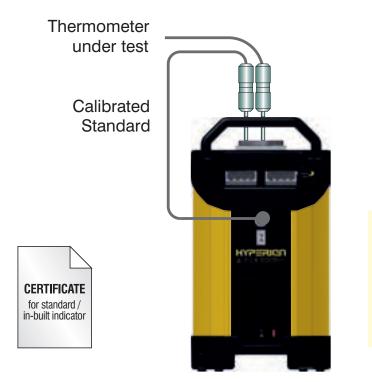


In addition to the temperature control system with its display the Isotech SITE model includes a separate temperature read out for the Standard (Reference) probe.

ISOTECH


Isotech ADVANCED Calibrator

Isotech ADVANCED models have inputs for both the test thermometer and a calibrated standard.

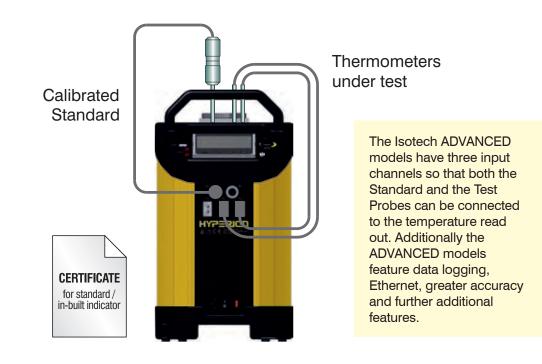

Using an External Indicator


Similar to the previous configuration but an external indicator is used - one instrument can be used with many calibration baths - the bath or baths do not need a calibration certificate, but they need an evaluation report.

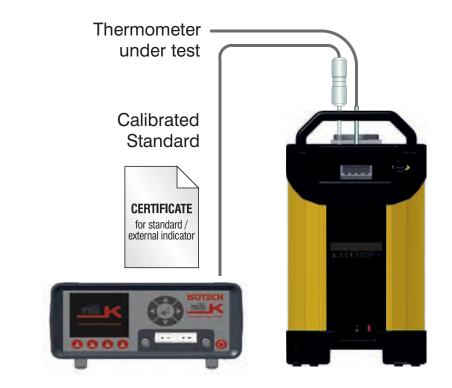

CALIBRATION SOLUTIONS

An example of Multiple Sensor Calibration

You can add a scanner for multi-probe calibration - the system can then be automated.

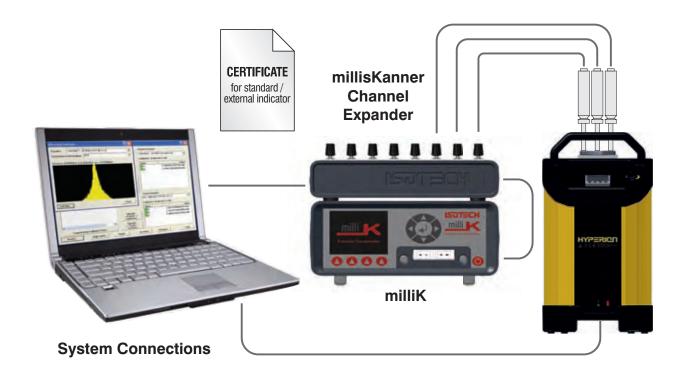

Isotech - award winning flexible solutions for all sensor types from -45°C to 1200°C By adding different accessories the 4000 Series can be used in up to six different modes - Dry Blocks, Liquid Baths, Ice Bath, Surface Sensor Calibrator, IR Thermometer Calibrator and even with ITS-90 Fixed Point

	Ø 65mm	Ø 65mm	Ø 35mm	@ 35mm	@ 35mm	Ø 65mm	@ 35mm	Ø 33.5mm
			lsocal-6	I	I		Dry Block Calibrators	
	HYPERION	DRAGO	EUROPA	VENUS	CALISTO	GEMINI	JUPITER	PEGASUS
Specifications								
Metal Block Bath	\checkmark	~	\checkmark	\checkmark	\checkmark	~	\checkmark	~
Stirred Liquid Bath	\checkmark	~	\checkmark	\checkmark	\checkmark			
Stirred Ice/Water Bath	\checkmark		\checkmark	\checkmark				
Blackbody Source	\checkmark	~	×	~	~	~	~	\checkmark
Surface Sensor	\checkmark	~	\checkmark	\checkmark	~		~	
ITS-90 Fixed Point	\checkmark	~	\checkmark	\checkmark	\checkmark			
Temperature Range (°C)		_						
1200°								150°C → 1200°C
800°						35°C → 700°C	35℃ → 660°C	
400°	-25°C → 140°C	30°C → 250°C	-45°C → 140°C	-35°C → 140°C	<u>30°C → 250°C</u>			
-100°	-20 0 -7 140°6							

SOTE

Isotech 4000 Series Available in three different versions

Isotech calibrators are available in three different versions; they all feature the same rugged casing with award winning calibration volumes. The range of features varies with the ADVANCED offering extra performance in terms of resolution and temperature stability. All models have common block sizes, interchangeable options and benefits providing calibration volumes with superior temperature uniformity and capacity for industrial probes.

Choose to suit your application:

SOTECH

ODEL

E

BASIC



Features

- Calibrate the entire loop by using a heat source rather than an electrical simulator, a test instrument and sensor can be calibrated as a system
- Simple To Use and Outstanding Value
- Rugged Case
- Calibrate all Sensor Types
- Ramp to Set Temperature
- Supply Voltage Power Correction and Digital Filtering
- PC Interface and Software

Features

As Basic but adds independent temperature indicator for reference probe

- Accepts Platinum Resistance Thermometers, Process Inputs Including 4 - 20mA, Thermocouples including Types K,N,R,S,L,PL2,T,J and E.
- Thermostat Testing, Stand Alone or with PC
- Configurable Units °C, °F or K

ISOTECH

Best performance

- Benefit from advanced temperature controller that provides best performance
- Resolution to 0.001°C Superior Stability
- Control parameters automatically optimised with temperature

View easily in all conditions

- Large Bright Colour Display
- Similar to Smartphone
- Crystal clear display with full colour graphics



Thermostat testing

Test Two Thermostats Simultaneously

View from anywhere

21st Century Connectivity with Ethernet and inbuilt webserver



Connect to the network and view the calibrator from anywhere, on your Notepad, Tablet or Smartphone

Save time and money

AUTOMATIC TEMPERATURE CYCLING

Program the ADVANCED with the temperature points you need, store commonly used ranges. The calibrator can then automatically follow these points and log your data

Safely store and secure all the data you need

- Massive Internal Memory can safely store all your data
- Store a lifetime of data on a USB Drive
- Choose Open CSV Files or Tamper Proof Data with Secure File Format for Data Security

Supports five languages

- English, French, German, Italian and Spanish
- Simple to use with clear user interface

Connect more probes

UP TO 3 INPUT CHANNELS

- Inbuilt reference system with two universal inputs for PRT, Thermocouple or Process Input and a third thermocouple input
- Input types: Process Inputs including 4-20mA, PRT, and Thermocouple Types B, C, D, E, J, K, L, N, R, S, T, U

Offset elimination

- Connect a reference probe to one of the external probe inputs and the calibrator can trim the block temperature to remove offsets
- Block adjusted to reference probe value
- Remove offsets
- Use in combination with automatic cycling

The software you need

PC tools for reviewing data with support for secure file format, managing temperature programs, easy configuration and data logging.

Introduction: 4000 Series

The 4000 series includes the Isocal-6 multifunction calibration system

The award winning ISOCAL-6 consists of a range of temperature calibrators designed to calibrate all temperature sensors. As a multi function temperature calibrator it can be used as a Dry Block with accessories added to allow use as a Stirred Liquid Bath, a Blackbody Source, a Surface Temperature Calibrator, an ITS-90 Fixed Point System and for low temperatures a Stirred Ice Bath.

ISOTECH

The ISOCAL-6 is a complete temperature calibration laboratory in a simple easy to use package.

A system designed to expand with you, to fulfill all future calibration needs. Giving the flexibility to add accessories when needed and meet current budgetary demands.

Metal Block Bath

A Metal Block Bath, (Dry Block Calibrator) provides fast and clean calibration of thermocouples, PRTs and other industrial temperature sensors. Isotech blocks use a combination of multi zone and advanced materials technology to ensure constant temperature zones to enable high accuracy calibration. Interchangeable 35mm diameter blocks allow several sensors to be calibrated simultaneously with fast heat up and cool down. For larger probes blocks are available up to 65mm diameter and with immersion depth of up to 300mm. An unmatched combination of leading performance and calibration capacity.

Stirred Liquid Bath

Remove the metal block and the Isocal-6 can be converted to a stirred liquid bath. Liquid bath operation allows angled or awkward shaped probes to be calibrated. Accuracies are improved over Dry Blocks alone and with a suitable reference probe performance of 0.005°C is achievable. In stirred liquid bath mode a reference probe should be used.

Stirred Ice / Water Bath

The ISOCAL-6 models that operate below 0°C can be used to provide a 0°C stirred ice / water bath. This provides a simple low cost way of checking that standards have not drifted in between calibrations.



Blackbody Source

Adding the blackbody target allows the testing of infrared thermometers. Low cost non-contact IR thermometers are increasingly being used in industry and the ISOCAL-6 is ideal to test and check these devices. The IR thermometer is focused on the target and compared to a reference probe in the block pocket.



Surface Sensor Calibrator

With the Surface Sensor Kit the test sensor is compared to a platinum resistance thermometer located just below the surface of the block. Again save the cost of buying additional equipment by adding accessories as required to expand the ISOCAL-6 for new calibration applications.

ITS-90 Fixed Point Apparatus

For the best possible performance with uncertainties to 0.0005°C (0.5mK) add an ITS-90 Fixed Point Cell. The most popular is the B8 Water Triple Point Cell, it is surprisingly affordable and simple to use - the triple point can be both created and maintained in the apparatus without the need for any other equipment or supplies.

4000 Series

ISOTECH

Calibration Options

Isotech's UKAS accredited calibration laboratory was established in 1980 and has grown to be a full scale laboratory providing calibration to the smallest of uncertainties. Isotech certificates are formally recognised in over 70 countries and carry the ILAC-MRA logo.

We can provide different certificates, reporting the temperature only (UKAS-TEMP), additionally we can include the inbuilt temperature indicator channels calibrating by electrical simulation, (UKAS-SYST) or we can just calibrate the input channels (UKAS-SIM).

CRUTHCATE OF CLARRATION Region Nondexect transcioner charge cargo and charge carg	
CLEANING CRUTCHING OR CLEANING CRUTCHING	
IDENTIFICATION Ga 454	
The second second second second and second s	-

The following Calibration Option	Code	
BASIC, SITE and ADVANCED	5 point calibration for block temperature; includes reference probe values when ordered with Site or ADVANCED	UKAS-TEMP
ADVANCED	5 point calibration for block temperature and reference probe (when ordered) and electrical simulation of indicator	UKAS-SYST
ADVANCED	Calibration of input channels, electrical simulation only	UKAS-SIM

Included Software

Model	Software	Connection
BASIC and SITE	Cal Notepad and I-Cal Easy LOG: Log and Monitor Calibrator	Serial
ADVANCED	I-Cal Easy LOG: Log and Monitor Calibrator	Ethernet
ADVANCED	Review Lite: Manage Log Data, supports secure data format	Ethernet
ADVANCED	Set Point Program Editor: Manage ADVANCED set point programs	Ethernet
ADVANCED	Isotech ADVANCED Block Configuration Utility: Fast configuration from a PC	Ethernet

Download: http://www.isotech.co.uk/downloads

Isotech Innovation

We have been providing calibration solutions for more than 30 years, from Primary Standards for National Metrology Institutes through to handheld calibrators for service engineers.

Our portable calibrators have set a number of world firsts, the first Dry Blocks to feature an independent inbuilt temperate indicator and the first Dry Blocks to reach -100°C. The award winning ISOCAL-6 was the first calibrator to offer both Dry Block and Liquid Bath use. It remains the only device to offer six modes of calibration.

Isotech Dry Block Range - expanded and upgraded

- New Rugged Lightweight Outer Cases
- Custom moulded components brings lighter and stronger Portable Calibrators
- New ADVANCED model more features than ever before

The ADVANCED Models include automatic temperature cycling which saves time and money with the calibrator automatically logging data over a series of calibration points

Temperature Cycling and **Logging** Features of the **ADVANCED** models

Automatic Temperature Cycling

- Save Time
- Save Money
- Store up to 100 Programs

Calibration temperatures and dwell times can be combined into 'programs' for the ADVANCED calibrator to follow, up to 100 programs can be stored, with each program having up to 25 segments.

Programs can be edited from the front panel or on a PC with Set Point Program Editor.

As the program runs all measurements are recorded to the large internal memory.

 123
 - Act
 C
 D
 F
 H
 H
 S
 M
 N
 P
 R

 I Oueredi Love
 3000 Hight
 800 Undar %
 500 Hight
 800 Hight

Automatic temperature cycling saves time and money with the calibrator automatically logging data over a series of calibration points

Data Logging

- Large Internal Memory to Store Your Data
- Creates CSV Files For Spreadsheets
- Creates Secure Tamper Proof Data to aid compliance to standards


The ADVANCED models log data to internal memory, the data is logged in a CSV format that can readily be opened and edited.

Data is also recorded in a secure tamper proof format as a binary 'UUH File'. These files can be opened with Review Lite software – download from our website. If the records

are found to have been tampered with then Review Lite will reject the record as corrupt.

Data files can be shared automatically across the network via FTP, transferred manually via USB or 'pulled' into Review's database.

For users working in regulated industries these features, when used in conjunction with company procedures and other equipment or process, aids compliance to standards such as FDA 21 CFR 11 and AMS2750E.

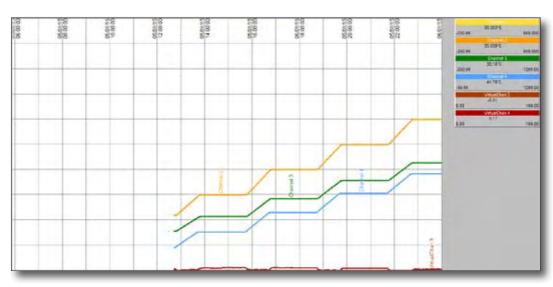





ISOCAL - 6 Range Hyperion • Drago

- Multi Function: Six Modes including Dry Block and Liquid Bath
- 65mm Volume: Ideal Liquid Bath
- Calibrate entire measurement loop using a heat source rather than an electrical simulator, a test instrument and sensor can be calibrated as a system

The Hyperion and Drago have large calibration volumes, 65mm x 160mm deep, which makes them ideal to use as portable liquid baths. Stirred liquid baths are suitable for temperature sensors of all types, sizes and shapes. Liquid Baths can provide smaller calibration uncertainties than dry blocks and when used with suitable reference thermometers, accuracies of up to 0.005°C can be achieved.

These models are part of the award winning Isocal-6 family and with a reference probe can be used with different accessories for Dry Block, Infrared, Surface Calibration and even with ITS-90 Fixed Point Cells for uncertainties to 0.001°C. In Dry Block Mode, the large 65mm diameter block allows for the calibration of either larger probes or for calibrating many sensors simultaneously.

As a Liquid bath the sensors can be placed directly into the stirred liquid thus avoiding the need for specially drilled blocks. If the liquid is directly in the block then the controller only model, or Basic (B) model, can be selected. However, instead of putting liquids directly in the block liquid containers can be used to facilitate rapid change of fluids. For greater accuracy, or when using a liquid container, Dry Block Insert, Blackbody Target or the Surface Sensor Kit a separate reference thermometer should be used to compensate for the varying offset between the controller and the accessory temperature.

An ideal arrangement would be to include the handheld Isotech TTI-10 or the bench top Isotech millik Precision Thermometer and Model 935-14-16 Semi Standard Platinum Resistance Thermometer.



Alternatively the SITE or ADVANCED model can be selected; the SITE includes a temperature indicator for a

selected; the SITE includes a temperature indicator for a reference probe. The ADVANCED also includes inputs for test thermometers, automatic temperature cycling, logging and additional sophisticated features.

All models include I-Cal Easy LOG software and the ADVANCED models additionally include software to manage logged data and configure the unit, see page 14 for more details.









PARAMETER	Model			
	Hyperion 4936	Drago 4934		
Temperature Range	-25°C to 140°C1	30°C to 250°C ²		
ADVANCED Range				
Stability: Dry Block / Liquid Bath	±0.005°C	±0.005°C		
Display Resolution	0.001°C over whole range	0.001°C over whole range		
Accuracy: RTD Input Channels	±0.05°C ±0	0.005% RDG		
Accuracy: Thermocouple Input Channels	E,J,K,N: ±0.2°C @ 660°C R: ±0.6°C S	S: $\pm 0.7^{\circ}$ C @ 660°C T $\pm 0.2^{\circ}$ C @ 150°C		
CJC Accuracy	±0.3	35°C		
BASIC / SITE Range				
Stability	±0.03°C	±0.03°C		
Display Resolution	0.01°C from -19.99 to 99.99°C the	en 0.1C: 0.01°C Over PC Interface		
COMMON Specifications				
Stability	Blackbody ±0.3°C Surface Sense	or ±0.5°C ITS-90 Cells ±0.0005°C		
Display Accuracy ³	0.15°C	0.15°C		
Uniformity - Radial, Liquid Bath Mode	<0.009°C	<0.007°C		
Uniformity - Axial, Liquid Bath Mode (40mm)	<0.011°C	<0.013°C		
Uniformity - Radial, Dry Block Mode (Between Wells)	<0.008°C	<0.008°C		
Uniformity - Axial, Dry Block Mode (40mm)	<0.040°C	<0.040°C		
Heating Time	-20°C to 140°C: 40 Mins	30°C to 250°C: 40 Mins		
Cooling Time	140°C to 20°C: 90 Mins 20°C to -25°C: 80 Mins	250°C to 30°C: 90 Mins		
Insert Size	65 x 160mm			
Insert Types	Standard 8 x 8mm + 2 x 4.5mm, Undrilled or Custom Drilled			
Power	115 or 230Vac 50/60Hz 200 Watts	115 or 230Vac 50/60Hz 1000 Watts		
Dimensions	384H (including handle) x 212W x 312D mm			
Weight	12kg	8kg		

(1) In ambient of 20°C: Minimum Temperature is 45°C Below Ambient, Absolute Minimum -35°C
 (2) In ambient of 20°C
 (3) Dry Block Mode only: Comparing 4.5mm Well to Display Value.

	ADVANCED	SITE	BASIC
		Sur?	
Digital Display of Set and Nominal Block Temperature	Yes	Yes	Yes
PC Interface	Ethernet + USB Host	Serial	Serial
Test Thermostats	Yes - Two Inputs	Yes - Single Input	No
Independent Temperature Indicator for Reference Probe	Yes	Yes	No
Additional Inputs for Units Under Test	Up to 3: Two universal inputs for PRT, Thermocouple or Process inputs and a further Thermocouple input	No	No
Automatic Temperature Cycling	Yes	No	No
Data Logging	Yes - Export to USB	No	No
Offset Elimination	Yes - block can follow reference input	No	No
Choose English, French, Italian or Spanish Language	Yes - on full colour display	No	No
In Built Web Server	Yes	No	No
Tamper Proof Data	Yes - Suitable for life science, automotive and aerospace applications	No	No








Metal Block Bath

Dry Block Calibrator provides fast and clean calibration of thermocouples, PRTs and other industrial sensors. Isotect blocks use a combination of multi zone and advanced materials technology to ensure constant temperature zones for high accuracy calibration.



Stirred Liquid Bath

Remove the metal block to convert to a stirred liquid bath. Liquid bath operation allows angled or awkward shaped probes to be calibrated. Accuracies are improved over Dry Blocks alone and with a suitable reference probe performance of 0.005°C is achievable.

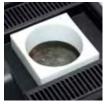

Stirred Ice / Water Bath

The ISOCAL-6 models that operate below 0°C can be used to provide a 0°C stirred ice / water bath. This provides a simple low cost way of checking that standards have not drifted in between calibrations.



Blackbody Source

Adding the blackbody target allows the testing of infrared thermometers. Low cost non-contact IR thermometers are increasingly being used in industry and the ISOCAL-6 is ideal to test and check these devices. The IR thermometer is focused on the target and compared to a reference probe in the block pocket.



Surface Sensor Calibrator

With the Surface Sensor Kit the test sensor is compared to a platinum resistance thermometer located just below the surface of the block. Again save the cost of buying additional equipment by adding accessories as required to expand the ISOCAL-6 for new calibration applications.



ITS-90 Fixed Point Apparatus

For the best possible performance with uncertainties to 0.0005°C (0.5mK) add an ITS-90 Fixed Point Cell. The most popular is the B8 Water Triple Point Cell, it is surprisingly affordable and simple to use - the triple point can be both created and maintained in the apparatus without the need for any other equipment or supplies.



World's First Multi-Functional Baths SIX FUNCTIONS

Hyperion / Drago Accessories



Dry Block Mode with Inserts

936-06-01a Standard Insert is: 8 x 8mm + 2 x 4.5mm all 157mm Deep. All Inserts have a 4mm tapped hole to suit supplied extractor tool.

Alternative Inserts 936-06-01b Blank Insert 936-06-01c Special Insert.



Stirred Liquid Mode with Liquid

Container Kit 936-06-02 Allows liquid bath use, includes container, magnetic stirrer, probe guide and sealing cap.

Stirred Ice Bath Mode with Liquid Container Kit

Uses same liquid kit to provide 0°C reference as a stirred ice bath.

Thermometer Support Kit 936-06-08 Supports up to eight thermometers into liquid. Suits probes 5mm - 8mm in

liqu dia 936 580 915

ımeter.			
6-06-07	C10 Oil	-35°C – 140°C	1L
0-06-09	C20 Oil	20°C – 200°C	1L
5/09	VH Oil	150°C – 250°C	1L



Infrared Calibration Mode with Blackbody Target 936-06-03

Use optional Probe 936-14-61DB.



Surface Sensor Calibration with Surface Sensor Kit 936-06-04 Includes an Insert and an angled platinum resistance thermometer.



ITS-90 Fixed Point Cells B8 Water Triple Point Cell (Hyperion) 17401 Slim Gallium Slim Cell 936-06-09 Cell Holder Assembly

Calibration



Includes three point traceable calibration certificate for block temperature **UKAS Calibration** UKAS Calibration available to order, legally traceable in more than 70 countries.



Standard Probe 935-14-61/DB Platinum Resistance Thermometer. 4mm diameter.




Current Loop Interface 935-06-161 24VDC Power Supply and Terminal Box. Powers 4-20mA Current Transmitters with 4mm terminal posts for easy connection.




Carrying Case 931-22-112 Sturdy case with room for accessories. Features wheels and pull out handle.



Isocal-6 Performance and Use

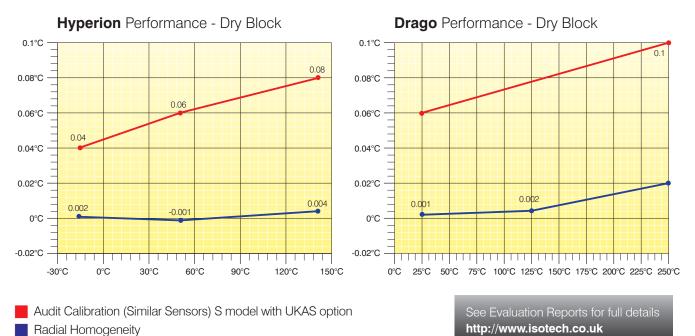


Alternative Methods of Calibrating with an Isocal-6

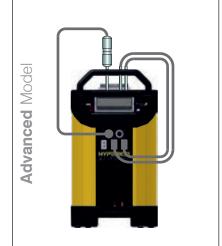


ADVANCED Model

- Digital Display of Set and Nominal Block Temperature
- Inbuilt three channel indicator for reference probe and units under test
- Advanced features including automatic Temperature Cycling and Logging
- Best Practice calibration with established traceability and uncertainty

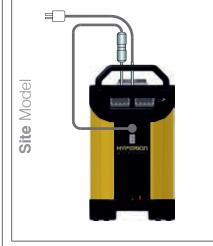

SITE Model

- Digital Display of Set and Nominal Block Temperature
- Inbuilt single channel indicator for reference probe
- Best Practice calibration with established traceability and uncertainty

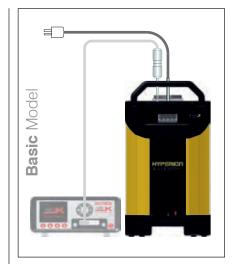

BASIC Model

- For Quick and Easy Testing
- Digital Display of Set and Nominal Block Temperature
- Use with a separate external indicator to compensate for gradients and loading



UKAS Calibration available for these systems - International Traceability - Best Practice See page 14



45 to 250°C



- Multi Function: Six Modes including Dry Block and Liquid Bath
- Fast Response 35mm x 160mm Calibration Volume
- Calibrate Whole Measurement Loop

These models will calibrate temperature probes from -45°C to 250°C with unrivalled flexibility. As a traditional Dry Block, several thermometers can be quickly calibrated.

Accessories are available to convert to a stirred liquid bath, for surface sensor calibration, to calibrate infrared thermometers and even to use as an ITS-90 Fixed Point System with calibration uncertainties as small as 0.0005°C. With excellent stability and distributed heating cooling zones for good uniformity these calibrators offer proven thermal performance.

These award winning calibrators are easy to use and are available in three versions – the Basic, the Site and the ADVANCED. The Basic has a digital display of set and nominal temperature, the Site additionally includes an inbuilt independent temperature indicator for a reference probe. The ADVANCED controller has inputs for reference and test thermometers with a further range of sophisticated features including automatic temperature cycling, secure data logging and full colour high resolution display.

All models include I-Cal Easy LOG software and the ADVANCED models additionally include software to manage logged data and configure the unit, see page 14 for more details.

These models meet the calibration capacity requirements of EURAMET/cg-13/v.01, "EA Guidelines on the Calibration of Temperature Block Calibrators, formerly EA10/13. accuracy and best performance.



http://www.isotechna.com









Parameter	Model				
	Europa 4520	Venus 4951	Calisto 4953		
Temperature Range	-45°C to 140°C (1)	-35°C to 140°C (2)	30°C to 250°C (3)		
ADVANCED Range					
Stability: Dry Block / Liquid Bath	±0.01°C	±0.01°C	±0.02°C		
Display Resolution	0.001°C over whole range	0.001°C over whole range	0.01°C over whole range		
Accuracy: RTD Input Channel		±0.05°C ±0.005% RDG			
Accuracy: Thermocouple Input Channel	E,J,K,N: ±0.2°C @ 660	0°C R: ±0.6°C S: ±0.7°C @ 660°	C T ±0.2°C @ 150°C		
CJC Accuracy		±0.35°C			
BASIC/SITE Range					
Stability	±0.03°C	±0.03°C	±0.03°C		
Display Resolution	0.01°C from -19	0.99 to 99.99C then 0.1C: 0.01C Ov	ver PC Interface		
COMMON Specification					
Stability	Blackbody ±0.3	°C Surface Sensor ±0.5°C ITS-90	Cells ±0.0005°C		
Display Accuracy ⁴	0.15°C	0.15°C	0.25°C		
Uniformity - Between Wells Dry Block Mode (Radial)	<0.008°C	<0.008°C	<0.02°C at 250°C		
Uniformity - Lower 40mm (Axial) Dry Block Mode	<0.040°C	<0.040°C	<0.25°C		
Uniformity - Radial Liquid Bath Mode	<0.02°C	<0.02°C	<0.011°C at 250°C		
Uniformity - Lower 40mm (Axial) As Liquid Bath	<0.026°C	<0.026°C	<0.02°C at 250°C		
Heating Time	-30°C to 140°C: 15 Mins	-30°C to 140°C: 15 Mins	25°C to 250°C: 15 Mins		
Cooling Time	140°C to 0°C: 15 Mins	140°C to 0°C: 15 Mins	250°C to 30°C: 25 Mins		
Calibration Volume		35 x 160mm			
Standard Insert	6 pockets, 2 x 4.5mm, 2	x 6.4mm, 1 x 8.0mm, 1 x 9.5mm dia	ameter, all 157mm deep		
Insert Types	Choice of Three - See Accessories				
	Ether	rnet - supporting software and USB	Host		
CJC Accuracy:		0.35°C			
Dimensions	384	l (including handle) x 212W x 312D) mm		
Power	300 Watts	150 Watts	300 Watts		
Voltage		115Vac or 230 Vac 50/60Hz			
Weight	14kg	10.2kg	8kg		

(1) In ambient of 20°C: Minimum Temperature is 65°C Below Ambient, Absolute Minimum -55°C (2) In ambient of 20°C: Minimum Temperature is 55°C Below Ambient, Absolute Minimum -45°C

(3) In ambient of 20°C
 (4) Dry Block Mode only: Comparing 4.5mm Well to Display Value.

	ADVANCED	SITE	BASIC
		SUN.	
Digital Display of Set and Nominal Block Temperature	Yes	Yes	Yes
PC Interface	Ethernet + USB Host	Serial	Serial
Test Thermostats	Yes - Two Inputs	Yes - Single Input	No
Independent Temperature Indicator for Reference Probe	Yes	Yes	No
Additional Inputs for Units Under Test	Up to 3: Two universal inputs for PRT, Thermocouple or Process inputs and a further Thermocouple input	No	No
Automatic Temperature Cycling	Yes	No	No
Data Logging	Yes - Export to USB	No	No
Offset Elimination	Yes - block can follow reference input	No	No
Choose English, French, Italian or Spanish Language	Yes - on full colour display	No	No
In Built Web Server	Yes	No	No
Tamper Proof Data	Yes - Suitable for life science, automotive and aerospace applications	No	No








Metal Block Bath

Dry Block Calibrator provides fast and clean calibration of thermocouples, PRTs and other industrial sensors. Isotech blocks use a combination of multi zone and advanced materials technology to ensure constant temperature zones for high accuracy calibration.



Stirred Liquid Bath

Remove the metal block to convert to a stirred liquid bath. Liquid bath operation allows angled or awkward shaped probes to be calibrated. Accuracies are improved over Dry Blocks alone and with a suitable reference probe performance of 0.005°C is achievable.

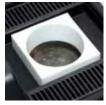


Stirred Ice / Water Bath

The ISOCAL-6 models that operate below 0°C can be used to provide a 0°C stirred ice / water bath. This provides a simple low cost way of checking that standards have not drifted in between calibrations.



Blackbody Source

Adding the blackbody target allows the testing of infrared thermometers. Low cost non-contact IR thermometers are increasingly being used in industry and the ISOCAL-6 is ideal to test and check these devices. The IR thermometer is focused on the target and compared to a reference probe in the block pocket.

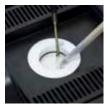

Surface Sensor Calibrator

With the Surface Sensor Kit the test sensor is compared to a platinum resistance thermometer located just below the surface of the block. Again save the cost of buying additional equipment by adding accessories as required to expand the ISOCAL-6 for new calibration applications.



ITS-90 Fixed Point Apparatus

For the best possible performance with uncertainties to 0.0005°C (0.5mK) add an ITS-90 Fixed Point Cell. The most popular is the B8 Water Triple Point Cell, it is surprisingly affordable and simple to use - the triple point can be both created and maintained in the apparatus without the need for any other equipment or supplies.



World's First Multi-Functional Baths SIX FUNCTIONS

Europa Venus Calisto Accessories





Dry Block Mode with Inserts

951-02-15 An Insert is included: (2 x 4.5mm, 2 x 6.4mm, 1 x 8mm & 1 x 9.5mm) x 157mm Deep. All Inserts have a 4mm tapped hole to suit supplied extractor tool.

Alternative Inserts

951-06-07 Alternative Insert type B 13mm, 10mm, 8mm, 5mm and 3.5mm dia. holes, all 157mm deep

951-06-08 Alternative Insert type C 8mm, 6 x 6.5mm dia. holes, all 157mm deep 951-02-15a Blank Insert without pockets for local machining. Includes M4 tapped hole

for supplied extractor tool. **951-02-15c** Custom insert. Isotech can provide custom drilled pockets, minimum of 3mm separation between holes.

Contact Isotech with your requirements



Stirred Liquid Mode with Liquid Container Kit 951-06-01

Allows liquid bath use, includes container, magnetic stirrer, probe guide and sealing cap. Stirred Ice Bath Mode with Liquid Container Kit

Uses same liquid kit to provide 0°C reference as a stirred ice bath (Not Calisto)



Thermometer Support Kit 921-02-06 Allows three thermometers to be suppended in the bath, including liquid in

glass types. **520-05-01** C10 Oil -35°C – 140°C 0.1L **951-06-06** C20 Oil 20°C – 200°C 0.1L **953-04-01** VH Oil 150°C – 250°C 0.1L



Infrared Calibration Mode with Blackbody Target 951-06-04 Use optional Probe 934-14-82/DB placed in the auxiliary block pocket for use as a reference.



Surface Sensor Calibration with Surface Sensor Kit 951-06-02 Includes angled platinum resistance thermometer.



ITS-90 Fixed Point Cells 17724M Slim Mercury Cell (Europa Only). B8 Water Triple Point Cell (Venus and Europa) 17401M Slim Gallium Cell (Europa, Venus and Calisto)



Standard Probe 935-14-82/DB Platinum Resistance Thermometer. Probe diameter 4mm, recommended pocket size 4.5mm. Angled head feature avoids sensors in block.




Current Loop Interface 935-06-161 24VDC Power Supply and Terminal Box. Powers 4-20mA Current Transmitters with 4mm terminal posts for easy connection.




Carrying Case 931-22-111 Sturdy case with room for accessories. Features wheels and pull out handle.



Isocal-6 Performance and Use

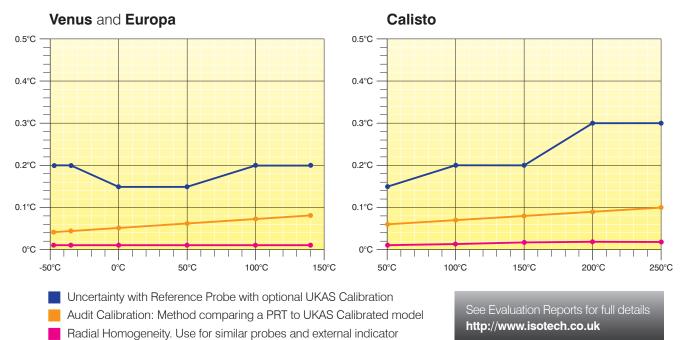


Alternative Methods of Calibrating with an Isocal-6

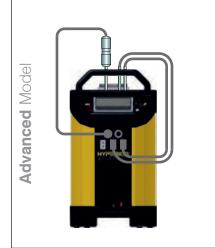


ADVANCED Model

- Digital Display of Set and Nominal Block Temperature
- Inbuilt three channel indicator for reference probe and units under test
- Advanced features including automatic Temperature Cycling and Logging
- Best Practice calibration with established traceability and uncertainty

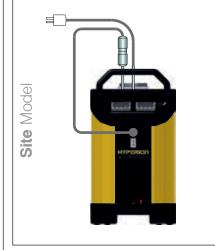

SITE Model

- Digital Display of Set and Nominal Block Temperature
- Inbuilt single channel indicator for reference probe
- Best Practice calibration with established traceability and uncertainty

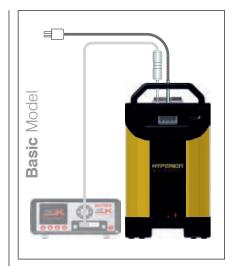


BASIC Model

- For Quick and Easy Testing
- Digital Display of Set and Nominal Block Temperature
- Use with a separate external indicator to compensate for gradients and loading



UKAS Calibration available for these systems - International Traceability - Best Practice See page 14



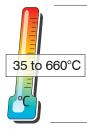

Dry Block Calibrator Jupiter

- Wide Operating Range to 660°C
- Fast Response
- Calibrate Whole Measurement Loop

The Jupiter Dry Block range offers industry-leading performance in an easy to use portable package - ideal for the calibration of thermocouples and platinum resistance thermometers. It has been designed for fast heating and cooling for convenient field use. For flexibility surface sensor and infrared thermometer accessories can be added.

The standard insert can hold up to six thermometers. For larger blocks see the Gemini range.

These award winning calibrators are easy to use and are available in three versions – the Basic, the Site and the ADVANCED. The Basic has a digital display of set and nominal temperature, the Site additionally includes an inbuilt independent temperature indicator for a reference probe. The ADVANCED controller has inputs for reference and test thermometers with a further range of sophisticated features including automatic temperature cycling, secure data logging and full colour high resolution display.

Isotech is a world leader in temperature calibration, providing many nations with their Primary Standards and operates a full scale UKAS accredited calibration laboratory. We can offer a range of calibration options to meet your requirements.

Benefit from our experience and understanding in calibration at all levels, our evaluation reports, our tutorials and uncertainty calculations.

These models meet the calibration capacity requirements of EURAMET/cg-13/v.01, "EA Guidelines on the Calibration of Temperature Block Calibrators, formerly EA10/13.



http://www.isotechna.com

ADVANCED models additionally include software to manage logged data and configure the unit, see page 14 for more details.







All models include I-Cal Easy LOG software and the



Parameter	Model		
	Jupiter 4852		
Temperature Range	35°C to 660°C		
ADVANCED Range			
Stability	±0.015°C @ 100°C ±0.025°C @ 650°C		
Display Resolution	0.01°C over whole range		
Accuracy: RTD Input Channel	±0.05°C ±0.005% RDG		
Accuracy: Thermocouple Input Channel	$E,J,K,N:\ \pm 0.2^\circC\ @\ 660^\circC R:\ \pm 0.6^\circC S:\ \pm 0.7^\circC\ @\ 660^\circC T\ \pm 0.2^\circC\ @\ 150^\circC$		
CJC Accuracy	±0.35°C		
BASIC/SITE Range			
Stability	±0.02°C @ 100°C ±0.03°C @ 650°C		
Display Resolution	0.01°C from 30.00 to 99.99°C then 0.1°C: 0.01°C Over PC Interface		
COMMON Specifications			
Display Accuracy ¹	0.5°C		
Blackbody Source	±0.3°C		
Surface Sensor Calibrator	±0.5°C		
Cools from 650°C to 150°C	in 60 minutes		
Heats from 30°C to 650°C	in 20 minutes		
Best Performance	See Graph		
Calibration volume	35mm diameter by 148mm deep		
Standard Insert	6 pockets, 2 x 4.5mm, 2 x 6.4mm, 1 x 8.0mm, 1 x 9.5mm diameter, all 140mm deep		
Indicator units	°С, °F, К		
Power	115Vac or 230Vac 50/60Hz 1000 Watts		
Dimensions	384H (including handle) x 212W x 312D mm		
Weight	8.5kg		

(1) Dry Block Mode only: Comparing 6.5mm Well to Display Value.

	ADVANCED	SITE	BASIC
Digital Display of Set and Nominal Block Temperature	Yes	Yes	Yes
PC Interface	Ethernet + USB Host	Serial	Serial
Test Thermostats	Yes - Two Inputs	Yes - Single Input	No
Independent Temperature Indicator for Reference Probe	Yes	Yes	No
Additional Inputs for Units Under Test	Up to 3: Two universal inputs for PRT, Thermocouple or Process inputs and a further Thermocouple input	No	No
Automatic Temperature Cycling	Yes	No	No
Data Logging	Yes - Export to USB	No	No
Offset Elimination	Yes - block can follow reference input	No	No
Choose English, French, Italian or Spanish Language	Yes - on full colour display	No	No
In Built Web Server	Yes	No	No
Tamper Proof Data	Yes - Suitable for life science, automotive and aerospace applications	No	No



UKAS Calibration available for these systems - International Traceability - Best Practice See page 14








Metal Block Bath

The Jupiter is supplied with an insert suitable for a wide range of sensors as standard.



Blackbody Source

Add the Blackbody accessory to allow calibration of infrared thermometers.

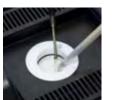

Surface Sensor Calibrator

The Jupiter can calibrate surface sensors by adding the surface sensor kit.



Standard Insert included. Size: 2 x 4.5mm, 2 x 6.4mm, 1 x 8mm and 1 x 9.5mm all 140mm deep

Alternative Inserts

Jupiter Accessories

852-09-03 Alternative Insert type B 13mm, 10mm, 8mm, 5mm and 3.5mm diameter



holes, all 140mm deep 852-09-04 Alternative Insert type C 8mm, 6 x 6.5mm diameter holes, all 140mm deep 852-07-07 Blank Insert without pockets for local machining. Includes M4 tapped hole for supplied extractor tool.
852-07-07C Custom Insert. Isotech can provide custom drilled pockets, minimum of

3mm separation between holes. Contact Isotech with your requirements.



Blackbody Kit 852-09-05 Includes a Blackbody target and Sensor.



Surface Sensor Kit 852-07-15 Includes angled thermocouple.



Calibration

Includes three point traceable calibration certificate for block temperature

UKAS Calibration Recommended: Options for block temperature and reference thermometer inputs (simulation). Legally traceable in more than 70 countries.



Air Cooling 853-04-02 For use with an air supply this accessory allows air to be blown into the block for rapid cooling.



Standard Probe 935-14-72/DB Platinum Resistance Thermometer for use up to 660°C. Probe diameter 6mm, recommended pocket size 6.5mm.



Current Loop Interface 935-06-161 24VDC Power Supply and Terminal Box. Powers 4-20mA Current Transmitters with 4mm terminal posts for easy connection.



Carrying Case 931-22-111 Sturdy case with room for accessories. Features wheels and pull out handle.





The world's leading National Metrology Institutes choose Isotech - shouldn't you?

Isotech manufacture the widest range of temperature calibration equipment from hand held thermometers to Primary Standards. With Isotech solutions you can expand your equipment no matter what the requirement.


Isotech have been pioneering the latest developments in Temperature Metrology for more than 30 years, benefit from our know how, experience and global network.

Jupiter Benefits


The Jupiter calibration block features uniform heating with a custom wound heater over an extended length of the block. The block itself is made from copper which has a very high thermal conductivity; much superior to the aluminium bronze alloys used elsewhere. A propriety process is used to protect the copper from oxidising. This combination of materials and expert knowledge delivers superior performance.

660°C Operation

The Jupiter ADVANCED operates to a maximum of 660°C; matched to the upper limit of the high temperature Isotech Semi Standard Platinum Resistance Thermometers. This allows maximum accuracy with no risk to exceeding the temperature limit of the PRT. This gives greater accuracy than extending the range beyond 660°C and having to use an inferior thermometer or thermocouple.

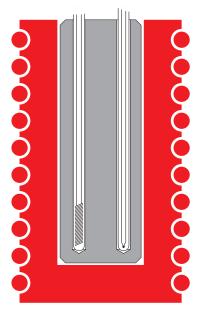


Superior uniformity by using copper block with extended length heating

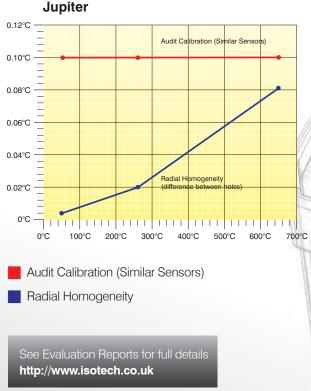


How To Order

- 1 Select Desired Options and Accessories
- 2 Supply Voltage:
 - Specify either 115Vac 50/60hz or 230Vac 50/60Hz

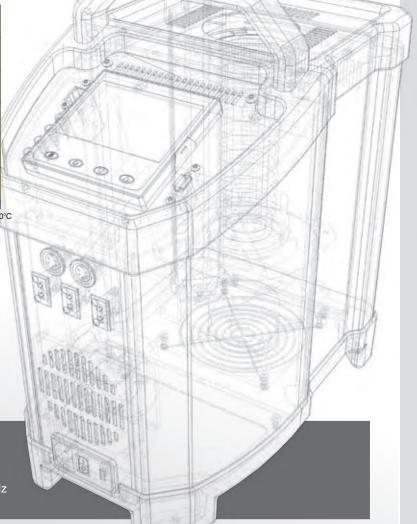



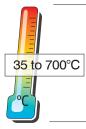

Dry Block Calibrator **Gemini**

- Calibrate Large Diameter Probes
- High Capacity Blocks 64 x 160mm
- Calibrate Whole Measurement Loop

The Gemini range of Dry Blocks have high capacity allowing a large number of probes to be calibrated together. They are also suitable to accept large diameter probes with the block volume of a nominal 64 x 160mm.

Whilst the large block takes longer to heat and cool than the Jupiter it can calibrate thermocouples, resistance thermometers, thermostats and sensors that are too large for the smaller blocks.

These award winning calibrators are easy to use and are available in three versions – the Basic, the Site and the ADVANCED. The Basic has a digital display of set and nominal temperature, the Site additionally includes an inbuilt independent temperature indicator for a reference probe. The ADVANCED controller has inputs for reference and test thermometers with a further range of sophisticated features including automatic temperature cycling, secure data logging and full colour high resolution display.

All models include I-Cal Easy LOG software and the ADVANCED models additionally include software to manage logged data and configure the unit, see page 14 for more details.

Available with a fixed block with four 8mm and four 19.5mm pockets or the LRI version which has a removable block. With the LRI model, blocks can be drilled to custom configurations.

Isotech is a world leader in temperature calibration, providing many nations with their Primary Standards and operates a full scale UKAS accredited calibration laboratory. We can offer a range of calibration options to meet your requirements.



http://www.isotechna.com

These models meet the calibration capacity requirements of EURAMET/cg-13/v.01, "EA Guidelines on the Calibration of Temperature Block Calibrators, formerly EA10/13.









Parameter	Model: Gemini 4857			
	550	700		
Temperature Range	35°C to 550°C	50°C to 700°C		
ADVANCED Range				
Absolute Stability over 30 mins	±0.01°C @ 100°C ±0.015°C	@ 300°C ±0.03°C @ 550°C		
Display Resolution	0.01°C over	whole range		
Accuracy: RTD Input Channel	±0.05°C ±0	±0.05°C ±0.005% RDG		
Accuracy: Thermocouple Input Channel	E,J,K,N: ±0.2°C @ 660°C R: ±0.6°C S: ±0.7°C @ 660°C T ±0.2°C @ 150°			
CJC Accuracy	±0.35°C			
BASIC/SITE Range				
Absolute Stability over 30 mins	±0.02°C @ 50°C ±0.03°C @ 250°C ±0.04°C @ 550°C			
Display Resolution	0.01°C from 30.00 to 99.99°C then 0.1°C: 0.01°C Over PC Interface			
COMMON Specifications				
Display Accuracy ¹	0.5°C			
Cools from 550°C to 275°C from 550°C to 60°C	35 mins (LRI: 132 mins) 345 mins (LRI: 420 mins)	-		
Heats from 30°C to 550°C from 50°C to 700°C	35 mins (LRI: 60 mins) -	- 110 mins (LRI: 120 mins)		
Best Performance	See Graph			
Calibration volume	65mm diameter	65mm diameter x 160mm deep		
Indicator units	°С, °F, К			
Voltage	115Vac or 230Vac 50/60Hz			
Power	600 Watts (LRI: 1000 Watts)			
Dimensions	384H (including handle) x 212W x 312D mm			
Weight	8.5kg	14kg		

(1) Dry Block Mode only: Comparing 6.5mm Well to Display Value.

	ADVANCED	SITE	BASIC
Digital Display of Set and Nominal Block Temperature	Yes	Yes	Yes
PC Interface	Ethernet + USB Host	Serial	Serial
Test Thermostats	Yes - Two Inputs	Yes - Single Input	No
Independent Temperature Indicator for Reference Probe	Yes	Yes	No
Additional Inputs for Units Under Test	Up to 3: Two universal inputs for PRT, Thermocouple or Process inputs and a further Thermocouple input	No	No
Automatic Temperature Cycling	Yes	No	No
Data Logging	Yes - Export to USB	No	No
Offset Elimination	Yes - block can follow reference input	No	No
Choose English, French, Italian or Spanish Language	Yes - on full colour display	No	No
In Built Web Server	Yes	No	No
Tamper Proof Data	Yes - Suitable for life science, automotive and aerospace applications	No	No



UKAS Calibration available for these systems - International Traceability - Best Practice See page 14





Gemini Accessories



Metal Block Sleeves

Gemini 550 Set of four Sleeves to suit the block. Optional single hole sizes 4, 6, 8, 10, 12, 14mm diameter all 150mm deep. 857-07-01 Undrilled sleeves for local machining. **857-07-03** 1 sleeve with 2 holes 4.5mm x 150mm deep. Gemini 700

Set of four Sleeves to suit the block. Optional single hole sizes 4, 6, 8, 10, 12, 14mm diameter all 150mm deep. 857-07-02 Blank sleeves for local machining. **857-07-04** 1 sleeve with 2 holes 4.5mm x 150mm deep.

Note: The use of sleeves will introduce an additional thermal gradient into the block. This can be avoided by using the LRI model with a block drilled for specific probes.







Gemini 700 LRI 976-07-02a Included as Standard Removable block with eight 8mm pockets 976-07-02b Blank Insert Insert without pockets for local machining 976-07-02c Custom Insert Contact Isotech with your requirements



Includes three point traceable calibration certificate for block temperature

UKAS Calibration Recommended: Options for block temperature and reference thermometer inputs (simulation). Legally traceable in more than 70 countries.



Current Loop Interface

935-06-161 24VDC Power Supply and Terminal Box. Powers 4-20mA Current Transmitters with 4mm terminal posts for easy connection.



Standard Probe

935-14-72/DB Platinum Resistance Thermometer for use up to 650°C

935-14-63 Type N Thermocouple for use up to 700°C



Carrying Case 931-22-111 - Gemini 550 / 700 931-22-112 - Gemini 550 LRI / 700 LRI Sturdy case with room for accessories. Features wheels and pull out handle.




Gemini LRI **Removable Block** Eight 8mm Pockets Can be custom drilled

Gemini

Fixed Block

Four 8mm Pockets Four 19.5mm Pockets

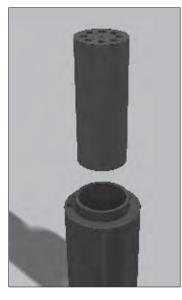





The world's leading National Metrology Institutes choose Isotech - shouldn't you?

Isotech manufacture the widest range of temperature calibration equipment from hand held thermometers to Primary Standards. With Isotech solutions you can expand your equipment no matter what the requirement.

Isotech have been pioneering the latest developments in Temperature Metrology for more than 30 years, benefit from our know how, experience and global network.


Gemini Benefits


- The Gemini has a large block with sufficient mass to accommodate larger sensors or a larger number of sensors. The fixed block has four 19.5mm pockets and four 8mm pockets.
- The LRI model with its removable block can be custom drilled. Whilst the larger block takes longer to heat than the fast response models they do allow larger probes to be accommodated and can be used with custom blocks to simulate applications, allowing probes that would otherwise be unsuitable for Dry Block calibration to be evaluated.




How To Order

- 1 Select Desired Options and Accessories
- 2 Supply Voltage: Specify either 115Vac 50/60hz or 230Vac 50/60Hz

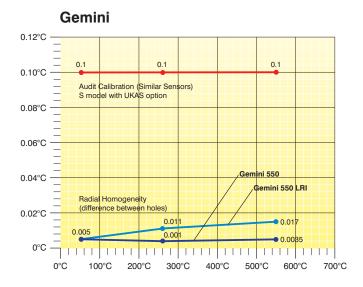

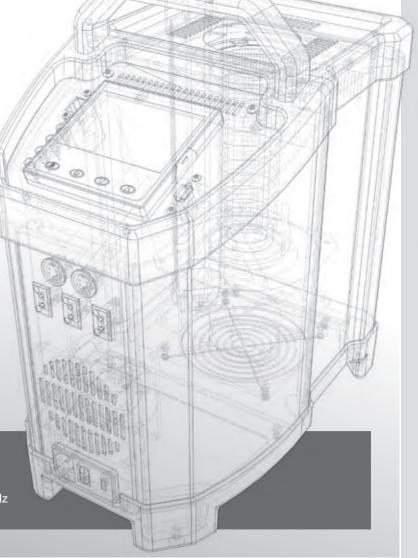



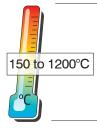

Dry Block Calibrator **Pegasus**

High Temperature Thermocouple Calibration Furnace

- Custom Furnace Design with Optimised Profile
- Calibrate Whole Measurement Loop

The Pegasus range offers extreme high temperature calibration in an easy to use portable package - ideal for the calibration of high temperature thermocouples. It has been designed for fast heating and finds applications in the glass, electrical power, automotive and material processing industries.

A Blackbody target can be added for the calibration of infrared thermometers.

The standard insert has four 8mm pockets 80mm deep. The metal insert is strategically placed beneath 50mm of insulation to provide optimal performance over the radiant temperature range.

The optional Blackbody target is used with a specially angled Type R thermocouple that sits immediately behind the target area.

These award winning calibrators are easy to use and are available in three versions – the Basic, the Site and the ADVANCED. The Basic has a digital display of set and nominal temperature, the Site additionally includes an inbuilt independent temperature indicator for a reference probe. The ADVANCED controller has inputs for reference and test thermometers with a further range of sophisticated features including automatic temperature cycling, secure data logging and full colour high resolution display.

The B model should be used with an external reference probe and indicator, such as the milliK. The thermocouples



under test should be calibrated by comparison to the external probe.

All models include I-Cal Easy LOG software and the ADVANCED models additionally include software to manage logged data and configure the unit, see page 14 for more details.



UKAS Calibration available for these systems - International Traceability - Best Practice See page 14









Parameter	Model	
	Pegasus 4853	
Temperature Range	150°C to 1200°C	
ADVANCED Range		
Stability	±0.05°C @ 150°C ±0.08°C @ 1200°C	
Display Resolution	0.01°C over whole range	
Input Channel Accuracy: Thermocouple	${\sf E}, {\sf J}, {\sf K}, {\sf N}: \ \pm 0.2^{\circ} {\sf C} \ @ \ 660^{\circ} {\sf C} {\sf R}: \ \pm 0.6^{\circ} {\sf C} {\sf S}: \ \pm 0.7^{\circ} {\sf C} \ @ \ 660^{\circ} {\sf C} {\sf T} \ \pm 0.2^{\circ} {\sf C} \ @ \ 150^{\circ} {\sf C}$	
CJC Accuracy	±0.35°C	
Input Channel Accuracy: RTD	±0.05°C ±0.005% RDG	
BASIC / SITE Range		
Stability	±0.1°C @ 150°C ±0.2°C @ 1200°C	
Display Resolution	0.1°C from 150°C to 999.9°C then 1°C: 0.01°C Over PC Interface	
COMMON Specifications		
Blackbody Source	±0.3°C	
Cools from 1200°C to 800°C 1200°C to 200°C	in 50 minutes* in 180 minutes* *substantially reduced by the cooling adaptor	
Heating Rate	25°C / minute	
Best Performance	See Graph	
Calibration volume	33.5mm diameter by 130mm deep	
Standard Insert	4 x 8mm Pockets all 80mm deep + 50mm top insulator	
Indicator units	°C, °F, K	
Power	115Vac or 230Vac (50 / 60 Hz) 800 Watts	
Dimensions	384H (including handle) x 212W x 312D mm	
Weight	13kg	

	ADVANCED	SITE	BASIC
		SQ1	
Digital Display of Set and Nominal Block Temperature	Yes	Yes	Yes
PC Interface	Ethernet + USB Host	Serial	Serial
Test Thermostats	Yes - Two Inputs	Yes - Single Input	No
Independent Temperature Indicator for Reference Probe	Yes	Yes	No
Additional Inputs for Units Under Test	Up to 3: Two universal inputs for PRT, Thermocouple or Process inputs and a further Thermocouple input	No	No
Automatic Temperature Cycling	Yes	No	No
Data Logging	Yes - Export to USB	No	No
Offset Elimination	Yes - block can follow reference input	No	No
Choose English, French, Italian or Spanish Language	Yes - on full colour display	No	No
In Built Web Server	Yes	No	No
Tamper Proof Data	Yes - Suitable for life science, automotive and aerospace applications	No	No







Metal Block Bath

The Pegasus includes an insert suitable for high temperature calibration of thermocouples.



Blackbody Source

Add the Blackbody accessory to allow calibration of infrared thermometers.

Pegasus Accessories



Metal Block Insert

Standard Insert Included Four 8mm pockets. Pocket depth 80mm + 50mm insulator. Effective depth 130mm. 853-06-02 Blank Insert Insert without pockets for local machining 853-06-02b Custom Insert Contact Isotech with your requirements



Blackbody Kit 853-06-03 Includes a Blackbody target and Sensor



Calibration

Includes three point traceable calibration certificate for block temperature

UKAS Calibration Recommended: Options for block temperature and reference thermometer inputs (simulation). Legally traceable in more than 70 countries.



Standard Probe 935-14-91 Type R Platinum Thermocouple for use up to 1200°C.



Air Cooling 853-04-02 For use with a compressor this accessory allows air to be blown into the block for rapid cooling.



Ceramic Insulators 853-06-04 Spare insulation pack Includes 2 x standard tops and 2 x standard bottoms.



Current Loop Interface 935-06-161 24VDC Power Supply and Terminal Box. Powers 4-20mA Current Transmitters with 4mm terminal posts for easy connection.



Carrying Case 931-22-111 Sturdy case with room for accessories. Features wheels and pull out handle.





The world's leading National Metrology Institutes choose Isotech - shouldn't you?

Pegasus Benefits

The Pegasus features a small tube furnace to allow operation to 1200°C in a portable case. With a ceramic furnace construction temperature gradients are larger than with lower temperature metal blocks.


The Isotech furnace benefits by using a specially wound furnace tube assembly. They are manufactured in our factory with the turns concentrated at the ends of the furnace, where the heat losses are greatest. This gives an

improved temperature profile and lower uncertainty.

Insulators are provided for the top and bottom of the furnace which further improve temperature uniformity.

The effective immersion depth is 130mm, 80mm in the metal insert and then a further 50mm in the furnace tube.

Benefit from Isotech's design and experience

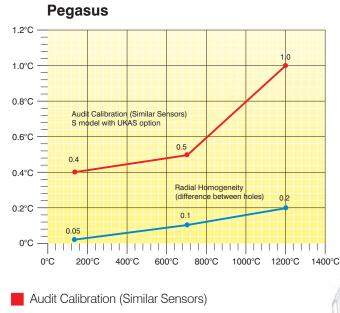


Radial Homogeneity

See Evaluation Reports for full details http://www.isotech.co.uk

How To Order

- 1 Select Desired Options and Accessories
- 2 Supply Voltage:
 - Specify either 115Vac 50/60hz or 230Vac 50/60Hz

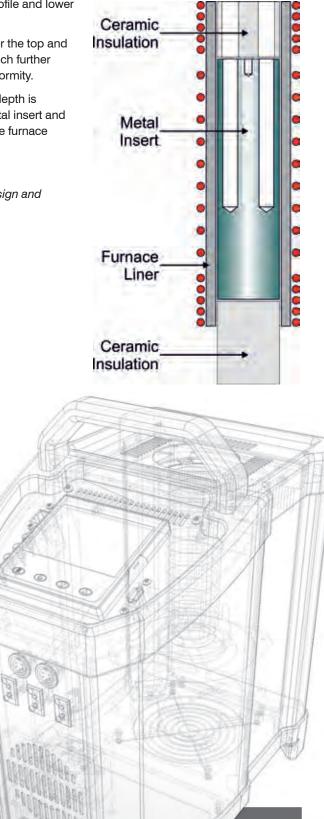



An Introduction to **Fast Calibrators**

This section focuses on the equipment needed for the rapid checking, testing and calibration of instrumentation and temperature sensors.

Service engineers and those working on site will appreciate the benefits of simple and fast temperature calibration. An engineer forced to carry a calibrator up a ladder or into a confined space will value the handheld QuickCal.

The Products Featured in this section have:

- Outstanding Value
- Compact Size with true handheld models
- Wide Operating Ranges
- Fast Response

Quick-Cals

There are two Quick-Cal models, handheld, portable and capable of operating from -12°C to 350° C

Fast-Cals

Fast-Cals work from -35°C to 650°C in three ranges, -35°C to 140°C, 30°C to 350°C and 35°C to 650°C.

During 2004, 20 experienced engineers from many parts of the world specified their ideal products for Industrial Calibration. Fast-Cal realizes their top ten requirements of:

- 1 Rugged
- 2 Lightweight
- 3 Easy to use on site
- 4 Low cost/high benefit ratio,
- 5 Fast response, high stability
- 6 Time saving features
- 7 Multiple sensor testing
- 8 Software
- 9 Modern design
- **10** Compliant with latest regulations.


One model is ideal for the Validation of Washer Disinfectors, Steam Sterilisers and Autoclaves. In place of a removable insert it has a fixed block with pockets for a reference probe and the type of test sensor commonly used in validation applications.



Immersion depth is very important

When selecting a Dry Block Calibrator depth of immersion is very important. The chart below provides guidance in selecting a bath for immersion depth. Note that sensors with

a long sensing length will require greater immersion. The chart is general and applies equally to all dry blocks - not to a particular model or manufacturer.

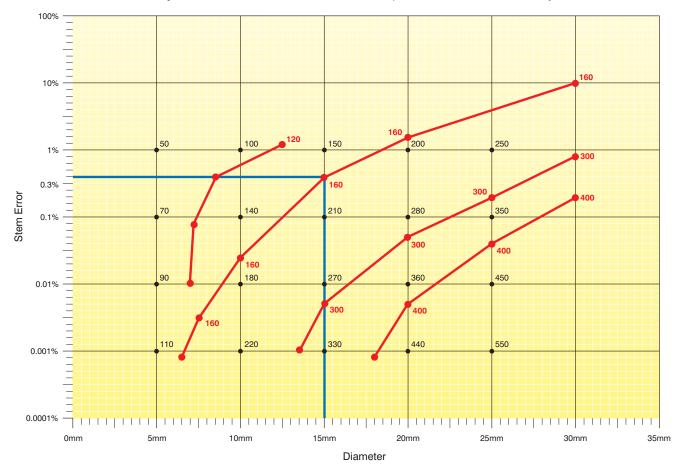

Immersion Depths for various diameter thermocouples or thermistors in a dry block bath

Example shows 0.3% stem error for a 15mm diameter thermocouple immersed 160mm in a dry block.

Note 1 For sensors immersed in stirred liquids the diameter of sensor can be doubled, or the minimum depth halved.

- **Note 2** The sensing length must be added to the above immersion depth calculation
- **N.B.** The above gives a good guide, however each sensor will be slightly different.

For full information on immersion depth visit the technical library on the Isotech website.

http://www.isotechna.com





Temperature Test Unit Quick-Cal

Fast
 Compact
 Two models -12°C to 350°C

Isotech's Quick Cal range of testers are ideally suited to the less exacting applications where our larger, better specified baths are unnecessary.

The Low Temperature Quick Cal covers the temperature range of -12° C to $+140^{\circ}$ C, making it ideal for the medical, sterilisation and food industries.

The High Temperature Quick Cal covers the temperature range from 30°C to 350°C.

Both have interchangeable inserts, so you don't have to buy another calibrator each time you find a new diameter sensor to calibrate. The Low Temperature Quick Cal is particularly exciting since it can work from a small 12 volt battery, making it truly portable, alternatively a mains supply is available with 100V - 240V, 50/60Hz input.

The Low Temperature Quick Cal comes complete with two inserts, one blank and one drilled 6.8mm diameter 120mm deep. Alternatively we can provide pre-drilled inserts to suit your needs. Each calibration well is 13mm diameter and 120mm deep.

The High Temperature Quick Cal has a single larger well and interchangeable inserts, 25mm in diameter by 120mm deep. Blank and other special drillings are available to order. The High Temperature Quick Cal is supplied 100-130 or 200-250V AC, 50/60Hz.

The manual which accompanies the Quick Cal contains helpful advice on how to get the best from your calibrator. Accessories for the above include a carry case, alternative inserts and a rechargeable battery for the Low Temperature Quick Cal.



Quick-Cal High

PERATURE RANGE 30-35

UTCK-CAL

(E

DANGER HOT



Specification		Specification		
Model	Low Temperature Quick-Cal	Model	High Temperature Quick-Cal	
Temperature Range	-12°C to +140°C In an ambient of 20°C The Quick-Cal will achieve lower temperature if used in a lower	Temperature Range	+30°C to +350°C In an ambient of 20°C	
		Set Point Resolution	0.1°C over range	
	ambient temperature	Accuracy	±0.1 to ±0.4°C using Comparison Techniques	
Set Point Resolution	0.1°C over range		$<100^{\circ}C \pm 0.4^{\circ}C$ against display	
Accuracy	± 0.1 to $\pm 0.4^{\circ}$ C using Comparison Techniques $\pm 0.4^{\circ}$ C against display temperature with a single probe in the 6.8mm drilled sleeve		temperature with a single probe in the 4.5mm pocket >100°C \pm 1.5°C against display temperature with a single probe in the 4.5mm pocket	
Stability	±0.05°C	Stability	±0.05°C	
Time to Set Point	9 minutes from 0°C to 100°C	Time to Set Point	9 minutes from Ambient to 350°C	
Additional Time for best stability	Typically 5 minutes	Additional Time for best stability	Typically 5 minutes	
Calibration volume	Two 13mm diameter pockets 120mm deep. Supplied with two sleeves, one blank and one drilled 6.8mm diameter 120mm deep	Calibration volume	25mm diameter 120mm deep Supplied with insert, 5 pockets, 8mm, 6.5mm, 6.5mm, 4.5mm, 4.5mm, all 115mm deep	
Voltage	15Vdc or 100-240V, Switch mode power supply (50 / 60 Hz)	Voltage	100-130V or 200-240V, (50 / 60 Hz)	
Power	70 Watts	Power	350 Watts	
Dimensions	Height 65mm Width 152mm Depth 175mm Weight 1.5kg	Dimensions	Height 65mm Width 152mm Depth 175mm Weight 1.5kg	

How To Order

Specify Model, Low or High Temp, Supply Voltage, Accessories and if UKAS Calibration is required.

Extra Sleeves

560-06-01A Standard Sleeve
560-06-01B Blank Sleeve
560-06-01C Special Sleeve (maximum bore 10mm)



Carrying Case

931-22-100



UKAS Calibration

Fahrenheit Operation

UKAS Calibration available to order, legally traceable in more than 70 countries.





Rechargeable Battery & Charger Battery is 12V 12Ah.

Unit can be supplied for °F Operation

On full charge battery: 60 minutes when cooling, 90 minutes when heating

Longer life available with larger batteries - consult lsotech

Temperature Range is -10°C to 123°C with battery operation





931-22-100

Extra Inserts

550-06-01B



UKAS Calibration UKAS Calibration available to order,

550-06-01A Standard Insert

550-06-01C Special Insert

Blank Insert

UKAS Calibration available to order, legally traceable in more than 70 countries.








Semi Standard Platinum Resistance Thermometer

935-14-98

Note: Also suits Quick-Cal Low-temp



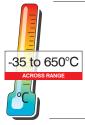

Industrial Temperature Calibrators Fast-Cal

- 25 x 148mm Calibration Volume
- Rugged
- Lightweight
- Easy to use on Site

Fast-Cal - the product range you designed !

During early 2004, 20 experienced engineers from many parts of the world were asked to produce their ideal specification for industrial calibration. Their top ten requirements, in order of priority, are: *Rugged, Lightweight, Easy to use on site, Low cost/high benefit ratio, Fast response, high stability, Time saving features, Multiple sensor testing, Software, Modern design, Compliant with the latest regulations.*

One model is ideal for the validation of Washer Disinfectors, Steam Sterilisers and Autoclaves. In place of a removable insert it has a fixed block with pockets for a reference probe and the type of test sensor commonly used in validation applications.

Fast-Cals are available in Basic or Complete models


A Basic Dry Block Calibrator

The thermometer under test is compared to the dry block controller value. Useful for moderate temperature ranges and quick testing.



An ISO 9000 Calibration System

A thermometer under test is compared to a calibrated standard, for true traceability and clearly meets the requirements of ISO9000. Complete models include a semi-standard probe.

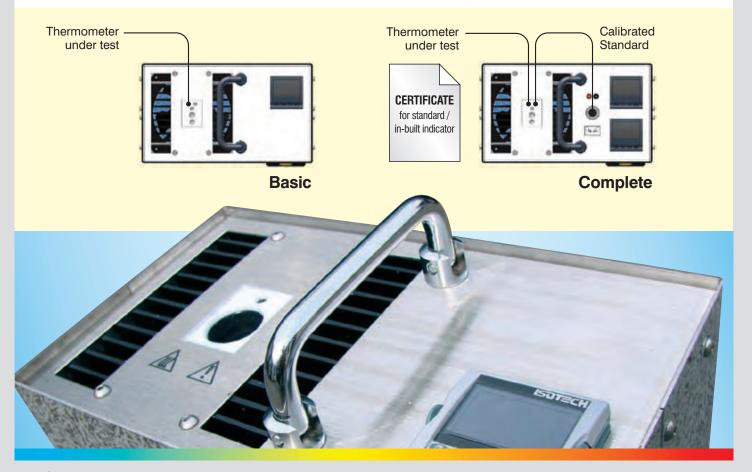



Model No.	HTM 2010	Low	Medium	High
Range at 20°C	-35 to +140°C (1)	-35 to +140°C (1)	30 to 350°C	35 to 650°C
Stability	0.02°C	0.02°C	0.03°C	0.03 to 0.05°C
Accuracy - Basic	±0.2°C	±0.2°C	±0.3°C	±1°C (500°C) ±2°C (650°C)
Accuracy - Complete (2)	±0.15°C	±0.15°C	±0.2°C	±0.5°C
Heating Time	-30 to +140°C In 15 minutes	-30 to +140°C In 15 minutes	50 to 350°C In 15 minutes	50 to 650°C In 20 minutes
Cooling Time	+140 to 0°C In 15 minutes	+140 to 0°C In 15 minutes	350 to 100°C In 40 minutes	650 to 300°C In 20 minutes
Calibration Capacity	145mm depth 8Ø, 6.5Ø, 3 x 4.5Ø	148mm depth 25mm Ø	148mm depth 25mm Ø	148mm depth 25mm Ø
Power	150W	150W	750W	750W
Traceable Certificate	Included	Included	Included	Included
UKAS Certificate	Extra	Extra	Extra	Extra
PC Communications	Included	Included	Included	Included
Software	Included	Included	Included	Included
4-20mA Input	Complete Models	Complete Models	Complete Models	Complete Models
Ramps and Dwells	Yes with Software	Yes with Software	Yes with Software	Yes with Software
Dimensions	228 x 248 x 143mm	228 x 248 x 143mm	228 x 248 x 143mm	228 x 248 x 143mm
Weight	6.60kg	6.60kg	6.35kg	6.35kg
(1) In ampliant of 20°C: Minimum To	non a materia in 550 Dalana Amatian	+ Abaabata Minimuma 4500		

(1) In ambient of 20°C: Minimum Temperature is 55°C Below Ambient, Absolute Minimum -45°C.

How To Order

Specify Model, Basic or Complete, Supply Voltage, Accessories and if UKAS Calibration is required.

Features (Basic & Complete)	HTM 2010	LOW	MEDIUM	HIGH
Rugged, hard-wearing stainless steel case	✓	✓	\checkmark	✓
Hygienic with no paint to chip or flake	\checkmark	\checkmark	\checkmark	\checkmark
No membrane keypads	✓	\checkmark	\checkmark	\checkmark
Lightweight - less than 6.6kg	✓	\checkmark	\checkmark	\checkmark
Thermal Switch Test - Complete Model with Software	✓	\checkmark	\checkmark	\checkmark
Easy to use on site, effortless transporting	✓	\checkmark	\checkmark	\checkmark
Fast response - high stability	✓	\checkmark	\checkmark	\checkmark
Built-in Insert holder for ultra fast changes and restarts	-	-	\checkmark	\checkmark
Removable inserts, drilled to your requirements	-	\checkmark	\checkmark	\checkmark
PC interface	\checkmark	\checkmark	\checkmark	\checkmark
Additional software available for automatic calibration	\checkmark	\checkmark	\checkmark	\checkmark
Certificate traceable to National Standards	✓	\checkmark	\checkmark	\checkmark
Five point UKAS Certification available	✓	✓	\checkmark	\checkmark
Mains power variance immunity	✓	\checkmark	\checkmark	✓

Extra Inserts

Blank

Note: The Complete model can capture a thermostat trip temperature either as a stand alone unit or with a PC and the included software.



Low, Medium & High

Standard Insert is:

 $2 \times 4.5 \text{mm} + 1 \times 6.5 \text{mm} +$ 1 x 8mm x 140mm Deep. All Inserts have a 4mm tapped

hole to suit supplied extractor tool.

HTM2010 Fixed Block

3 x 4.5mm + 1 x 6.5mm + 1 x 8mm x 145mm Deep.

	Low & Medium	Hi
Standard	907-02-03	90
Custom	907-02-03c	90

In addition to the standard insert, specially

drilled or blank inserts are also available.


907-02-03b





Carrying Case included



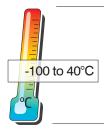

-100°C Dry Block Isis

- 35 x 160mm Calibration Volume
- 20,000 hours testing, equal to 10 years use
- Minimum Temperature -100°C (-148°F)
- No expensive liquids

The Isis Dry Block offers operation to temperatures as low as -100°C, and is the only block bath working to such a low temperature. Now it is possible to calibrate temperature sensors such as PRTs, Thermocouples and Thermistors at ultra low temperatures without the need for a liquid bath.

Portability and Safety

Unlike a liquid bath the Isis requires no costly, or hazardous fluids and offers greater portability. This will be of particular value to calibration engineers working on site with low temperature freezers as encountered in pharmaceutical, aeronautical and food environments.

The minimum operating temperature is less than stirred liquid laboratory calibration baths and users in laboratories will also benefit by avoiding the ongoing need for expensive fluids.

The maximum operating temperature is 40°C, a little higher than the minimum operating temperature of Isotech Hot Blocks. This permits covering the range from -100°C to 650°C or higher with just two Isotech blocks. By limiting the maximum temperature the reliability and operating life of the cooling engine is maximised, and has been agreed in consultation with the licensee of the cooling technology, see below.

Cooling Technology

The Isis makes use of a Free Piston Stirling Cooler (FPSC) which provides a massive 80 Watts of cooling power to the calibration block. Specialist materials, patent applied for, are used for the heat transfer from the FPSC to the block.

Operating Life

Reliability is a prime attribute of this revolutionary new product. Testing at 20,000 hours (nominally equivalent to 10 years at 40 hours use each week) shows that -100°C is still possible, with an increase in cooling time <10%.



Benefits

Isotech can offer full support with options for UKAS / ILAC calibration, tutorial on getting the best calibration uncertainties and a full range of supporting reference thermometers, indicators and software.

The Isis has a large insert 35mm diameter by 160mm deep. This allows for calibration of multiple sensors. For thermal validation applications there is an insert with pockets for a reference probe (6.5mm) and 20 x 3.5mm pockets for thermocouples. This allows a single calibration cycle to validate up to 20 probes at a time.







METAL BLOCK INSERTS

Standard Insert A



9.5mm, 8.0mm, 6.4mm, 6.4mm, 4.5mm, 4.5mm All 157 deep

ALTERNATIVE INSERTS

Standard Insert B



13mm, 10mm, 8mm, 5mm and 3.5mm All 157 deep

Standard Insert C



8mm and 6 x 6.5mm All 157 deep

Validation Insert



6.5mm and 20 x 3.5mm All 157 deep

Model	525 Isis	
Temperature Range	-100°C to +40°C	
Approximate time to Temperature from Ambient	-40°C 3 -60°C 4 -80°C 6	20 minutes 30 minutes 40 minutes 60 minutes 90 minutes
Absolute Stability	At 0°C ±0.03°C (30 m At -90°C ±0.02°C (30	
Radial Homogeneity (similar pockets)	0.01°C	
Vertical gradients (over bottom 40mm)	0°C 0.1°C -90°C 0.2°C	
Calibration Volume	35mm diameter x 160	mm deep (Excludes Insulating Cap)
Standard Insert	6 thermometer wells a All 157 deep)	as standard (9.5mm, 8.0mm, 6.4mm, 6.4mm, 4.5mm and 4.5mm
Power	200W	
Voltage	100-240Vac, 50/60Hz	
Dimensions	215mm(W) x 420mm	(D) x 640mm (H)
Weight	22.7kg	

CALIBRATION SOLUTIONS 43



Introducing Absolute Calibration with ITS-90 Fixed Point Cells

For calibration to the smallest of uncertainties thermometers are calibrated by placing them into a series of Fixed Point Cells. For example pure aluminum freezes at 660.323°C so by first melting a cell containing pure aluminum, then placing a thermometer into it as the metal changes state, from a liquid to a solid, a very precise calibration point is realized.

This absolute or fixed point calibration is performed by National Metrology Institutes providing primary standards and directly realizing the International Temperature Scale, ITS-90. Isotech's solutions for Primary Standards are found in a separate publication, *"Volume 1: Solutions for Primary & Secondary Laboratories."*

Isotech also offer a range of ITS-90 Fixed Point systems that are less expensive, easier to use and more robust than the larger cells used by the international NMIs.

For some countries, where the local industry needs are less demanding Slim Cells are used by NMIs and Isotech can offer UKAS calibration with uncertainties from 0.5mk to 2mK over the range -38°C to 660°C.

Users in industrial and secondary laboratories benefit from using Slim Cells to calibrate to smaller uncertainties than is possible with dry blocks or liquid baths. The Isotech Slim Water Triple Point Cell is comparable in cost to a specially drilled metal insert, putting it in the reach of all calibration engineers. Using a Water Triple Point Cell allows cost effective checking of standards between calibrations, and to help determine when a thermometer needs recalibration. Water triple point cells have uncertainties less than 0.001°C at a very modest cost.

In order to use an ITS-90 Fixed Point Cell, apparatus is needed, it must create a zone of constant temperature around the cell so that the cell can melt of freeze uniformly. Isotech equipment uses multi zone heating or for optimal performance a heat pipe or heat siphon. To calibrate a thermometer it must be sufficiently immersed that further immersion would make no temperature change to the thermometer.

A new innovation from Isotech is the Isothermal Tower, which combines apparatus, a heat siphon, fixed point cell and an immersion compensation device.* The Isothermal Towers are simple to use integrated devices providing optimal performance.

There are also ranges of apparatus that can accept a range of cells, see table opposite. These models can also be used without cells, including use as Dry Blocks

Fixed Point State Temperature °C Argon **Triple Point** -189.3442 Mercury **Triple Point** -38.8344 Water **Triple Point** 0.010 Gallium Melt Point 29.7646 Indium Freeze Point 156.5985 Tin Freeze Point 231.928 Zinc Freeze Point 419.527 Aluminium Freeze Point 660.323 Silver Freeze Point 961.78

for immersion depths of up to 300mm

- ideal for larger sensors.

* Patents applied for



Equipment

The ITS-90 Cell needs equipment in order to melt, freeze or maintain it. Many combinations of Isotech equipment can be used; liquid baths, dry blocks and furnaces. Isotech have a range of Dry Blocks that allow the fixed points to be easily realized at an affordable price. The combination allows you to calibrate at a "Point On The Temperature Scale" hence the name POTTS.

Point	Temperature	Suitable Apparatus	
Mercury	-38.8344°C	Europa	P
B8 Water Triple Point	0.01°C	Europa Venus	
B12 Water Triple Point	0.01°C	Oceanus	
Gallium	29.7646°C	Europa Venus Calisto	
Indium	156.5985°C	Medusa 510 Medusa 511 ISOTower 490	
Tin	213.928°C	Medusa 510 Medusa 511 ISOTower 491	
Zinc	419.527°C	Medusa 510 Medusa 511 ISOTower 492	
Aluminium	660.323°C	Medusa 510 Medusa 511 ISOTower 493	
Silver	961.78°C	Oberon	
Copper	1084.62°C	Oberon	

ISOTowers are intergrated devices including the ITS-90 Fixed Point Cell





ITS-90 Isothermal Towers ISOTower

- High accuracy
- Fast to temperature
- Simple to use
- Patented

The most accurately defined temperatures are those defining points (fixed points) of ITS-90.

The leading Primary Laboratories use large fixed point cells in deep calibration furnaces that utilise Heat Pipes to eliminate temperature gradients. This combination of cell and furnace gives the smallest of uncertainties.

In the patented Isothermal Towers the fixed point cell and heat pipe (or heat siphon) have been combined to produce the ideal realisations for calibrating standard thermometers.

Thermometers can only be calibrated accurately if they are immersed sufficiently.

In Isothermal Towers a heated block (Immersion Compensator, patent applied for) sits on top the heat siphon/cell to fully compensate for the immersion characteristics of the unit under test.

The Isothermal Towers performance has been fully evaluated against the most detailed and demanding requirements ever written: CCT/2000-13.

All Isothermal Towers; Indium, Tin, Zinc and Aluminium meet all the requirements of CCT/2000-13 allowing laboratories to realise the smallest uncertainties, at a fraction of the cost of conventional Metrology Furnaces with Primary Standard Cells

You can purchase three Isothermal Towers; Tin, Zinc and Aluminium for a similar price as one conventional cell and heat pipe apparatus!



Isothermal Towers are simple to use, and very robust. Operation is risk free, as a combined apparatus there is no need to handle a fragile cell. No need for specialist training courses. Isothermal Towers remove the mystery from fixed point calibration.

Easily set to provide a melt or freeze of 24 hours or more, lending themselves for automatic calibration and providing your lab with an all day long plateau.





Perfect Audit Item

As an audit item, an accreditation authority can send the device to laboratories for intercomparison.


Because the cell, apparatus and immersion compensator are a single entity, the performance is unambiguous unlike existing systems where cell and apparatus are often separated during intercomparison. Accreditation authorities love them.

Available to Hire

Additionally Isothermal Towers are available to hire from Isotech and a growing number of Isotech Distributors to allow laboratories to audit themselves by intercomparing their cells and standard thermometers to a UKAS calibrated Isothermal Tower.

Transportable


ITS-90 Isothermal Towers are transportable by carrier; there are no fragile glass parts!

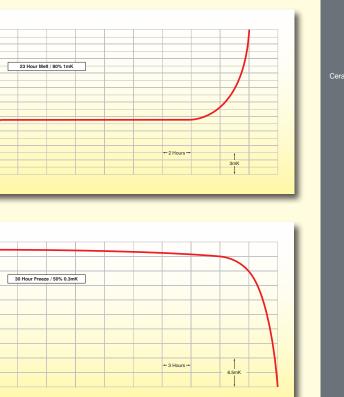

■ Isothermal Towers include a traceable calibration certificate.

This includes a graph of one freeze, one melt plateau and a certificate of purity for the metal inside the siphonic cell. As an option, UKAS calibration is available to one of two services, see table over for the uncertainties.

Full data available at www.isotech.co.uk/isotower

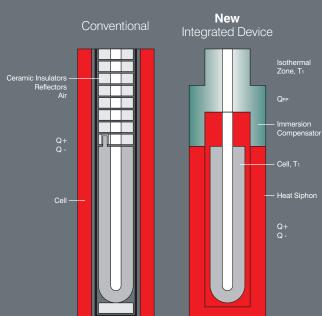

A fixed point cell is not long enough to eliminate heat conductance along the thermometer calibrated in it. Currently, using long furnaces, heat shunts and reflective baffles an attempt is made to reduce these losses.

The ISOTower uses a combined metal clad fixed point cell and heat siphon, which when heated provides an isothermal environment for the metal within to change state. The outer wall of the cell becomes the inner wall of the heat siphon with cost as well as performance benefits.

Additionally an Immersion Compensator is used to compensate for the stem conduction problems caused when a thermometer under test is not sufficiently immersed into a fixed point cell.

Benefits of the ISOTower over a conventional Quartz Cell and Apparatus

ISOTower

- Robust no glass parts
- Easily Transported
- Integrated Device known immersion characteristics
- Uniquely integrated cell, apparatus and correction for thermometer stem conduction
- Simple and safe to use with increased confidence in results

Conventional Quartz Cell and Apparatus

- Fragile and Risk of Breakage
- Difficult and expensive to Transport
- Cell certified separate from apparatus, stem conduction unknown

ISOTECH

Specification

Model	490	
ITS-90 Point	Indium	
Temperature	156.5985°C	
Metal Purity	6N	
Plateau Duration		
UKAS Uncertainty: Premium Service*	±0.7mK	
UKAS Uncertainty: Standard Service*	±2mK	
Heating Time	2 hrs	
Pocket Diameter		
Total Immersion Depth		
Depth of metal surface to bottom of ree	entrant tube	
PC Interface		S
Power		
Voltage		
Dimensions		
Weight		

* UKAS Calibration is Optional, Uncertainties apply to whole system

491	492	493
Tin	Zinc	Aluminium
231.928°C	419.527°C	660.323°C
6N	6N	6N
Up to	30 hrs	
±0.8mK	±1mK	±2mK
±2mK	±2mK	±6mK
2 hrs	2 hrs	2 hrs
8 r	nm	
290	mm	
180	mm	
Supplied with PC (Cable and Software	
900 \	Natts	
110 Vac or 230 Vac 50/60Hz		
H 430 mm x W 31	0 mm x D 300 mm	
15	kg	

The ISOTowers are protected by two patents.

Immersion Compensator


Fully compensates for the immersion characteristics of the thermometer under test.

Siphonic Cell

A combined Fixed Point Cell and Heat Pipe

The outer wall of the fixed point cell is the inner wall of the heat pipe or siphon, thus integrating the cell and its apparatus, giving a guaranteed performance from the cell.



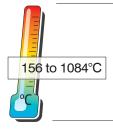

Slim Fixed Point Cells Sealed

■ Ultra Pure >99.9999% 6N

- 35 Year Plus History
- For Optimal Realisations

Water

The Isotech B8 30 130 Cell is small enough to fit into portable Dry Blocks such as the Venus and Europa models. For the larger blocks like the Oceanus Isotech recommend the B 12 46 210 with its increased immersion depth

Gallium and Mercury Cells

Like the B8 Cell the Slim Gallium Cell can be used in portable blocks like the Venus and Europa, or in stirred liquid baths. The Hydra and Orion have accessories available to support the cells.

Indium to Copper

Isotech's Slim Cells have been in constant use since their introduction in 1990. The cells have always been made from the highest quality graphite and 6N (99.9999%) pure metals.

After further investment in the lab, and gaining smaller uncertainties from UKAS, we reviewed and further refined our range of metal clad cells to give better accuracy and performance. The new professional ranges of cells have more metal inside providing an active immersion depth in the metal of 160mm.

Metal Clad

Isotech produced the first metal clad cells in 1990 and have much experience in the manufacture and calibration of high quality proven metal clad cells.

Metal clad cells are recommend for all points from Indium to Aluminium.

Quartz Clad

These cells are recommended for Silver and Copper points, whilst available for the lower temperature points the metal clad versions are recommended as they are more robust, have the same performance and are more cost effective.

Equipment for Slim Cells

The Slim Cells can be used in the same apparatus as the larger cells, and the greater immersion depth will give the lowest uncertainties.



Cost effective dedicated desktop apparatus like the POTTS, "Points on the Temperature Scale" can be used to automatically bring the cell to the plateau. These simple to use systems conveniently provide long flat plateau for low uncertainty calibration of thermometers.

UKAS certification of our Slim Cells

The metal clad fixed point cells are intercompared to our reference cells for smallest uncertainties. Isotech now offer two UKAS services depending on the amount of measurements we make on the cell under test.

In our standard and recommended service we perform one melt, one freeze and one intercomparison. In our premium service, in order to reduce uncertainties we perform two or three melts, two or three freezes and two intercomparisons.

The two optional UKAS services with the uncertainties are tabulated below:

Isotech UKAS Calibration Uncertainties (k=2)

Cell	Premium Service UKAS Schedule Note 4	Standard Service UKAS Schedule Note 5	
Mercury	±0.5mK	±1mK	
Water	±0.1mK (B12)	±0.5mK	
Gallium	±0.5mK	±1mK	
Indium	±0.7mK	±2mK	
Tin	±0.8mK	±2mK	The la
Zinc	±1mK	±2mK	be found on th
Aluminium	±2mK	±6mK	01



be found on the Isotech website or at www.ukas.org



Available Types			Additional	Model	
Cells	Temperature	Uncertainty ¹	Uncertainty ²	Metal Clad	Model
Water	0.01°C	±0.5mk	±0.3mk	N/A	B8 30 130
Water	0.01°C	±0.1mk	±0.3mk	N/A	B12 40 210
Water	0.01°C	±0.1mk	±0.3mk	N/A	B12 46 210
Gallium	29.7646°C	±0.5mk	±0.3mk	17401M	N/A
Mercury	-38.8344°C	±0.5mk	±0.1mk	17724M	N/A
Higher Temperature					Quartz Clad
Indium	156.5985°C	±0.7mk	±0.7mk	17668ML	17668QS
Tin	231.928°C	±0.8mk	±0.8mk	17669ML	17669QS
Zinc	419.527°C	±1mk	±1.5mk	17671ML	17671QS
Aluminium	660.323°C	±2mk	±3mk	17672ML	17672QS
Silver	961.78°C	±30mk		N/A	17673QS
Copper				N/A	17674QS

Isotech cells are of the highest purity available. Open cells conform to CCT/2000-13. Sealed cells are sealed to one atmosphere with 6N pure argon at the freeze temperature.




A free report is available, www.isotech.co.uk/pdfs/SlimCells.pdf

N.B. Other SPRTs may give different results depending on the stem conduction properties.

1. The uncertainty applies when the cells are sufficiently immersed in deep apparatus.

2. When these cells are used in bench-top apparatus the additional uncertainty should be included for stem conduction effects. The value is typical for the 670 SPRT, others SPRTs may give different results depending on the stem conduction properties.



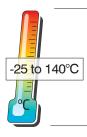

ISOCAL - 6 Range Oceanus-6

- 50 x 300mm Calibration Volume
- Use with Water and Gallium Fixed Points
- Can also be used as a Dry Block, Liquid Bath...
- -25°C to 140°C

The Oceanus-6^{PLUS} has all the advantages of the Isocal-6 models but with a substantially larger calibration volume, 50mm diameter by 300mm deep. The Oceanus-6 can be used as a Dry Block, a Liquid Bath, a Blackbody Source for infrared thermometers, a Surface Sensor Calibrator and for performance to a few mK (0.001°C) ITS-90 Fixed Points. The Oceanus-6^{PLUS} is available in two models, the BASIC (B) and the SITE (S). The B model includes a sophisticated temperature controller with a dual display for Set Temperature and Dry Block Temperature.

The S model includes a built in digital thermometer to which an external standard thermometer can be connected, for Dry Block use this will give greater accuracy eliminating temperature gradient and loading errors.

For Liquid Bath, Surface Sensor or Blackbody use an external thermometer should always be used - either with the S model or the B model and a separate stand-alone indicator. For Lab use the Oceanus- 6^{PLUS} can be used with a laboratory performance temperature indicator such as one of the Isotech True Temperature Indicators with performance, for similar sensors, down to hundredths of a degree. The Oceanus- 6^{PLUS} offers unprecedented accuracies of $\pm 0.0002^{\circ}$ C (2 Sigma) at the Water triple point and the Gallium melt temperature of 29.7646°C and up to $\pm 0.005^{\circ}$ C in the stirred liquid bath option (by comparison).








Specification			
Model Temperature Range	Oceanus-6 45°C below ambient to +140°C	Display Resolution	0.01 -19.99 to 99.99 0.1 100.0 to 140.0 PC can display 0.01 across
Absolute stability over 30 minutes	Dry Block Bath ±0.03°C Stirred Liquid Bath ±0.025°C		whole range with the software included
	Ice/Water Bath ±0.001°C Blackbody Source ±0.3°C	Indicator units	°C, °F, K
	Surface Sensor Calibrator ±0.5°C ITS-90 Fixed Point ±0.0002°C	Power	108 to 130V (50 / 60 Hz) or 208 to 240V (50 / 60 Hz) 300 Watts
Computer Interface	Included with Software	D	
Cools from	20°C to -10°C in 90 minutes	Dimensions	Height 430mm Width 310mm
Heats from	-10°C to 80°C in 60 minutes		Depth 300mm
Calibration volume	50mm diameter by 300mm deep	Weight	17kg
Standard Insert	6 pockets, all 8.0mm diameter and 250mm deep		

Isocal-6 Performance and Use

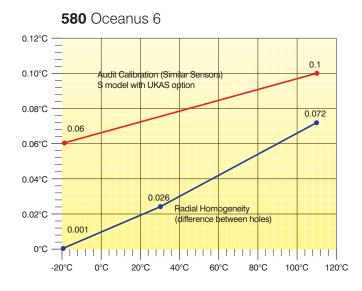

Calibration and Uncertainty

A certificate, traceable to National Standards, is included as standard. Recommended is an optional UKAS five-point calibration.

The accuracy will depend very much on the mode of use and the types of sensor to be used. Please contact Isotech for tutorials and uncertainty calculations and comprehensive evaluation reports. The 580 Oceanus-6^{PLUS} meet the Calibration Capacity requirements of EURAMET/cg-13/v.01, "Guidelines on the Calibration of Temperature Block Calibrators".

Features (Basic & Site)

Dry Block	\checkmark
Stirred Liquid bath Option	✓
Stirred Ice Bath Operation	✓
Surface Sensor Option	✓
Infrared Calibration Option	✓
ITS-90 Fixed Point Cells	Water, Gallium
Additional 8mm Pre-heat Pocket	✓
Configurable Units: °C, °F and K	\checkmark
Supply Voltage Power Correction	✓

Oceanus 6+

Additional Features (Site)

Oceanus 6+

✓
✓
✓
✓
✓
✓
\checkmark
\checkmark



Isocal-6 Ultimate Flexibility - Oceanus

Calibrate all sensor types - Thermocouples, PRT's, Thermistors, Thermostats, Infrared, Surface Sensors...



1 Metal Block Insert 580-06-03
Standard Insert included

Alternative Inserts


580-06-04 Blank Insert without pockets for local machining. Includes M4 tapped hole for supplied extractor tool.

580-06-05 Custom insert. Isotech can provide custom drilled pockets, minimum of 3mm separation between holes. Contact

580-06-06



with your requirements. Adjustable Equalising Block. Six pockets 8mm x 157mm deep.

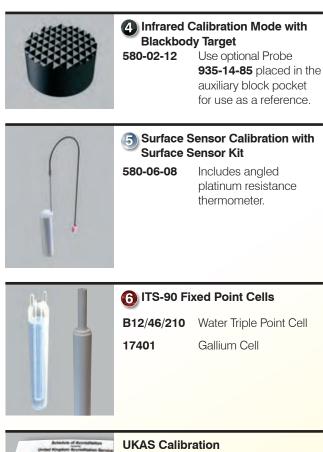




2 Stirred Liquid Mode with Liquid Container Kit 580-06-07 Allows liquid bath use, includes

container, magnetic stirrer, probe guide and sealing cap.

3 Stirred Ice Bath Mode with Liquid **Container Kit**

Uses same liquid kit to provide 0°C reference as a stirred ice bath

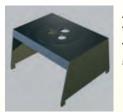

Additional Accessories for use with the above kit

Thermometer Support Kit 580-06-00

Allows three thermometers to be suspended in the bath, including liquid in glass types

C10 Oil 936-06-07

1 Litre (-35°C to +140°C)



UKAS Calibration available to order, legally traceable in more than 70 countries.



Standard Probe

935-14-85/DB Platinum Resistance Thermometer for use up to 250°C.

4.5mm diameter Angled head feature avoids sensors in block.



Carrying Case 931-22-58

Sturdy case accommodates the unit with room for accessories

How To Order Specify Model, Basic or Site, Supply Voltage, Accessories and if UKAS Calibration is required.





45 x 285mm Calibration Volume
 Use for Comparison and Fixed Point Calibration
 Use with very long thermometers

Isotech have a wide range of Dry Blocks to suit probes requiring a large immersion depth. These products feature large and deep calibration volumes. As such they are less portable than the earlier Dry Blocks, but have higher capacities and retain outstanding temperature uniformity, this uniformity is so good that these larger products are also apparatus for Secondary Laboratories to realize the Fixed Points of ITS-90.

Medusa 510 has a maximum operating temperature of 550°C. The Medusa Model 511 can be used to 700°C and features three zone control. In addition to the main heating zone there are additional top and bottom heaters which compensate for the end losses creating a constant temperature zone across the well.

For Comparison Calibration the Medusa should be used with an insert, the standard insert has six 8mm pockets 250mm deep. Also available is an insert 44mm diameter x 170mm deep which is suspended from the top of the block so that the height is user adjustable. For flexibility the Medusa can also be used with accessories for infrared thermometers and surface sensors. The Medusa is available in two models, the BASIC (B) and the SITE (S). The B model includes a sophisticated temperature controller with a dual display for Set Temperature and Dry Block Temperature.

The S model includes a built-in digital thermometer to which an external standard thermometer can be connected giving greater accuracy, eliminating temperature gradient and loading errors. Also included in the site model is a timer which can set the bath between two temperatures, and automate ITS-90 fixed point operation. For Surface Sensor and Blackbody use an external thermometer is recommended. For laboratory accuracy the Medusa can be used with a high-end temperature indicator such as an Isotech TTI model.

Includes as standard: Windows Software, Computer Interface and a Ramp to Set Point Feature. Increased resolution of ± 0.01 available throughout the range via the PC interface and from 0.01 to +99.99 locally on the autoranging front display. The controller features multi-point block to display correction giving good absolute accuracy.

The S model has universal sensor input allowing Platinum Resistance Thermometers, Thermocouples (types K, N, R, S, L, B, PL2, T, J and E) along with Linear Process Inputs including 4-20mA current transmitters to be displayed on the in-built indicator. The indicator can be programmed with up to five calibration points to provide high accuracy digital probe matching. The indicator and controller are both addressable over the communications link.




Fixed Point Cells Available

Material	Temperature
Gallium	29.7646°C
Indium	156.5985°C
Tin	231.928°C
Zinc	419.527°C
Aluminium	660.323°C

LOT	ECH

Specification			
Model	510 Medusa	511 Medusa 3	
Temperature Range	30°C to 550°C	50°C to 700°C	
Absolute stability over 30 minutes	Metal Block Bath Blackbody Source Surface Sensor Calibrator ITS-90 Fixed Point	±0.03°C ±0.1°C ±0.5°C ±0.001°C	
Computer Interface	Included with	Software	
Cools from	550°C to 30°C i	in 5 hours	
Heats from	30°C to 550°C in 90 minutes		
Uncertainties	Refer to Uncertainties Graph		
Calibration volume	45mm diameter by 285mm deep		
Standard Insert	Six 8mm pockets all 250mm deep		
Display Resolution	(0.01) to 9 (0.1) 100.0 to PC can display 0.01 across whole ra	o 650.0	
Indicator units	°C, °F,	K	
Power	108 to 130V or 208 to 1000 Watts	240V 50 / 60Hz 1800 Watts	
Overall dimensions	Height 430 Width 310 Depth 300	Omm	
Weight	17kg	25kg	

Performance and Use

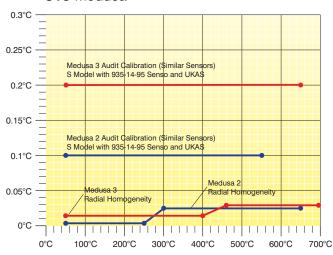

Calibration and Uncertainty

A certificate, traceable to National Standards, is included as standard. Recommended is an optional UKAS five-point calibration.

The accuracy of the Medusa will depend very much on the mode of use, see the Uncertainty Graph for typical uncertainties. NTPL calculate the uncertainties to UKAS requirements. The Medusa meets the Calibration Capacity requirements of EA-10/13, "EA Guidelines on the Calibration of Temperature Block Calibrators."

Features (Basic & Site)	Medusa
Dry Block	✓
Surface Sensor Option	\checkmark
Infrared Calibration Option	✓
ITS-90 Fixed Point Cells	✓
Additional 8mm Pre-heat Pocket	✓
Configurable Units: °C, °F and K	\checkmark
Supply Voltage Power Correction	\checkmark

Additional Features (Site)	Medusa
Independant Temperature indicator	\checkmark
Universal Input Types PT100	✓
Thermocouples Types K,N,R,S,L,PL2,T,J,E	✓
Linear Process Inputs Including 4-20 mA	✓
Stand Alone Thermostat Testing	✓
Thermostat Testing With PC	✓
Five Point Digital Probe Matching	✓
Configurable Units: °C, °F and K	\checkmark

510 Medusa



510 Medusa & 511 Medusa 3

Calibrate all sensor types - Thermocouples, PRT's, Thermistors, Thermostats, Infrared, Surface Sensors...



510 Metal Block Insert			
510-06-01	Standard Insert included		
510-06-02	Blank Insert without pockets for local machining		
510-06-03	Special Insert. Contact Isotech with your requirements		
510-06-04	Adjustable Equalising Block		
511 Metal B	lock Insert		
511-06-01	Standard Insert Included		
511-06-02	Blank Insert without pockets for local machining		
511-06-03	Special Insert. Contact Isotech with your requirements		
511-06-04	Adjustable Equalising Block		



ITS-90 Fixed Point Cells

 ITL17401M
 Galium Slim Cell (510)

 ITL17668ML
 Indium Slim Cell

 ITL17669ML
 Tin Slim Cell

 ITL17671ML
 Zinc Slim Cell

 ITL17672ML
 Slim Aluminium Cell (511)

 510-05-00
 Cell Basket for 510

 510-05-01
 Cell Basket for 511



UKAS Calibration

UKAS Calibration available to order, legally traceable in more than 70 countries.



Blackbody Kit

511-06-05

510-06-05 For 510. Includes a Blackbody target and Sensor.

> For 511. Includes a Blackbody target and Sensor.



510 Surface Sensor Calibration with Surface Sensor Kit

510-06-06 Includes an insert and an angled thermocouple.

511 Surface Sensor Calibration with Surface Sensor Kit

511-06-06 Includes an insert and an angled thermocouple.



Standard Probe

935-14-95H/DB Platinum Resistance Thermometer for use up to 650°C.




Carrying Case

931-22-58

Sturdy case accommodates the unit with room for accessories

How To Order Specify Model, Basic or Site, Supply Voltage, Accessories and if UKAS Calibration is required.



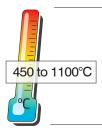

POTTS Oberon

- 50 x 300mm Calibration Volume
- Compact Heatpipe Furnace
- Suits Aluminium, Silver or Copper Fixed Points
- Can be used for comparison and infrared calibration

Model 426 is for Aluminium, Silver or Copper slim fixed point cells as well as for comparison calibration. Heatpipes provide the ideal conditions for the creation and maintenance of slim ITS-90 cells.

The furnace core is a specially-designed stress-free isothermal heat pipe, which provides a very low thermal gradient along the core working length.

The heatpipe is designed so that the inner wall is not subject to thermal expansion stresses from the outer wall before the heat pipe reaches conduction temperature. The working fluid is permanently and safely sealed within the plasma-arc welded enclosure.

The Oberon can be used with Blackbody Fixed Point Cells.



Oberon for Aluminium or Copper Slim Fixed Point Cells

Model	426	Accessories		
Temperature Range	450°C to 1100°C	Metal Block Bath	426-06-05 Adjustable Equalising Block. Six pockets 8mm x 160mm	
Stability	±0.05°C		deep.	
Display resolution	0.1°C	ITS-90 Fixed Point	ITL M 17672QS Aluminium Quartz Clad Slim Cell	
Cavity size	50mm diameter			
Time to temperature	300mm deep 4 hours	Apparatus	ITL M 17673QS Silver Quartz Clad Slim Cell	
			ITL M 17674QS Copper Quartz	
	Supplied as standard with serial interface, PC adaptor cable and		Clad Slim Cell	
	Cal NotePad	Inconel Basket	426-04-00	
Power	110 Vac, 1.5kW, 50/60Hz	including insulators		
	(230 Vac Transformer available)	230/110V	935-19-43	
Dimensions	Height 410mm Width 415mm	Transformer		
	Depth 280mm	How to order		
	Weight 30.5kgs	426 Oberon High Temperature Furnace Please specify voltage required		

ISOTECH

Precision Thermometers

Isotech have a range of innovative precision thermometers to match the calibration requirements of all labs, from the most demanding of National Metrology Institutes through to the needs of those calibrating industrial sensors.

Precision Thermometers

Isotech have a range of precision thermometers, from a two channel handheld thermometer, the TTI-10 to a bench mounting thermometer with performance to 1mK (0.001°C) at a ground breaking new price.

True Surface Temperature Measuring System

This is a true temperature indicator for use with surface temperature measurement, ideal for use with the Small Hot-Plate Model 983.

Semi Standard Resistance Thermometers and Thermocouples

This section includes a range of "Semi Standard" thermometers that can be used with the TTI's. These precision semi standards are more rugged and affordable than the standard thermometers and ideal for industrial applications.

Fixed Resistors

Isotech have a miniature resistor with an ultra low temperature coefficient, model 836, with accuracies of $\pm 0.005\%$ and temperature coefficients of less than 1ppm. These resistors offer an outstanding cost to performance ratio and will find use alongside a model from our TTI range and in other areas of industrial calibration.





In addition to Precision Thermometers a range of Thermometry Bridges are available. The microK range have accuracies of $< \pm 0.1$ ppm to suit Primary and Secondary Laboratories.

Details are in Catalogue 1: Solutions for Primary and Secondary Laboratories.

Thermometer Selection Guide						
Model	SPRTs	PRTs	Thermistors	Thermocouples	Accuracy at 0°C	Features
TTI-10					0.01°C	Handheld, two channel
TTI-22					0.001°C	Sets new Standard for Price to Performance Ratio
milliK					0.003°C	SPRTs, PRTS, Thermistors and Thermocouples
Model 954						8 Channel PRT Switch for TTI-22 and TTI-7 PLUS
Model 958						8 Channel Thermocouple Switch for TTI-22 and TTI-7 PLUS





True Surface

Temperature Measuring System

- Indicates True Surface Temperature
- 30°C to 350°C
- Resolution 0.1°C



Semi Standard

Resistance Thermometers and Thermocouples


- Ideal for Industrial Applications
- Can be supplied with UKAS
- Choice of temperature ranges and sizes



Fixed Resistors

- Wide range of Values
- Oil Filled
- Precise and Stable for Industrial References



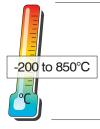

Handheld Thermometer **TTI-10**

- High Accuracy Handheld Thermometer
- High Resolution, to 0.001°C
- Perfect Standard for use with Isocal-6, Fast-Cal & Dry Blocks

The TTI-10 is a high accuracy handheld temperature indicator with two platinum resistance thermometer inputs. The high precision makes the instrument particularly suitable as a portable reference thermometer to use alongside Isotech temperature calibrators such as the Fast-Cal, Isocal-6 and Dry Block ranges. It is also suited for high accuracy measurements in industrial and scientific applications.

TTI-10 brings laboratory level performance of up to 10mK (0.01°C) and resolution up to 0.001°C in a portable handheld instrument. Battery life is typically 20 hours from a 9V PP3 battery and a protective rubber boot offers protection in field use.

The instrument can capture the minimum, maximum and average values over up to 4000 measurements with a logging rate selectable in the range of 1 second to 30 minutes.

The TTI-10 has an easy to use "learning calibration mode" that allows the TTI-10 to be system calibrated with a Platinum Resistance Thermometer simply by comparing it to a calibrated standard thermometer, no need to calculate coefficients or data, simply enter the reference probe temperature or temperatures and the TTI-10 does the work for you.

The USB interface allows connection to Isotech Cal Notepad software with its charting and logging features.

TTI-10 supports Isotech Semi Standard Platinum Resistance probes with system uncertainties (probe and instrument) as low as 20mK. We recommend the 935-14-61 and 935-14-16 probes detailed below and have special calibration deals available. Other probes and ranges are available, refer to Semi Standards – Platinum Resistance Thermometers in catalogue.





Input Connectors Highest quality latching metal 'Lemo' connectors.



Rubber Sleeve The TTI-10 Handheld Thermometer is supplied with a protective rubber boot.



Specifications

Input Channels	Two: 100 Ohm PRT, EN 60751	Connectors	High Quality Latching Metal: Lemo:
Range	(Pt100), Four Wire -200°C to +850°C	Working temperature	0°C to +40°C
Units	°C, °F and Ohms	Display	2-line LCD
Resolution	0.001°C from -199.999°C to +199.999°C remaining range 0.01°C		Display Single Channel or Dual Channels Simultaneously
Accuracy: Instrument	±0.012°C from -80°C to 199.999°C ±0.02°C ±0.0015% RDG from 200°C	Housing	Plastic (ABS) supplied with protective rubber boot
Only	to 660°C	Weight	300g
Logging	Record Average, Min and Max over 4000 measurements	Power Supply	9V battery PP3 (or via USB Cable)
		Battery Life	Typically 20 Hours
Measuring interval	Adjustable: 1 second to 30 minutes	Dimensions	200 x 85 x 40 mm (LxWxH)
PC Interface	USB - Cable Included		

Options

Semi Standard PRT 935-14-112-TTI	Isotech Semi Standard Platinum Resistance Thermometer: Fast Response, 2m Cable Length, four wire with Lemo plug fitted	
UKAS System Calibration TTI-10-14-112-SYST	Recommended: -50°C to 199.999°C Four Point System Calibration, Uncertainty across range 0.025°C (25mK)	
Semi Standard PRT 935-14-61-TTI	Isotech Semi Standard Platinum Resistance Thermometer: Fast Response, 2m Cable Length, four wire with Lemo plug fitted	
UKAS System Calibration TTI-10-14-61-SYST	Recommended: -50°C to 199.999°C Four Point System Calibration, Uncertainty across range 0.02°C (20mK)	

Semi Standard PRT 935-14-116-TTI

UKAS System

Calibration

Isotech Semi Standard Platinum Resistance Thermometer: General Purpose,2m Cable Length, four wire with Lemo plug fitted

Recommended: 0°C to 420°C Four Point System Calibration, TTI-10-14-116-SYST Uncertainty across range 0.04°C (40mK)

Carrying Case 931-22-101






Recommended Probes (Fit TTI-10 Carry Case)

Model	Maximum Range	Diameter	Length (A)	Sensing Length (B)	Handle (C)	Cable	Application
935-14-112/TTI	-50°C to 250°C	3mm	225mm	6mm	No Handle	2m PTFE	Fast Response, Low Stem Conduction
935-14-61/TTI	-50°C to 250°C	4mm	300mm	6mm	19 x 120mm	2m PTFE	Fast Response, Low Stem Conduction
935-14-116/TTI	-100°C to 450°C	6mm	350mm	25mm	19 x 120mm	2m PTFE	General Purpose

-B



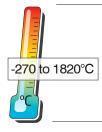

Precision Thermometer millik

- Wide Range of Sensors, SPRTs, PRTs, Thermistors, Thermocouples and 4-20mA Transmitters
- High Accuracy, < ±5ppm for PRTs, ±2µV for Thermocouples and ±1µA for Transmitters
- Logs and Controls Isotech Temperature Sources Massive logging capacity - supports Dry Blocks and Liquid Baths

The millik Precision Thermometer from Isotech sets a new standard for the high accuracy measurement and calibration of Platinum Resistance Thermometers, Thermistors, Thermocouples and Process Instrumentation (4-20mA) over the range -270°C to 1820°C.

In addition to low uncertainty measurements from Reference Standards and Industrial sensors, the milliK can control Isotech temperature sources, sequencing through a programmable list of temperature set points and log data to internal memory or a USB drive.

The milliK has two input channels for sensors and a third channel for current. It can be expanded to become a measuring system with up to 33 channels reading SPRTs, RTDs, Thermistors, or Thermocouples with the option to control calibration baths and log readings accurately.

Benefiting You

The milliK sets a new standard for value, versatility and accuracy - < \pm 5ppm over range for PRTs, $\pm 2\mu$ V for Thermocouples and $\pm 1\mu$ A for current transmitters, see table.

Supporting a wide range of sensors and functions it replaces individual devices making it a cost effective calibration solution.

A robust design and operation from AC or DC power allows the milliK to be used in the laboratory, test room or out in the field.

The milliK can display in °C, °F, K, Ohms, mV and mA with numeric and graphical display modes. The large back lit display makes configuring the instrument and setting the scrolling strip charts intuitive. The USB port allows for the use of a mouse, keyboard or USB Drive.

Built on World Leading Technology

In 2006 Isotech launched the microK range of thermometry bridges which quickly established themselves as the instrument of choice for National Metrology Institutes and Primary Laboratories with innovative features, accuracy and versatility.

In response to industry demands for greater accuracy, the milliK now brings the same design philosophy used in the microK to those outside the Primary Laboratory. Users calibrating industrial sensors in the laboratory, pharmaceutical plants, food and beverage plants, aerospace, power industries and service companies will welcome the milliK as a solution to increase measurement confidence, ensure high accuracy traceable calibration, improve quality as well as ensure safety and lower energy consumption.



No Compromise Design

The design team have considered industrial users and applications in order to avoid measurement errors and problems encountered in some instruments from other manufacturers:

- Eliminates Thermal EMF Errors in PRTS Fast current reversal technology and solid state switching eliminate thermal EMF effects avoiding the errors that occur with fixed DC instruments.
- Lead Wire Correction

PRT lead wire errors are eliminated for up to 30m of four core screened cable. Also supports lead wire correction for three wire PRTs.

Galvanic Isolation

Not only are the two sensor channels galvanically isolated, the 4 - 20mA input is also separately isolated. The benefits of the advanced design are no ground loops, improved safety and noise immunity.

High Resolution

The display resolution is 0.0001°C (0.1mK) made possible by using a powerful Sigma Delta Analogue to Digital converter to achieve a true measuring resolution of just $28\mu\Omega$ equivalent to 0.00007°C (0.07mK) for PRT inputs.

Expandable

The millisKanner adds eight channels, and each can be configured individually as a SPRT, PRT, Thermistor or Thermocouple input. A maximum of four millisKanners can be added, providing up to 32 channels - all controlled from the milliK touch screen or an RS232 connection.



Reliable

Like the award winning microK range, the milliK is all solid state. There are no mechanical relays, switches or potentiometers which would reduce reliability.

Input Connectors

No compromise design ruled out lower cost problematic connectors and the SPRT / PRT inputs are via the highest quality gold plated push / pull self latching circular connectors overcoming the problems seen elsewhere where thermometers have been designed to a budget.

Outstanding CJC Performance and Flexibility

Again, the no compromise design philosophy led to a specially developed rugged thermocouple connector made from alumina and incorporating a digital temperature sensor for optimal cold junction accuracy.

Three CJC modes allow thermocouple operation with internal automatic compensations, external 0°C reference systems or the milliK can measure the



1 The milliK can connect to Isotech temperature sources

Dry Blocks (Basic & Site only), Liquid Baths and Furnaces Can cycle the bath through a series of temperatures logging the data - all without a PC.



4 Logs

The milliK can record time stamped data to internal memory or a USB Memory Drive.

junction with a probe on an unused channel, useful for automated systems.

21st Century Design

Utilising a powerful internal operating system and fast 32 Bit processor the milliK has the power and capacity to overcome the memory limitations of older instruments.

Store Probe Data

There is sufficient memory for an almost unlimited number of standard probes, allowing the storing of calibration data for both resistance thermometers and thermocouples. The digital matching of probe data allows the instrument to show the true temperature. The instrument will warn if a probes calibration time has expired.

Data Logging

Older instruments are limited to a maximum number of logged data points, the milliK is limited only by storage space. The internal memory can store more than six months of data, and with a low cost USB Memory stick the milliK can log continuously for a lifetime

Data Management

Probe data and logged measurements can be exported to a USB Memory drive at the push of a button. Additionally, the instrument is future proof with future software updates applied from a USB drive.

Connectivity and Communications

With USB host, two serial interfaces and Ethernet it is easy to communicate with the milliK whether it is on the bench next to a PC or remotely using a LAN or WAN connection. These interfaces are fitted as standard.

The milliK includes a PC lead and software.


Open Calibration

The millik is readily calibrated against resistance and voltage standards. There are no internal adjustments and the calibration commands are simply sent via RS232 or from the front panel (password protected). The procedure is open and fully documented unlike some other instruments where there is no choice but to return to the manufacturer.



2 Wide range of sensors

The milliK can use Standard Reference probes and read from industrial sensors being calibrated, including 4 - 20mA transmitters - all to high accuracy.

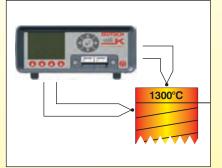


5 Safety

The milliK inputs are galvanically isolated, with the 4 - 20mA input separately isolated avoiding problems with high voltage pick up common when using thermocouples in high temperature furnaces.



3 Expandable


The milliK can be expanded to have a maximum of 33 high accuracy channels. The millisKanner has eight expansion channels, with each channel configurable for SPRT, PRT, Thermistor or Thermocouple input type.

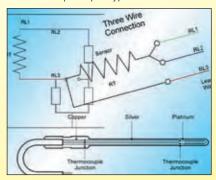

6 Designed to eliminate and protect against real world problems

The milliK eliminates thermal EMF errors, compensates for lead wire resistance and warns if a probe is out of calibration.



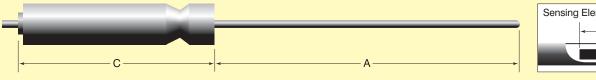

Recommended Probes (Fit millik Case)

Model	Maximum Range	Diameter	Length (A)	Sensing Length (B)	Handle (C)	Cable	Application
935-14-61/TTI	-50°C to 250°C	4mm	300mm	6mm	19 x 120mm	2m PTFE	Fast Response, Low Stem Conduction
935-14-116/TTI	-100°C to 450°C	6mm	350mm	25mm	19 x 120mm	2m PTFE	General Purpose

For further options and details, see Reference Probes - Semi Standards, pages 76-81.

For laboratory standard thermometers we recommend for SPRTs the Isotech Model 909Q and for thermocouples the Model 1600 Type R, see Catalogue 1: Solutions from Primary & Secondary Laboratories.



- The choice for high accuracy temperature measurement - expandable for multichannel operation
- Higher Accuracy than DAQ Systems
- Ideal for industrial sensor calibration alongside Dry Block and Liquid Baths.
- Expandable to 33 channels with no loss of accuracy


milliK shown with optional channel

expander

annel



UKAS Calibration available for these systems - International Traceability - Best Practice

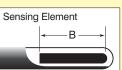




Specifications					
Input Channels Channels 1+2	3 SPRTs, PRTs, Therr Thermocouples		Input Connectors		LemoEPG.1B.306. HLN 6-pin gold plated contacts
Channel 3	Process Inputs 4 - 20mA Isolated 24VDC Power Supply Included			Thermocouples:	Miniature Thermocouple socket (ASTM E
Ranges	SPRTs: 0-115 PRTs: 0-460 Thermistors: 0-500	Ω		4-20mA:	1684-05) 4mm sockets
Display Units	Thermocouples: ± 4-20mA: 0 °C, °F, K, Ω, mV, mA	30mA	Interfaces	10/100MBit Ethernet (RJ45 socket) USB (2.0) host 2 x RS232 (9-pin D-type plug, 9600 Baud)	
Accuracy SPRTs/PRTs: Thermistors:	Initial 5ppn 50pp	n 7ppm	Display	89mm / 3.5" QVGA (TFT LCD with LED b	
Thermocouples: 4-20mA:	2μV 0.01%	4μV	Operating Conditions	Operating:	0-45°C / 32-113°F, 0-99% humidity
Temperature Accu	uracy Initial	Over 1 year		Full Specification:	15-30°C / 50-85 °F, 10-90% humidity
(over full		4mK 7mK	Statistics	In Addition to Instant	
Thermisto Thermocouples:	Ice Point Ref	Internal CJC		user can select mean measurements with	
Type B @ 1000°C Type E @ 600°C Type J @ 600°C Type K @ 600°C Type L @ 600°C Type N @ 600°C	$ \begin{array}{ c c c c c c c c c c c c c c c c c c c$	$\begin{array}{ c c c c c c c c c c c c c c c c c c c$	Measurement Time	PRTs (4-wire): 0.4s (3-wire): 0.7s Thermistors: 0.4s Thermocouples (ice (internal CJC): 0.7s (external CJC): 1.0s	point): 0.4s
Type R @ 1000°C Type S @ 1000°C Type T @ 200°C Au-Pt @ 600°C	±0.09°C ±0.12°C ±0.10°C ±0.14°C ±0.02°C ±0.03°C ±0.06°C ±0.08°C	±0.16°C ±0.24°C	Cable Length	Limited to 10Ω per c capacitance (equival typical 4-core screer	lent to 100m of
Resolution	Resistance (PRTs): (Thermistors): Voltage: Current:	0.00001Ω 0.001Ω 0.00001mV 0.0001mA	Logging	Capacity to store > stamped measureme memory	
	Temperature:	0.0001°	Recommended Probes	Isotech Semi Standa Isotech Model 909 S	
Temperature Conversions	PRTs:	IEC60751(2008), Callendar-van Dusen, ITS90	Power	88-264V (RMS), 47-6 6W maximum or 4 x	
	Thermocouples:	IEC584-1 1995 (B,E,J,K,N,R,S,T), L, Au-Pt	Dimensions	255mm x 255mm x 1 4.5" (W x D x H)	114mm / 10" x 10" x
	Thermistors:	Steinhart-Hart, polynomial	Weight	2.25kg / 5lb	
Sensor Currents	SPRTs/PRTs:	1mA and 1.428mA ±0.4% (reversing)	Optional Carring Case	931-22-102	
	Thermistors:	5μA (reversing)			
Keep-Warm Current	SPRTs/PRTs:	1mA and 1.428mA			

NOTE: Due to our program of continual development and improvement, we reserve the right to amend or alter characteristics and design without prior notice.



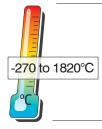

Channel Expander millisKanner

- Expands milliK to a maximum of 33 Channels
- Supports SPRTs, PRTs, Thermistors and Thermocouples
- Universal Inputs for Flexibility

The millisKanner channel expander has eight input channels, and each can be configured individually for SPRT, PRT, Thermistor or Thermocouple input. This gives ultimate flexibility with no need for separate devices for resistance or thermocouple inputs. A maximum of four devices can be added to the milliK providing 33 sensor inputs as well as the 4 - 20mA Process Input.

With no loss of accuracy and total flexibility a milliK system can be configured to suit a wide range of reference thermometers and units under test. This adaptable system saves on cost with no need for separate dedicated expansion modules and the flexibility maximises the usefulness of each channel.

A solid state design avoids mechanical relays and provides high reliability. The inputs are isolated with galvanic isolation between the contacts and the PSU and also from the control circuitry which allows for better measurements and lower noise.



The millisKanner is controlled from the milliK with plug and play operation.

For use as a standalone switch for PRTs, the device has UP / DOWN touch buttons or can be operated via RS232. The temperature of the input thermocouple connectors can be read directly over RS232 to facilitate reference junction compensation.

pply provided)



UKAS Calibration available for these systems - International Traceability - Best Practice

Model	millisKanner	Max devices per milliK	4
Channels Channel Indication	8 LED	Power	5V DC at 100mA (mains power sup
Input Connectors Termin plugs	Terminal Post, accepts 4mm	CJC Sensor	Digital
	plugs, spades and bare wires Miniature Thermocouple Sockets	Size Weight	255 x 255 x 80mn 2.3kg
Control	Automatic: Plug and Play from the milliK Manual: UP / DOWN buttons or RS232		





- Add up to 4 millisKanners
- Expand to a maximum of 33 Channels
- Same Accuracy no loss in performance
- No Need for Multiple Modules each channel can be selected for 3 or 4 Wire PRTs, SPRTs, Thermocouples or Thermistors

millisKanner shown here with milliK, order seperately


SOTECH

Thermometer

http://www.isotechna.com

R



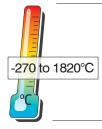

Channel Expander millisKanner

- Expands milliK to a maximum of 33 Channels
- Supports SPRTs, PRTs, Thermistors and Thermocouples
- Universal Inputs for Flexibility

The millisKanner channel expander has eight input channels, and each can be configured individually for SPRT, PRT, Thermistor or Thermocouple input. This gives ultimate flexibility with no need for separate devices for resistance or thermocouple inputs. A maximum of four devices can be added to the milliK providing 33 sensor inputs as well as the 4 - 20mA Process Input.

With no loss of accuracy and total flexibility a milliK system can be configured to suit a wide range of reference thermometers and units under test. This adaptable system saves on cost with no need for separate dedicated expansion modules and the flexibility maximises the usefulness of each channel.

A solid state design avoids mechanical relays and provides high reliability. The inputs are isolated with galvanic isolation between the contacts and the PSU and also from the control circuitry which allows for better measurements and lower noise.



The millisKanner is controlled from the milliK with plug and play operation.

For use as a standalone switch for PRTs, the device has UP / DOWN touch buttons or can be operated via RS232. The temperature of the input thermocouple connectors can be read directly over RS232 to facilitate reference junction compensation.

pply provided)



UKAS Calibration available for these systems - International Traceability - Best Practice

Model	millisKanner	Max devices per milliK	4
Channels Channel Indication	8 LED	Power	5V DC at 100mA (mains power sup
Input Connectors Termin plugs	Terminal Post, accepts 4mm	CJC Sensor	Digital
	plugs, spades and bare wires Miniature Thermocouple Sockets	Size Weight	255 x 255 x 80mn 2.3kg
Control	Automatic: Plug and Play from the milliK Manual: UP / DOWN buttons or RS232		





- Add up to 4 millisKanners
- Expand to a maximum of 33 Channels
- Same Accuracy no loss in performance
- No Need for Multiple Modules each channel can be selected for 3 or 4 Wire PRTs, SPRTs, Thermocouples or Thermistors

millisKanner shown here with milliK, order seperately

SOTECH

Thermometer

http://www.isotechna.com

R





True Temperature Indicator TTI - 22

- Accuracy to 0.001°C, 1mK
- Warns if calibration due date exceeded
- No mechanical relays, long life

Quite simply the Isotech TTI-22 High Accuracy Thermometer sets new standards in the price to performance ratio for industrial and secondary resistance thermometry. If you need high accuracy at an affordable price you have to look at the TTI-22.

The TTI-22 has an accuracy of 0.001°C and a resolution of 0.0001°C (0.00004 Ohms). It has two input channels, is lightweight (1.8kg) and will operate for more than 10 hours from two small AA cells. It has both RS232 and Ethernet ports.

Simple to use, supporting both Industrial 100 Ohm probe and SPRTs to ITS-90, 25.5 and 100 Ohm. Up to 30 probe calibrations can be stored along with the calibration expiry date so the instrument can warn when the calibration time has been exceeded.

Built in statistics calculation can show you both the measured and average values along with the standard deviation over previous measurements.

The Isotech TTI-22 is ideal as a reference standard alongside liquid calibration baths, for the smallest uncertainty calibration with Dry Blocks or for demanding stand alone measurement applications.

Previously this level of performance was confined to specialist laboratories with expensive thermometry bridges; TTI-22 delivers 5 to 10 times the performance of comparably priced instruments.

- The TTI-22 uses the same patented measurement technique as the earlier TTI-2.
- Each measurement performs a zero point and gain correction.
- The switched polarity DC measuring current (0.4mA) eliminates thermal EMFs.
- Surface mount construction ensures long term reliability.



Model	TTI-22
Inputs	2 channel Pt100 (BS EN 60751 / IEC 751) or 25.5/100Ω SPRT to ITS-90
Measuring Current	0.41mA
Self Heating Test Current	0.29mA (0.41mA / √2)
Measuring Time	1.44 seconds for both channels
Measuring Range	-250 to 960°C (0 to 440 Ohm)
Resolution	Temperature: 0.0001°C, 0.1mK Resistance: 0.00004Ω, 40 μ Ω
Uncertainty of Measurement	Temperature: 0.001°C, 1mK 100 Ohm PRT Resistance: 0.4m Ω @ 20°C Instrument only, uncertainty with sensor dependant on range and sensor type.
Reference Resistor	Internal 380Ω TCR ±0.3ppm / °C Stability ±5ppm / year
Interface	RS232, Ethernet, built-in web server provides simple temperature display
Ambient Temp. Range	10°C to 30°C
Power Supply	7.5VDC, 250mA power adaptor or 2 x AA batteries (typically >10 hours operating time)
Case Dimensions	Width: 190mm Height: 112mm Depth: 240mm Weight: 1.8kg



The TTI-22 continually compares the connected sensor to a highly stable precision internal reference resistor. For a Pt100 at 0°C the annual stability for absolute measurement is typically ± 1.3 mK (5ppm x 100 Ω = 0.5m Ω / 1.3mK).

For comparison calibration, when a reference probe is compared to a calibrated standard, the long term stability is not important as any change of value is cancelled in the comparison. The temperature coefficient is 0.3ppm / °C and the measuring time, for both channels, is just 1.44 seconds.

The instrument can be configured to measure ratio of the measured

resistance of the two input channels, a technique familiar to users of older style thermometry bridges.

The overall uncertainty of the instrument and probe together will be determined by the model of probe and the temperature range. For the majority of applications the contribution of the instrument uncertainty will be negligible compared to the uncertainty of the calibrated probe.

Recommended probes include the Isotech 909/100 and 670SQ /100, 935-14-16, 935-14-95L and H. The TTI-22 includes Cal Notepad software for easy monitoring and logging of data. It is fully compatible with Isotech I-Cal Easy which can automate comparison calibration.



931-22-106 Optional Carrying Case





Terminal Adaptor Model 956


- Time Saving
- Accepts Bare Wire, Spades or Banana Plugs
- Suits a wide range of Isotech instruments

Our TTI range use high quality 'Lemo' connectors for the Pt100 inputs. This simple accessory provides 4mm Terminal Posts for the connection of bare wires, spade terminals or 4mm plugs - useful if a lot of probes are going to be used with the instrument.

The adaptor connects to the TTI via two flexible cables, terminated with the appropriate Lemo connector.

New compact design.



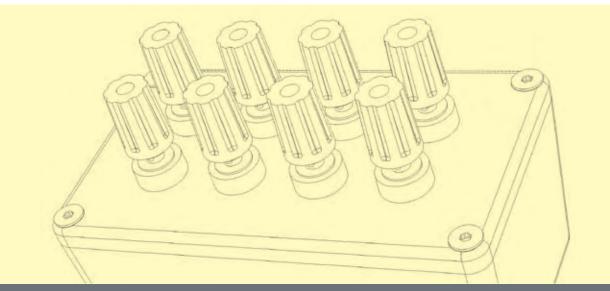

Specifications

Dimensions	Height 65mm Width 115mm Depth 105mm (including plugs)
Weight	0.4kg

How to Order

956 Terminal Adaptor suits milliK, TTI-6, TTI-7 or TTI-10 956/TTI-22 Terminal Adaptor for TTI-22



Miniature Fixed Resistor Model 836

- Wide range of Values
- Oil Filled
- Precise and Stable for Industrial References

Isotech produces a miniature resistor with ultra-low temperature coefficient and ultra-high stability.

This is achieved because the resistors are oil filled and hermetically sealed.

The function of hermetic sealing is to eliminate the ingress of moisture and oxygen both of which play a role in both short and long term degradation of unsealed resistors. A further enhancement in both short and long term stability is achieved by oil filling. The oil also acts as a thermal conductor allowing the device to accept short periods of overload without degradation.

With accuracies of $\pm 0.005\%$ and long term drift of less than 5ppm, these devices are virtually secondary standards that can be carried in sets for daily or periodic calibration of factory systems.

Resistance Values

We keep in stock the following standard values: $10\Omega,\,25\Omega,\,100\Omega,\,1000\Omega,\,10,000\Omega$

UKAS Calibration

For the highest quality traceability we recommend that the 836 be UKAS Certified.

Measured Quantity Instrument or Gauge Range	Frequency	Best measurement Capability expressed as an Expanded Uncertainty (<i>k</i> =2)
DC Resistance 0.1 Ω to 1000 Ω 1 KΩ to 100 MΩ		±10ppm ±20ppm
AC Resistance 2.5 Ω to 400 Ω 400 Ω to 1000 Ω	75 Hz 75 Hz	±15ppm ±100ppm
		the second se

The latest schedule can be found on the Isotech website or at www.ukas.org.



Please Note:

We offer other Resistor Ranges including the models 456, the SRA and the SRB ranges.

For more information please contact lsotech, or visit our website www.isotech.co.uk



Model	836 Miniature Fixed Resistor
Power Rating	0.5 watt
Nominal Temperature Coefficient of Resistance	+0.6ppm/°C (0°C to +25°C) -0.6ppm/°C (+25°C to +60°C)
Resistance Tolerance	(Initial Resistance Accuracy) ±0.005%
Resistance Range	5 ohms to 3.3 megaohms
Current Noise	${<}0.010 \mu\text{V}$ (RMS) / Volt of applied voltage
Thermal EMF	0.1µV/°C maximum 0.05µV/°C typical
Connections	Screw Terminal Posts
Stability	Typically 1ppm per year at 1mA
Dimensions	Height 30mm Width 89mm Depth 58mm (including terminals)
Weight	90g

How to Order

836 Miniature Fixed Resistor Please state Ohms Value Required Please state if UKAS Certification is required



Selector Switch 8 Way

- Eight Channel Local and RS232 Control
- PRT and Thermocouple Models
- Use with Isotech TTIs and Automation Software

Isotech produces two eight way selector switches, one for resistance thermometers Model 954 and Model 958 for thermocouples.

These switches have been designed for use in conjunction with our TTI range. The switches allow easy selection of connected sensors. They can be operated from either the front panel switch or from an RS232 interface that is provided as standard. Channel status is indicated via front panel LEDs. The Selector Switches can be located adjacent to the sensors being calibrated, giving more flexibility than a permanently connected or stacked system.

The PRT Switch has 4mm terminal posts that can accept bare wires or 4mm plugs. The thermocouple switch has eight miniature thermocouple connectors. These thermocouple connectors are thermally bonded to a platinum resistance thermometer that measures the temperature of the connector and hence the "cold junction".

The TTI range temperature indicators feature the ability to measure a remote cold junction and this permits a mixture of thermocouple types to be connected through the box. The I-Cal Easy Software supports Switchbox models 954 and 958 and, for automatic operation, two boxes can be connected together with a "master / slave" lead allowing them to be controlled from a single RS232 port and up to 16 sensors to be switched.

The software can automatically switch between the boxes and connect the appropriate output to the TTI. This 16 channel operation is not convenient without the software and manual operation of two boxes together is not recommended.


Advantages

- Use with TTI-6 and TTI-7 PLUS easily switch up to eight sensors manually or with RS232.
- RTD and Thermocouple Models.
- Use with I-Cal Easy Software for automatic switching and temperature calibration, add a second box (either type) to calibrate up to 16 sensors.
- Switches are stand-alone allowing them to be positioned anywhere in a laboratory for most efficient operation.








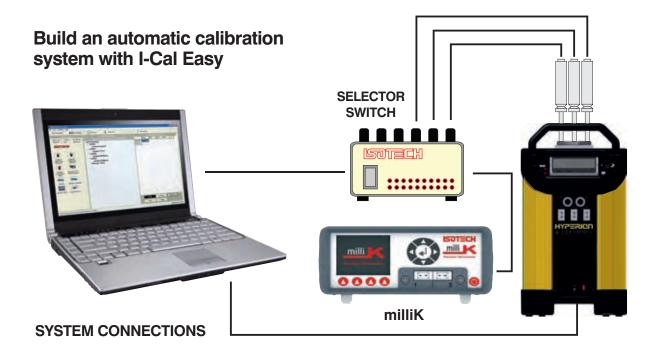

Model	954 RTD Selector Switch	Model	958 TC Selector Switch
Channels	Eight - four wire (four pole)	Channels	Eight - two wire (two pole)
Control	Front panel switch And RS232	Control	Front panel switch And RS232 (Also compatible with Isotech VLT
Connectors	4mm Terminal post		system)
Internal Circuit Resistance	<250mΩ	Connectors	Miniature Thermocouple Connectors
Thermal EMF, typical	2μ V after 1 minute of channel set 6 μ V after 30 minutes of channel set	Internal Circuit Resistance	<250mΩ
Power	5 VDC 100-250 VAC, 50 / 60Hz	Thermal EMF, typical	2μ V after 1 minute of channel set 6μ V after 30 minutes of channel set
	Power Supply Included	Reference Junction Measuring Device	100Ω 1/10 Din Pt100
Dimensions How to Order	Height 91mm Width 141mm Depth 165mm Weight 1kg	Thermal Coupling	<0.2°C* *Basis of test. At ambient 20°C ±2°C the internal Pt100 agreed with the connected thermocouples to ±0.2°C (including all measurement errors)
954 RTD Select	or Switch		using IEC584-1995 and IEC751-1995. The uncertainty of this test was $\pm 0.3^{\circ}$ C which includes the reproducibility of the test thermocouples.
		Power	5 VDC 100-250 VAC, 50 / 60Hz Power Supply Included
		Dimensions	Height 64mm Width 141mm Depth 165mm Weight 1kg
		How to Order	

How to Order 958 TC Selector Switch



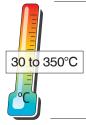

Surface Measurement Model 944

- Indicates True Surface Temperature
- 30°C to 350°C
- Resolution 0.1°C

The fundamental problem with surface temperature measurement is that it is subject to large stem conduction errors, also because heat conducted from the surface of the hot-plate causes a localised cold spot to be created which means that the temperature indicated by the hot plate is not necessarily the temperature at the point of measurement.

An ideal system would not disturb the heat-flux from the hot-plate.

During 1993 such a system was described (ref. "Progress in Contact Thermometry" 1993 B. D. Foulis) and Isotech have the inventors permission to make and market the device World-wide.

Principal of Operation

A fine wire type N thermocouple is used as the surface temperature sensor, a second junction 2 to 3mm along the thermocouple, senses the temperature difference due to heat flux along the sensor.

A heater heats the thermocouple stem until the temperature gradient is zero, thus creating a measurement without stem conduction, or disturbance of the hot-plate's surface.

The 944 can be used with the Isotech Small Hot Plate model 983. A traceable calibration certificate is available to order.



Model	944 True Surface Temperature Measurement System				
Temperature Range	30°C to 350°C				
Resolution of display	0.1°C or 0.1°F				
Stability	±1°C				
Accuracy	$\pm 2^{\circ}$ C with TRACEABLE Certification $\pm 5^{\circ}$ C without Certification				
Probe Assembly	Probe Diameter Probe Length Lead Length	150mm			
Power Supply	100V - 120V, 50 / 6 200V - 240V, 50 / 6				
Dimensions	Height 90mm Width 153mm Depth 265mm (exc	cluding plugs)			
Weight	4kg				
How to Order					

Model 944 & Probe 935-14-81 Please state supply voltage required Please state if Calibration is required





Surface Sensor Calibrator Small Hot Plate

- Low Cost Portable Hot Plate
- PC Interface and Software
- Stable to ± 0.1°C

The Isotech Small Hotplate is a lightweight portable calibration system purpose designed for surface mounted sensors. The flat surface plate is made from precisionmachined aluminum. The sensor to be tested is simply placed on the surface, for higher accuracy a calibrated surface sensor can be placed alongside and the two compared.

Good thermal contact is ensured by the flat disc that is recessed to allow the optional use of a heat transfer paste or fluid. Uniform heat distribution is achieved with a flat spiral heater clamped to an integrating block below the surface of the plate. The typical accuracy that can be achieved 1°C but this will be influenced by the type of sensor to be calibrated.

The internal control sensor is located immediately below the plate's surface.

A protective cover that can fit over the block is included along with a comprehensive handbook.

The temperature range is from 35°C to 350°C, which is set by an advanced, but easy to use temperature controller. The controller has 0.01 resolution below 100°C (0.1° above 100°). A PC interface is included as standard along with an RS232 converter lead and Windows software.




Notes:

A similar model but with a black high emissivity surface is available. Many of the dry block calibrators featured within this book have accessories available for surface sensor calibration.

Model	983 Small Hot Plate
Temperature Range	35°C to 350°C
Stabilisation Time	10 minutes
Cools from	350°C to 100°C in 125 minutes
Heats from	50°C to 350°C in 20 minutes
Uncertainties	Dependant on sensors and method of use 1°C typical
Calibration volume	Flat Plate 71mm diameter
Display Resolution	0.01 to 99.99 0.1 100 to 350.0 PC can display 0.01 across whole range with the software included
Units	°C, °F, K
Power	100 to 115V (50 / 60 Hz) or 200 to 230V (50 / 60 Hz) 200 Watts
Dimensions	Height 115mm Width 230mm Depth 225mm
Weight	3.9kg
How to Order 983 Small Hot Plate	

Please specify voltage required

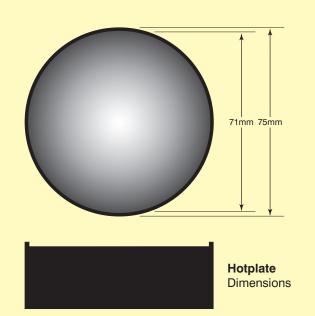



Reference Probes – Semi Standards Platinum Resistance Thermometers

- High Stability Reference Probes
- Wide Temperature Ranges
- High Stability Platinum Coil Elements

These industrial platinum resistance thermometers are ideal for field and lab use. Suitable for use as working standards in Dry Blocks and Liquid Baths or as high accuracy probes for our range of True Temperature Indicators.

All the thermometers are metal sheathed and both less fragile and more affordable than the Isotech range of true Standard Platinum Resistance Thermometers that are used in laboratories and are found in our publication *"Solutions for Primary and Secondary Laboratories"*.

All the thermometers use handmade coil wound platinum sensing elements to give high accuracy and low drift. Isotech's UKAS accredited lab can calibrate to the smallest of uncertainties.

Calibration should be specified to suit the particular operating range and application. Isotech can advise on which service is appropriate to match the temperature range and application.






Universal Specifications

Ro Alpha Standard Stability Recommended Current Self Heating at 1mA Calibration

Connection Max. Handle Temperature $\begin{array}{l} 100\Omega\pm0.05\ \Omega\\ 0.003850\pm0.000005\\ \text{IEC }60751\\ 0.010\ \Omega/\text{year}\\ 1\text{mA}\\ 0.004^\circ\text{C}\\ \text{Optional UKAS Calibration at extra cost.}\\ \text{See table for typical uncertainties}\\ \text{Four Wire}\\ 80^\circ\text{C}\\ \end{array}$


After manufacture all Isotech Semi Standard PRTs are thermally pre-conditioned to provide optimal stability.



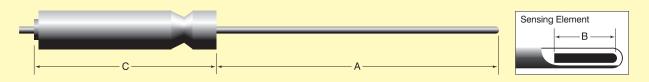


General Purpose Probes

Model	Range	Diameter	Length (A)	Sensing Length (B)	Handle (C)	Cable	Application
935-14-112	-50°C to 250°C	3mm	225mm	6mm	No Handle	2m PTFE	General Purpose/TTI-10
935-14-61	-50°C to 250°C	4mm	300mm	6mm	19 x 120mm	2m PTFE	Fast Response, Low Stem Conduction
935-14-13	-196°C to 250°C	6mm	350mm	25mm	25 x 115mm	2m PTFE	Low Temperature
935-14-113	-100°C to 250°C	6mm	350mm	25mm	19 x 120mm	2m PTFE	General Purpose
935-14-16	-100°C to 450°C	6mm	450mm	25mm	19 x 120mm	2m PTFE	General Purpose
935-14-116	-100°C to 450°C	6mm	350mm	25mm	19 x 120mm	2m PTFE	General Purpose/Fits milliK Case
935-14-72	-50°C to 670°C	6mm	375mm	25mm	No Handle	2m PTFE	Fits Jupiter / Gemini Carry Case
935-14-98	-50°C to 350°C	4mm	300mm	8mm	No Handle	2m PTFE	Low Stem Conduction

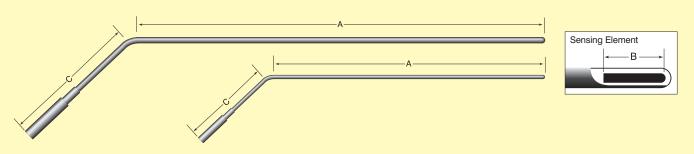

Angled Probes - angled head provides maximum clearence at top of calibration bath

Model	Range	Diameter	Length (A)	Sensing Length (B)	(C)	Cable	Application
935-14-82	-50°C to 250°C	4mm	165mm	6mm	50mm	1.5m PTFE	Europa - Venus - Calisto
935-14-85	-50°C to 250°C	6mm	385mm	25mm	35mm	0.54 m PTFE	Oceanus-6

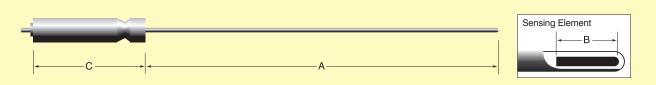

Working Industrial Standards

These thermometers use premium grade wire wound elements to IEC-751 and the same internal construction as our working Standard SPRTs. The 95L is optimised for low temperature with minimum stem conduction. The 95H is optimised for high temperature operation. Both models employ strain free construction.

Model	Range	Diameter	Length (A)	Sensing Length (B)	Handle (C)	Cable
935-14-95L	-200°C to 165°C	6mm	480mm	25mm	25 x 115mm	2m PTFE
935-14-95H	-80°C to 670°C	6mm	480mm	25mm	19 x 120mm	2m PTFE

Termination Options

Bare Wire (BW) TTI suits milliK and TTI-1 to TTI-7, TTI-b – suits TTI-22 DB Connector for Dry Block Calibrator Site Indicator

How to Order

Please Specify Model Type and Termination Option (for example 935-14-13/BW) Please state whether UKAS Certification is required

ISOTECH

Typical Uncertainties of PRT Semi Standards with Range

Temperature

Uncertainty mK

				935-14-95H*	
Model		935-14-61*		935-14-72	935-14-95H*
	935-14-95L*	935-14-13	935-14-13*	935-14-16	935-14-72
-196	25	N/A	25	N/A	N/A
-80	20	N/A	20	25	25
-50	15	15	15	20	20
0	10	10	10	15	15
50	10	10	10	15	15
156	10	10	10	15	20
232	N/A	15	15	20	25
420	N/A	N/A	N/A	40	40
550	N/A	N/A	N/A	N/A	50
660	N/A	N/A	N/A	N/A	50

*Preferred Models

The above uncertainties do not include long term drift Typical Stability of correctly used semi standard is 0.01°C/year at 0°C Actual uncertainty of a probe determined at time of calibration

Isotech have generated a long history of many of our semi-standards.

Here are a few documented facts:

The 935-14-95 model has the widest temperature range and in consequence is likely to suffer the largest changes in characteristics.

Guy Snelling sent the following email about the 935-14-95.

UKAS Calibration

available for these systems

International Traceability - Best Practice

ISOTECH

I though that you might like to see the calibration history of one of our probes from the past 12 years.


You may recall that we purchased this probe to use as a laboratory standard when our company was still young. This particular probe is still in daily use and is regularly taken to 600°C in our dry block calibrator. While we handle it with care, being in daily use for 12 year it has take the occasional mild knock and accidental abuse - I believe that it was even taken to close to 700°C once, although I wasn't involved so I can't testify to the temperature reached.

You'll see from the attached history of the calibration by our NMI that the probe has remained stable and accurate, and bearing in mind the daily variations in temperature that it has undergone, these results are testimony to the high quality of this product.


John, you are to be congratulated on developing and producing such a fine measuring instrument, and feel free to use us as a product reference any time.

Kind regards,

Guy Snelling Temperature Metrologist InterCal (South Africa)

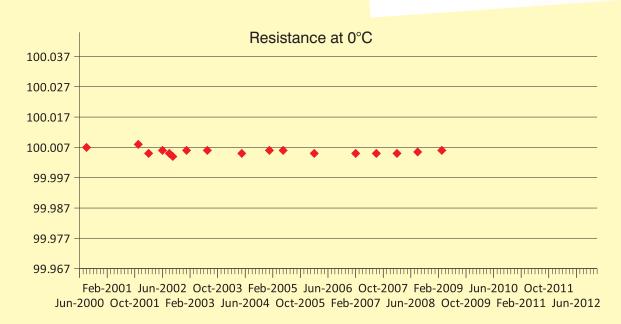



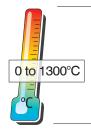

Reference Probes - Semi Standards **Thermocouples**

- Wide Temperature Ranges
- Noble Metal & Type N for best life, stability and reproducibility
- Can be supplied with UKAS calibration

These thermocouples are suitable for use as references in Isotech Dry Blocks and for use with temperature indicators. Details of our laboratory grade Standard Thermocouples with separate cold junctions can be found in our publication *"Solutions for Primary and Secondary Laboratories"*.

These semi standards are lower cost and suitable for a variety of industrial applications.

The 935-14-91 is constructed from Platinum and Platinum Rhodium alloys and can be used to 1300°C. Recommended for the Pegasus 1200 and general purpose applications. It has 1M of compensating cable terminated with a miniature thermocouple plug. The 935-14-88 is similar to the 14-91 but is made entirely from precious metals, with platinum wires all the way to the miniature plug.

There is a range of high quality mineral insulated metal sheathed (MIMS) Type N thermocouples. These devices are lower cost than the noble metal types and can be bent to a desired shape if required. They are suitable for use in Isotech Dry Blocks and for general purpose measurement and calibration applications.

The system accuracy or uncertainty will depend on the application and what instrument they are used with. The table shows the uncertainties that we can offer with optional UKAS calibration from our accredited laboratory.





Noble Metal Thermocouples

Platinum wire for best performance, ceramic sheath construction. Carry case included.

Model	Diameter	Length (A)	Range	Application	Туре	
935-14-91/R	5mm	300mm	0 to 1300°C	Pegasus	R	
935-14-91/S	5mm	300mm	0 to 1300°C	General Purpose	S	
Termination: 1M	extension cable t	to miniature plug				
935-14-88/R	5mm	300mm	0 to 1300°C	Working industrial standard	R	
935-14-88/S	5mm	300mm	0 to 1300°C		S	
Termination: 1M	platinum cable to	o miniature plug				

———— A —

Type N Thermocouples

Recommended base metal thermocouple, low cost metal sheathed.

Model	Diameter	Length(A)	Termination	Range	Application	Туре
935-14-63	3mm	300mm	1M Cable Miniature Plug	0 to 1200°C	Gemini 700 Jupiter 650	Ν
935-14-64	3mm	300mm	Miniature Plug	0 to 1200°C	General Purpose	Ν
935-14-65	3mm	300mm	Standard Plug	0 to 1200°C	General Purpose	Ν
935-14-66	3mm	500mm	1M Cable Miniature Plug	0 to 1200°C	General Purpose	Ν
935-14-67	3mm	500mm	Miniature Plug	0 to 1200°C	General Purpose	Ν
935-14-68	3mm	500mm	Standard Plug	0 to 1200°C	General Purpose	Ν



Isotech UKAS Calibration Uncertainties (k=2)

Item	Measured Quantity Instrument or Gauge	Temperature Range	(k=2) Best measurement capability expressed as an uncertainty (±)
1	Temperature Platinum Thermocouples	-50°C to 0°C 0°C to 50°C	0.5K 0.45K
		50°C to 660°C 660°C to 1100°C Above 1100°C to 1300°C	0.4K 0.7K 1.7K
2	Other Thermocouples	-196°C -80°C to 300°C Above 232°C to 420°C Above 420°C to 660°C Above 660°C to 1100°C Above 1100°C to 1300°C	0.3K 0.25K 0.3K 0.4K 0.8K 2.2K The latest schedule can be found on the Isotech website or at www.ukas.org.



UKAS Calibration available for these systems - International Traceability - Best Practice



Refer to Chart

Temperature Range Calibration Dimensions Refer to Chart A UKAS Calibration Certificate can be provided at extra cost Refer to Chart

How to Order Please Specify Model Type (for example 935-14-65) Please state whether UKAS Certification is required

ISOTECH

Introduction to Temperature Calibration Software

Software

Isothermal Technology's range of calibration software saves you time and lowers calibration costs. Isotech have a tested solution to calculate coefficients for industrial probes, for SPRTS, fit thermocouple error curves and fully automate the calibration of sensors.

ITS-90 Software

Icarus is software for the ITS-90 Laboratory to calculate between resistance and temperature for SPRTs. It allows for the calculation of coefficients and the printing of charts and certificates, see our publication "Solutions for Primary and Secondary Laboratories"

Cal NotePad

Isotech calibration equipment is supplied with Cal NotePad. This software allows equipment to be controlled remotely, monitored and the logging of data which can be imported into spreadsheet software. Cal NotePad supports the connection of both an Isotech Furnace, Bath, Block and an Isotech TTI Temperature Indicator.

I-Cal Easy Software

Both Cal Notepad and I-Cal Easy LOG are available at no cost to users of Isotech equipment.

I-Cal Easy LOG can be upgraded to the full version 'I-Cal Easy' to enable the automatic calibration of temperature sensors, from controlling the calibration run to printing certificates and calculating coefficients. When I-Cal Easy LOG is first installed it allows all the features of I-Cal Easy to be used over a 30 day evaluation period. After 30 days it can be used in LOG mode indefinitely.

Software Comparison Chart	Cal Notepad	I-Cal Easy LOG	I-Cal Easy
Connect to Isotech Equipment	✓	✓	\checkmark
Included with Isotech Equipment	✓	✓	Demo
Monitor and Record Data	✓	✓	✓
Supports Isotech Scanner and Switchboxes	×	✓	✓
Automatic Sensor Calibration	×	×	\checkmark
Print Certificates	×	×	\checkmark
Design Custom Certificates	×	×	\checkmark
Calculate Coefficients for IEC-751 (Cvd) and IST-90	×	×	\checkmark
Produce Thermocouple Tables	×	×	\checkmark



Demonstration versions of our software are available to download for evaluation:

http://www.isotechna.com



Calibration Software Cal NotePad

Easy to use

C. Note Pad

- Interface to Isotech Block Baths & Temperature Indicators
- Log Chart and Export data Control Calibration Bath - Read Standard

The purpose of Cal NotePad (CNP) is to automatically log and display the temperature of an Isotech calibration bath together with the unit under test. Cal NotePad can be used with baths (or indicators) without PC interfaces by the user typing in values from the keyboard. The Cal NotePad can be used to identify the operator and the unit under test. With the click of a button data is logged with time information, it is also possible to log continually.

The calibration bath temperature can be changed from the PC or from the calibration bath's controller - Cal NotePad will display the temperature changes as they occur on the re-scalable chart display.


For traceable calibration the unit under test should be compared to a calibrated standard thermometer. Cal NotePad can record the actual temperature of the bath from either the in-built indicator of an Isotech SITE model or from a variety of external instruments see list. If the external instrument has two channels e.g. Isotech TTI then the unit under test may be connected to channel B for logging with CNP. Alternatively the value can be typed in from the keyboard. Similarly the calibration bath controller value, actual temperature, SITE indicator value or unit under test value may also be entered manually.

Then the manually entered data is combined with that gathered automatically and the resultant file can be opened in an external application such as Excel for the preparation of reports, certificates etc.

Cal NotePad is designed for ease of use, it will give a chart of the system. When the operator determines the system is stable - easily seen from the chart, then clicking a button will record time, operator, serial numbers of unit under test along with controller and indicator values.

Cal NotePad can be used for semi-automatic calibration, see I-Cal Easy for a fully automated calibration solution.

http://www.isotechna.com

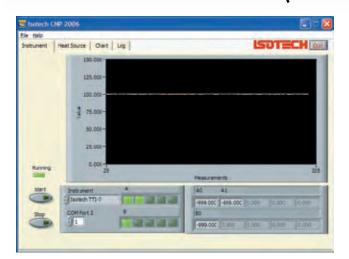

System Requirements Laptop or Desktop Windows 98/2000/XP/7 Serial Ports: RS232 or Adapaters

CNP is compatible with the following Isotech calibration equipment:-Calibration Baths, Furnaces, Dry Block and POTT models with a serial interface. Temperature Indicators: Isotech milliK, TTI-1, TTI-2, TTI-5, TTI-6, TTI-7, TTI-10 and TTI 22



I-CAL Easy

Calibration Software I-Cal Easy

- Fully Automatic Calibration
- Design and Print Certificates
- Calculate Coefficients
- Prints PRT and TC tables
- Supports more equipment
- Try the full version free for 30 days

Use I-Cal Easy to automate sensor calibration, enter up to 20 calibration points and let the software set the bath, wait for stability and log the data automatically. Choose the stability criteria and how many points to record at each calibration temperature. Automatic temperature calibration the easy way.

I-Cal Easy lets you use a built-in template or design your own certificate. Add text, data fields and graphics on single or multiple pages, then publish the calibration data to the certificate. Do you want to include or calculate coefficients? Then drag your data to the ITS-90 or Callendar Van Dusen calculators. For thermocouples use the powerful regression calculator to fit error curves.


Other systems have limited the user with a built-in template and the need to pay extra for any changes, with I-Cal Easy just build in your own certificate in minutes!

> Evaluate I-Cal Easy free for 30 Days. After the evaluation period the software continues to function as I-Cal Easy LOG which allows remote monitoring and logging. The advanced features including automatic calibration and certificate printing are disabled until a license is purchased but the remote monitoring and logging features remain available indefinitely.



http://www.isotechna.com

	SC 64 (E) Enos	e 15 kela	alle 10 Car	n Teachtie	MI YestEuks	Tel Hard Eastern	C Santa Mar	Plane 164
wite N	wints destroyment that	1						
	(Descelling	Chartest Times	Statement .	Med Merer	Sergeant.	Selgount Tol		No.
-	2211825	111,502,611	- IL BOUND	31.45.160	21.4014	1868	1008	Post
	++++ A28 -	111.020017	21.0003000	22-051200	22,9425	75.8	-	Post
-		THE STORES.	(1 Minute)	21.60206	20.0014	100	75.6	No. of
	+111828	10.821000	38 (#13963)	28.001000	20 1021	96,2	State -	
	11111037	101 230225	10.00200012	25 erroal.	21403	764	19.6	P.F.
-	et ti GI	111.822144	28 MIC2002	28.867380	24 4403	76.8	76.8	24
-	+4118328	101322005	21.000000	21.861.061	214618	794.5	75,8	Ph
-	-cingit.	111.022000	23 0000000	24 MILTER	20 4820	94.0	76.2	
-	+1116A88.	101.014010	21-0010001	21.011000	0.403	- 10	96.0	10
-	41111022	101.002000	25.0000000	23 (75,76)	20 48225	78.2	76.2	24
-	-110800	THE ADDRESS.	14-04030-001	25.824981	214029	12.8	10.0	-
-	+115827	111.023410.	ja entrete.	28.671085	34.8028	18.0	1.182	.pe
	+117850	101 820800	Distances .	23 512,001	12 4024	14.8	16.0	175
	arres all	TTO NUMBER	1110001	31.074.0	Casteria	147	74.0	
_	-CTMD	111.428222	11.07200000	23.40100	224828	71.8	75.8	Pe
	10000	10.24622	14 Arrester	21.462.00	20.6656	784	76.6	
-	1111812	111.000000	22.07030000	22.002701	22,0220	16.2	75.5	ma .
-	inen ditt.	111.224625	21-04230443	21.652541	(25-8615)	162	164	Pet
	10111424	111.874883	28.0700000	28.80 (98)	24 6028	76.6	764	24
_	1411828	101.540053	29-8053600	26 (07 00)	254004	796.0	196.0	-

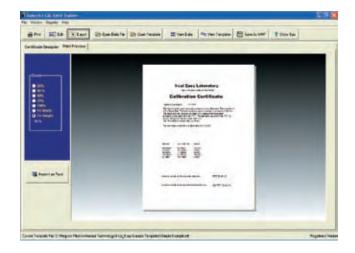




I-Cal Easy supports the TTI-7, TTI-22, microK and Isotech Dry Blocks, Liquid Baths and Calibration Furnaces. Additional support for other and third party instruments is available, contact Isotech for details.

I-Cal Easy provides a powerful but easy to use automatic calibration system. A graphical setup lets you drag and drop instruments and equipment onto the appropriate PC port - no need to create config files. In addition to the comprehensive manual balloon tips guide you as to the operation of each control. Once familiar with the system this balloon help feature can be turned off. The criteria for stability can be set to suit all types of equipment; Dry Blocks, Liquid Baths and High Temperature Furnaces. Once the system is stable choose how many measurements to take at each calibration point and have the average value appear on the certificate. Create one or multiple page certificates, as many as required to suit different customers and different types of calibration, Thermocouple, Industrial PRTs and SPRTs. Drag and drop data and text fields onto the certificate, link to logos and other graphic elements.

The in-built calculator will calculate coefficients for both IEC 751, ITS-90

and for thermocouples you can choose what order of regression to fit an error curve. Try the demo version and see how easy it is to drag data to the calculator and export the results straight to a certificate.

Judge for yourself how this compares to any other software. The demonstration version will run without restriction for 30 days and enable you to learn how to use I-Cal Easy and save time by rapidly producing certificates to your own requirements. Try it and see why we are confident that I-Cal Easy is the market leader.

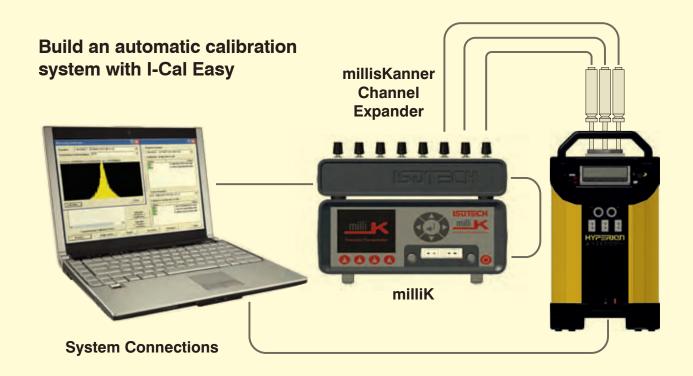

The number of sensors that can be calibrated depends upon the equipment used. Calibrate up to 16 sensors with two switch boxes, expandable to 32.

Svetom	Requiremen	ate
System	nequiremen	115

Computer
Operating System
Display
Optional Camera:
Serial Ports

Laptop or Desktop
Windows 98/2000/XP/Vista/7/8
1024 x 768 or greater
Check for supported types

Serial Ports A maximum of three ports are required. One for the Dry Block, one for the TTI and one for the switchboxes. (Two switchboxes can be operatede from a single port) **Note:** Available in different languages -Spanish and Chinese currently available.

http://www.isotechna.com

ISOTECH

Introduction to **Blackbody Sources**


This section contains dedicated blackbody sources for low uncertainty calibration of infrared thermometers. A range of portable primary blackbody sources combine high emissivity with excellent temperature uniformity. The cylindrical cavity design minimises the effects of air movement and ambient changes.

Many of the sources can be used with high purity ITS-90 Fixed Point cells where the thermometer is calibrated against the freezing temperature of a pure metal.

How To Calibrate Infrared Thermometers


1 With a Primary Standard Source

The temperature source has an emissivity approaching unity and sufficient uniformity so that the test thermometer can be compared to a traceable contact thermometer.

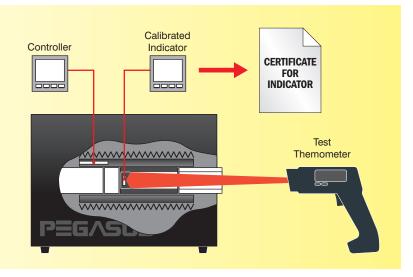


2 Calibration using a secondary standard sourceWith this method the blackbody is calibrated with a standard thermometer and the test thermometer is compared to the source. **Controller Controller**



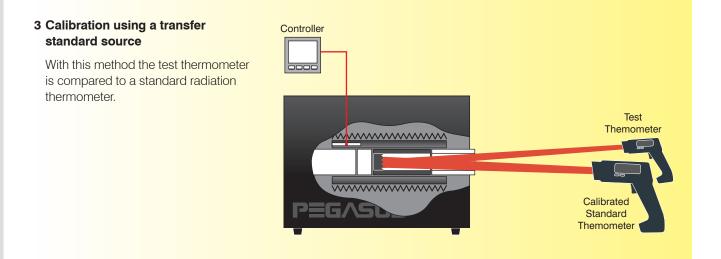

4 Using ITS-90 Fixed Point Cells

The test thermometer is calibrated, not against a source or other thermometer, but against a fixed temperature from an ITS-90 Fixed Point Cell. For example by melting a quantity of pure Gallium to obtain a fixed temperature of 29.7646°C

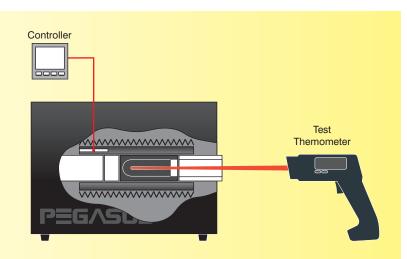

Blackbody

A blackbody has been defined as either a source with zero reflectivity or a source emitting the maximum possible radiation (at all wave lengths) for its temperature.

Emissivity

Emissivity is the ratio of the radiation emitted by a surface to that emitted by a black body at the same temperature.

Isotech has a range of sources having a high emissivity combined with thermal uniformity for use as Primary Standard Sources for low uncertainty calibration wavelength independent calibration.





Types of Equipment

Primary Black Body Sources


■ Hyperion R, Gemini R, Medusa R, Oberon R.

Black body sources covering the range from -10°C to 1100°C. Can also be used with fixed point cells. Aperture sizes ranging from 20 to 65mm.



Cyclops

Temperature range from 100°C to 1300°C.

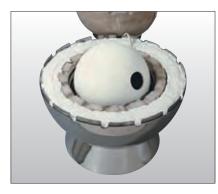

Special Applications


Model 988 Useful with Thermal Imaging Systems Covers Human Body Temperature

Model 975 Secondary Source Emissivity 0.95






Model 550-02 Low Cost 50 - 350°C

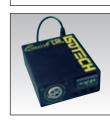


Calibration of Simple Low Cost Infrared Thermometers.

Many of Isotech's calibrators for PRTs and Thermocouples can also be used for testing simple Infrared devices, see Isocal-6 or Dry Block sections.

Cyclops Blackbody Source

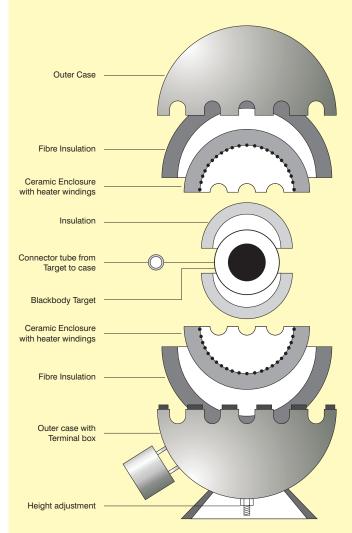

Effective emissivity of cavity in relation to temperature indicated by the thermocouple (TC)

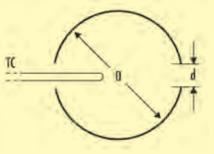




Isotech Blackbody Fixed Point Cells


Primary Standard Cells

Point	Part Number	Temperature	Apparatus
Indium	998-06-00A	156.60°C	Medusa R
Tin	998-06-00B	231.93°C	Medusa R
Zinc	998-06-00C	419.53°C	Medusa R
Aluminium*	998-06-00D	660.32°C	Oberon R
Silver*	998-06-00E	961.78°C	Oberon R
Copper*	998-06-00G	1084.62°C	Oberon R

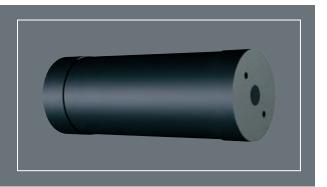


■ Medium Temperature "Hockey Puck" Cells

Point	Part Number	Temperature	Apparatus
Gallium	431-03-00	29.7646°C	Gemini R
Indium	976-05-00A	156.60°C	Gemini R
Tin	976-05-00B	231.93°C	Gemini R
Zinc	976-05-00C	419.53°C	Gemini R

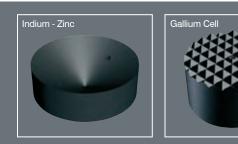


High Temperature Cells for the Pegasus R


Point	Part Number	Temperature	Apparatus
Indium	970-06-00A	156.60°C	Pegasus R
Tin	970-06-00B	231.93°C	Pegasus R
Zinc	970-06-00C	419.53°C	Pegasus R
Aluminium*	970-06-00D	660.32°C	Pegasus R
Silver*	970-06-00E	961.78°C	Pegasus R

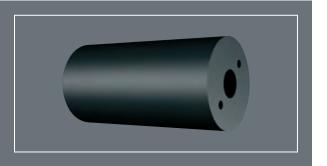


* These cells are required to be surrounded by gas for protection at high temperatures (See model 984)

Gallium Cells for Medusa R or Hyperion R to special order. Cells are provided with Certificate of Metal Analysis.

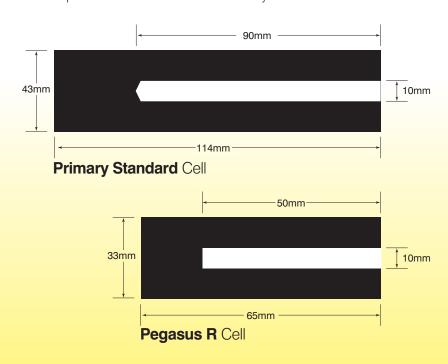

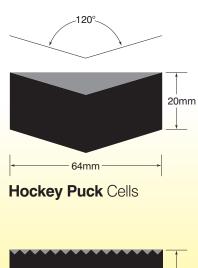

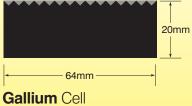



Gas Flow System Model 984

- Designed for Isotech Blackbody Cells
- Protect High Temperature Cells
- Gas flow interuption alarm

The higher temperature Isotech Blackbody Fixed Point Cells consist of high purity metals within a graphite body. Graphite reacts with air to form Carbon Dioxide. The rate of the reaction is temperature dependant. The effect is small at low temperatures but increases at higher temperatures.

For Indium and Tin cells the effect is small and for Indium, Tin and Zinc cells in general no attention needs to be made. Yet for Aluminium, Silver and Copper Cells the oxygen must be excluded or the cells would be quickly damaged.

The Model 984 Gas Flow System connects between an inert gas supply, such as Argon or Nitrogen and the Cell in its apparatus.

The Model 984 has a regulator and a flow meter optimised to easily set the flow to 12 litres/hour (0.2 l / min) and features an audible alarm should the gas flow be interrupted.





Model

Input - Output Connectors Maximum Flow Alarm Power Voltage Dimensions


984

s Genevac 16KF 16 litres/hour Audible 15W 12Vdc Height 240mm Width 120mm Depth 220m (excluding connection pipes) 2.5kg

Weight

How To Order Model 984 Gas Flow System Supplied with external power supply, 2 x connecting pipes and fittings



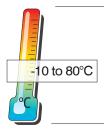

Blackbody Source Hyperion R

- Low Temperature Radiation Pyrometer Primary Source
- 50mm Diameter Cavity
- 0.995 Emissivity

The Hyperion R Portable Blackbody Calibration Source allows for calibration of noncontact infrared thermometers over the temperature range -10°C to 80°C.

It is suitable for use as a primary radiation source for infrared thermometers from sub zero to 80°C.


Laboratory performance and low uncertainty calibrations are ensured by the combination by high emissivity and excellent temperature uniformity.

The digital temperature controller allows the block temperature to bet set to any value from -10°C to 80°C.

Traceability of the radiance temperature is established by a separate, built-in temperature indicator and included platinum resistance thermometer.

A three point traceable calibration certificate is included. UKAS calibration of the resistance thermometer is available, as is radiometric calibration.

Uniformity of the block is ensured by using distributed thermoelectric heat pumps with the benefit of solid state, vibration free cooling.

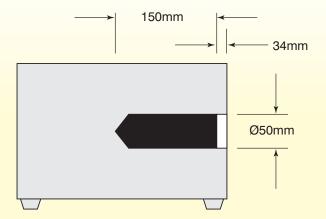

Hyperion R



10kg

100-130 or 208-240 Vac

Model	982
Temperature Range	-10°C to 80.0°C
Emissivity	Greater than 0.995
Stability	±0.1°C
Display Resolution	0.01°C
Heating Time	40 minutes to 80°C
Cooling Time	45 minutes to -10°C
Aperture Diameter	50mm
Cavity Depth	150mm
PC Interface	included
Power	200 Watts typical
Voltage	100-130 or 208-240
Dimensions	H 310mm
	W 265mm
	D 200mm

Weight

Options

Orifice Plates 10, 20, 30, 40mm 812-01-06 (Restricts Cavity Aperture) **Carrying Case** 931-22-64

How To Order

Model 982 Hyperion R State Supply Voltage Please state any special calibration requirement





Blackbody Source Gemini R

- 30°C to 550°C
- Emissivity > 0.995
- 65mm Diameter Cavity

The Gemini R 550 Portable Blackbody Calibration Source allows for calibration of noncontact infrared thermometers over the temperature range 30°C to 550°C.

It is suitable for use as a primary radiation source for infrared thermometers.

Laboratory performance and low uncertainty calibrations are ensured by the combination by high emissivity and excellent temperature uniformity.


The digital temperature controller allows the block temperature to bet set to any value from 30°C to 550°C.

Traceability of the radiance temperature is established by a separate, built-in temperature indicator and included platinum resistance thermometer.

A three point traceable calibration certificate is included. UKAS calibration of the resistance thermometer is available, as is radiometric calibration.

Uniformity of the block is ensured by using distributed heating technology.

For the smallest of uncertainties the Gemini R may be used with Isotech ITS-90 Fixed Point Cells, Gallium 29,7646°C. Indium 156.5985°C, Tin 231.928°C and Zinc 419.527°C. The cells are provided with a certificate of metal purity.

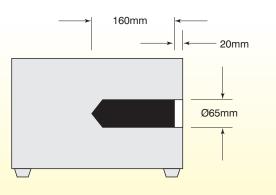

Gemini R

Test Report

The variation seen on the controller's temperature indication over a 5 minute period was $\pm 0.2^{\circ}$ C. Similar variations were detected by a radiation thermometer looking into the cavity.

Using a portable radiation thermometer having a target diameter of 13mm,

the 65mm target was surveyed. Maximum temperature differences of $\pm 1^{\circ}$ C were found for set temperatures in the range 100°C to 500°C.

Temperatures along the inner 100mm of the 160mm long cavity were measured at 400°C and 500°C, using a hand held fibre-optic radiation probe. Maximum temperature differences of $\pm 4^{\circ}$ C were found.


The temperature, as shown on the controller, agrees with the cavity temperature as measured by a radiation thermometer, where calibration is traceable to National Standards, to within $\pm 2^{\circ}$ C.



Model	976
Temperature Range	30°C to 550.0°C
Emissivity	Greater than 0.995
Stability	±0.1°C
Display Resolution	0.01°C to 99.99; 0.1°C from 100 to 550
Heating Time	45 minutes
Aperture Diameter	65mm
Cavity Depth	160mm
PC Interface	included
Power	1000 Watts typical
Voltage	100-130 or 208-240 Vac 50/60Hz
Dimensions	H 310mm W 265mm D 200mm
Weight	10kg
Options Fixed Point Cell	
Gallium Hockey Puck Cell	431-03-00
Indium Hockey Puck Cell	976-05-00A 976-05-00B
Tin Hockey Puck Cell Zinc Hockey Puck Cell	976-05-00D 976-05-00C
Orifice Plates 10, 20, 30, 40, 50mm (Restricts Cavity Aperture)	976-01-05
Carrying Case	931-22-64

How To Order Model 976 Gemini R State Supply Voltage Please state any special calibration requirement



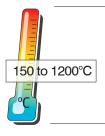

Blackbody Source Pegasus R

150°C to 1200°C

- Compact 20mm Diameter Cavity
- Emissivity > 0.995 Cavity 20 x 65mm

The Pegasus R is a compact furnace for calibrating radiation pyrometers.

The temperature of the furnace is set on a controller, whilst an independent indicator, whose sensor fits into the cavity, indicates the actual radiance temperature. The sensor can be removed for external calibration or the complete system can be calibrated.


The cavity diameter is 20mm, the depth 65mm. Emissivity is 0.995. The cavity is removable and a fixed point cell may be put in its place. The cavity inside the fixed point cell is 10mm in diameter by 65mm deep to the tip of a 120° cone.


Blackbody target radiation source for use with Pegasus R.

For calibration radiation thermometers in the wavelength range 0.9 to 14 micrometres. A blackbody radiation source 'based on a design study by England's National Physical Laboratory (NPL)' for Isothermal Technology Ltd is housed, with suitable insulation, in the Pegasus tube furnace. The aperture is 20mm in diameter.

A Pegasus blackbody source has been calibrated at NPL (Nat. Physical Laboratory) with an uncertainty of $\pm 2^{\circ}$ C and the calibration was found to be reproducible after a period of about 2 months. A scan across the aperture at 444°C showed that the source was uniform to better than 0.3°C.

Traceability may be established with a UKAS certificate for the in-built indicator and supplied probe (935-14-40).

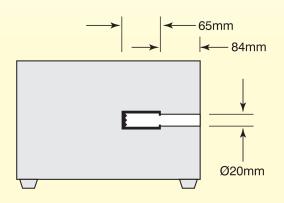

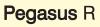



Model

Temperature Range Emissivity Stability Display resolution

Cavity size

Heating Rate PC Interface Power

Dimensions

Weight

Options

Indium Blackbody Cell Tin Blackbody Cell Zinc Blackbody Cell Aluminium Blackbody Cell Silver Blackbody Cell Probe Carrying Case Gas Flow System


How to order

Model 970 Pegasus R Please state supply voltage required Please state any special calibration requirement

970 150°C to 1200°C 0.995 ±0.1°C 0.1°C to 999.9; 1°C from 1000 to 1200 20mm diameter 65mm deep 25°C/minute Included 800W typical 100-130 or 208-240 Vac 50/60Hz Height 310mm Width 265mm Depth 200mm 13 kg

970-06-00A 970-06-00B 970-06-00C 970-06-00D 970-06-00E 935-14-40 931-22-64 984-00-00



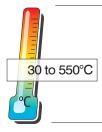

Blackbody Source Medusa R

- 30°C to 550°C
- Emissivity > 0.995 Cavity 45 x 285mm
- Accepts Isotech Primary Blackbody Cells

The Medusa R Blackbody Calibration Source allows for calibration of noncontact infrared thermometers over the temperature range 30° C to 550° C.


It is suitable for use as a radiation source for infrared thermometers. The cavity is 45 x 285mm deep and suitable for use with the larger Isotech fixed point cells.

Laboratory performance and low uncertainty calibrations are ensured by the combination by high emissivity and excellent temperature uniformity.

The digital temperature controller allows the block temperature to bet set to any value from 30°C to 550°C.

Uniformity of the large block is ensured by using distributed heating technology.

For the smallest of uncertainties the Medusa R may be used with Isotech ITS-90 Fixed Point Cells, Gallium 29.7646°C, Indium 156.5985°C, Tin 231.928°C and Zinc 419.527°C. The cells are provided with a certificate of metal purity.

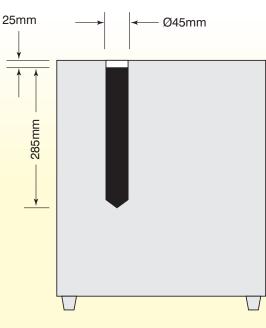

Medusa R



999

Model

Temperature Range Emissivity Stability Display Resolution

Heating Time Aperture Diameter Cavity Depth PC Interface Power Voltage

Dimensions

Weight


Options Fixed Point Cells Indium Large Primary Cell Tin Hockey Large Primary Cell Zinc Large Primary Cell

How To Order Model 999 Medusa R State Supply Voltage Please state any special calibration requirement

30°C to 550.0°C Greater than 0.995 ±0.1°C 0.01°C to 99.99; 0.1°C from 100 to 550 45 minutes 45mm 285mm included 1000 Watts typical 100-130 or 208-240 Vac 50/60Hz H 480mm W 425mm D 260mm 17kg

998-06-00A 998-06-00B 998-06-00C



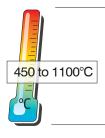

Blackbody Source **Oberon R**


- 450°C to 1100°C
- For High Temperature Blackbody Fixed Points
- Utilises a Sodium Heatpipe

The Oberon R uses a Sodium Heat Pipe to ensure an exceptionally low temperature gradient along the furnace core. It is ideal for the realization of Aluminium, Silver or Copper ITS-90 Fixed Points.

It may be used as a blackbody source over the range 450° C to 1100° C.

The furnace heater is of the non-inductive bird-cage design insulated by twin bore alumina tubes. The heatpipe is designed so that the inner wall is not subject to thermal expansion stresses from the outer wall before the heat pipe reaches conduction temperature. The working fluid is permanently and safely sealed within the plasma-arc welded enclosure.



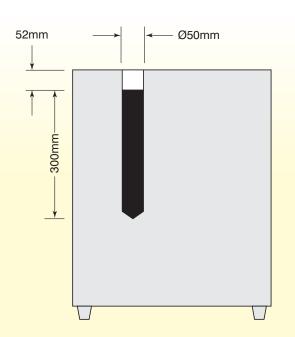

Oberon R


Nodel	426	
emperature Range	450°C to 1100°C	
Emissivity	greater than 0.995	5
Stability	±0.05°C	
Display resolution	0.1°C to 999.9; 1°C from 1000 to	1090
Cavity size	50mm diameter 300mm deep	
ime to temperature	4 hours	
PC Interface	Included	
Supply	110Vac, 3kW, 50/6 (230Vac and 110\ Isolating Transform	/ac to 110Vac
Dimensions	Height 410mm Width 415mm Depth 280mm	
Veight	30.5kgs	
Options		
Aluminium Primary Blackl Silver Primary Blackbody Copper Primary Blackboo Gas Flow System	Fixed Point Cell	998-06-00D 998-06-00E 998-06-00G 984-00-00
230v/110v Transformer		935-19-43
10v/110v Transformer		935-19-48

How to order

Model 426 Oberon R Please state voltage required

Please state any special calibration requirement



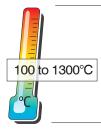


Blackbody Source Cyclops

- Spherical Blackbody Source
- Wide Temperature Range 100°C to 1300°C
- Can be adapted for Thermocouple Calibration

The Cyclops Model 878 is a spherical blackbody source. It consists of an inner black sphere that sits inside a spherical furnace and is suitable for use as a radiation source for infrared thermometers.


The inner sphere has a nominal diameter 230mm and is accessed by an optical sighting tube. The furnace can be supplied in one of two constructions, one providing an aperture size of 17mm and one of 45 mm.

The furnace can be adapted for thermocouple operation by replacing the inner sphere with an equalizing block and using a different control sensor.

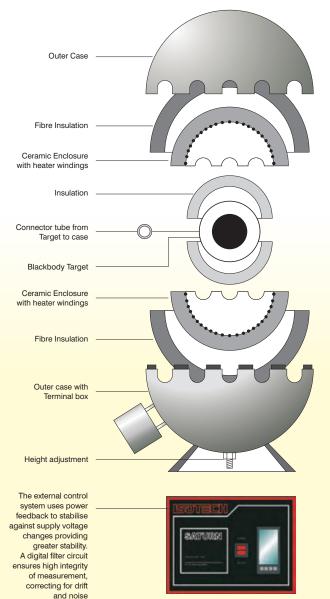



Cyclops Assembly Diagram Plan View (shown in section)

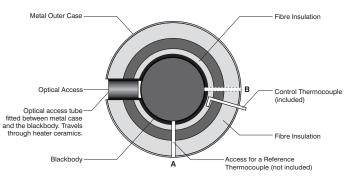

Model Temperature Range Emissivity Stability Display resolution

Time to temperature

PC Interface Power Voltage

Dimensions Weight **878** 100°C to 1300°C Greater than 0.999 ±0.1°C 0.1°C to 999.9; 1°C from 1000 to 1300 90mins hour to 700°C 4 hours to 1300°C included 3kW typical 100-130 or 208-240 Vac 50/60Hz 425mm Diameter 25 kg

Options

Ceramic Equalising Block to accept up to 8 thermocouples 878-02-08.

How to Order

Model 878 Cyclops. Please state supply voltage required Please state target diameter either 17mm standard or 45mm to special order.

Reference Thermocouple Access Point is at position A as Standard, specify if required at point B. State Target Diameter 17 or 45mm





Blackbody Source Model 988


20°C to 45°C


- Emissivity better than 0.97 ±0.02
- Controller Resolution 0.01°C

This blackbody source has been introduced to meet the demand for a simple, cost effective but high accuracy calibrator for the calibration of thermal imagers and infrared thermometers used at temperatures around ambient.

A 70mm diameter ridged plate is heated or cooled with an internal solid state thermoelectric heat pump. The temperature of the plate can be set from 20° C to 45° C to a resolution of 0.01°C.

Evaluation showed the advantages of fitting a stainless steel tube around the plate to give better uniformity and less sensitivity to draughts and ambient temperature effects.

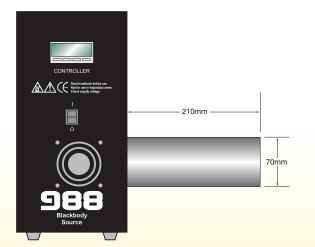

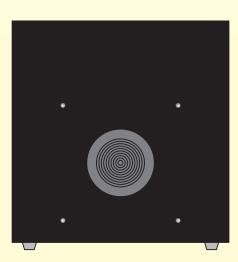



http://www.isotechna.com

Model Temperature Range Resolution Target Size Emissivity Combined Accuracy / Stability Power Voltage Dimensions

Weight

Optional PC Interface

Accessories Switch Mode Power Supply Supplied as Standard 100 - 240 VAC


How to Order Model 988 20°C to 45°C ±0.01°C 70mm Diameter 0.97 ± 0.02 ±0.2°C (±0.3°F)

60 Watts

988

12 Vdc H 230mm W 225mm D 115mm 4kg



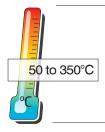


Greybody Source Model 975

- 50°C to 350°C
- Emissivity > 0.95
- 70mm Ridged Plate Target

When the high accuracy of the Gemini R is not necessary this product offers a cost effective solution for the calibration and testing of infrared thermometers.

The Greybody Model 975 gives fast accurate results with a larger sensing area. A temperature sensor sits just under the target surface and controls the temperature of the source. A custom designed surface sensor is used to set the controller calibration and a traceable certificate is supplied with each source.

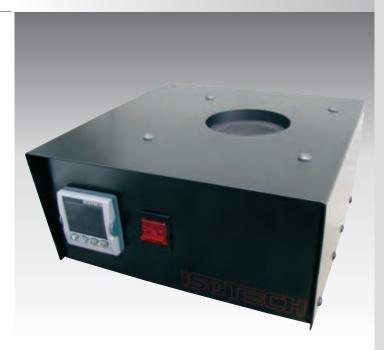

http://www.isotechna.com



Model

Temperature Range Display Resolution

Heating time Target Size Stability Accuracy Emissivity Power Voltage PC Interface

Dimensions

Weight

How to Order Model 975 Greybody Source Please state voltage required

975 Greybody 50°C to 350°C 0.01°C 50 to 99.99 0.1°C 100 to 350 35 minutes Ridged Plate, 70mm Diameter ±0.2°C ±2 >0.95 180 Watts 100-130 or 208-240 Vac Included H 115mm W 230mm D 225mm 3.9kg






Blackbody Source QuickCal Blackbody

30°C to 350 °C

- Emissivity > 0.99
- Controller Resolution 0.1°C

When the high accuracy of the Gemini R is not necessary this product offers a cost effective solution for the calibration and testing of infrared thermometers. The Quick Cal Blackbody is ideal for rapid and portable checking of infrared thermometers.

It features a cylindrical cavity 25mm diameter by 115mm deep.

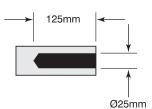

Quick-Cal



931-22-100 Optional Carrying Case

DANGER HOT 1



http://www.isotechna.com

Model

Temperature Range Display Resolution Heating time Target Size Stability Accuracy Emissivity

Power Voltage Dimensions

(F

TEMPERATUR

Weight

±0.2°C ±0.5 >0.99 Surface coating 0.98 - cavity gives overall emissivity of >0.99 300 Watts 100-130 or 208-240 Vac H 65mm W 152mm D 175mm 1.5kg

550 QuickCal Blackbody

25 x 115mm Cavity with end cone

30°C to 350°C

0.1°C

9 minutes

How to Order

Model 550-02 Blackbody Source Please state voltage required Optional Carry Case 931-22-71



Thermocouple Referencing Techniques

Isotech has a world leading range of thermocouple referencing equipment. From laboratory models for standard thermocouples through to large scale installations used in power stations, aeronautical, industrial and research institutes.

Types of Equipment

Automatic Ice Point

The junctions are maintained at a fixed temperature of 0°C

Features

- Provides 0°C reference as adopted by thermocouple tables
- Reference from one to a 100 Junctions
- Will reference any type or combination of thermocouple types
- Can be bench or rack mounted.

Constant Temperature Ovens

The junctions are maintained at a fixed temperature typically in the range 40 to 75°C

Features

- The junctions are maintained at a fixed elevated temperature.
- Reference up to 100 Junctions
- Will reference any type or combination of thermocouple types
- Can be provided as bench, rack mounted or in a weatherproof wall mounting case to IP66

Isothermal Systems

The junctions are not maintained at a fixed temperature, but are held in a thermal reservoir with an output signal which is fed to the instrumentation system.

Features

- The junctions float at ambient temperature in an isothermal block - no loading errors
- Can be fully passive with no power requirement
- Will reference any type or combination of thermocouple types
- Can be provided as bench, rack mounted or in a weatherproof wall mounting case to IP66

Types of Housing





Bench Mounting

Wall Mounted




Rack Mounted

Types of Junction


Single Junction

This is particularly useful with Type T thermocouples. With Type T one leg is copper, so only the Copper / Constantan wire needs to be processed thermally.

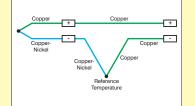

Double Junction


This is the most commonly used method. The input signal being connected to the "double junction" with the output on copper wires.

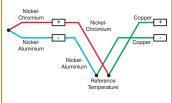



Isotech Custom Engineering

Isotech have more than 30 years experience in consulting and designing referencing systems. With a reputation for outstanding reliability, long term support and excellence Isotech can provide systems from a single, to several thousand junctions.

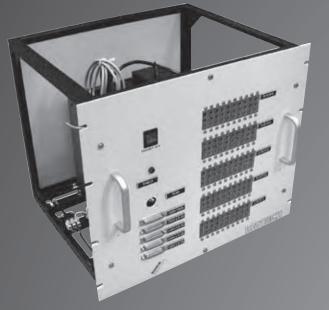

Solutions for Aeronautical, Power Stations, Environmental Monitoring, Space, Boiler Rooms, Maritime...












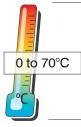

Thermocouple Probes Cold Junction

- Looms for Compensation Units
- Saves on On Site Wiring
- Maintains Accuracy of Unit

Probes are suitable for use with all lsotech thermocouple reference units or may be used with any other equipment including ice flasks and ice point reference units.

They can be supplied to suit a single thermocouple, or up to ten double junctions in a single assembly.

These probes are normally specified to order, or are made to match the reference equipment to which they will be fitted. Normal lead length, L2, is 1M but can be specified at the time of order along with the probe length, L1.

All wire material is to the highest grade available and PTFE insulated, numbered and colour coded for the appropriate thermocouple type.

Double junctions are most commonly supplied (four wire connections per junction). Single junctions (two wire connections per junction) can also be supplied.


For Types R & S cold junctions are most commonly made from low cost compensating cable but can be supplied in platinum material to special order.

Standard thermocouple Types are, K, E, J, T, N, U, S and R.

Other materials available on request

For Isotech reference units simply advise which model the junctions are for, and the types and number required. For custom junctions the following information is required.

Number of Junctions Required	
Thermocouple Type	
Single or Double Junction	
Length of Probe, L1	
Length of Lead, L2	
Maximum Probe Diameter	

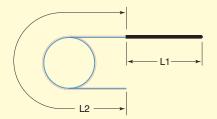



Standard Combinations Available

Code	Single Junction Probes
К	Nickel Chromium vs Nickel Aluminium
Е	Nickel Chromium vs Copper Nickel (Constanta
J	Iron vs Constantan
Т	Copper vs Constantan
N	Nicrosil vs Nisil
U	Copper vs Cupronic
S	Platinum vs Platinum 10% Rhodium
R	Platinum vs Platinum 13% Rhodium
Code	Double Junction Combinations examples
K	Nickel Chromium vs Copper
V	Niekel Aluminum ve Conner

an)

- K Nickel Aluminum vs Copper
- J Iron vs Copper
- J Constantan vs Copper
- U Cupronic vs Copper
- S/R Platinum vs Copper
- S Platinum 10% Rhodium vs CopperR Platinum 13% Rhodium vs Copper

("U" is a substitute metal alloy combination for Pt/Pt Rh types in the range 0 to 50° C).

Other materials are available on request.

How to order

Model 880 Cold Junction Probes are normally specified for each order. Please discuss your exact requirements with us before ordering.





Thermocouple Reference Unit TRU Model 938

- Suitable for Laboratory or High Capacity Applications
- Works in high ambients up to 65°C
- Reliable Solid State Design


The TRU (Thermocouple Reference Unit) supplies a stable and accurate 0°C, or elevated reference temperature.

It is a self-contained all solid-state unit using Peltier technology which provides maintenance free operation.

The TRU features rapid temperature shift even from high ambient temperatures and is stable within 10 minutes from switch on.

An alarm will be activated should the reference temperature deviate by more than 0.2°C.

The 938 can be ordered with a choice of one of two block types. B1 is recommended for laboratory use and has 6 pockets 6.2mm x 130mm deep and a 4.2mm pocket and can accommodate up to 36 junctions. Block B2 is for higher capacity applications and can accommodate up to 100 junctions with 8 8.2mm pockets x 76mm deep and a 4.2 pocket for an optional monitoring PRT.

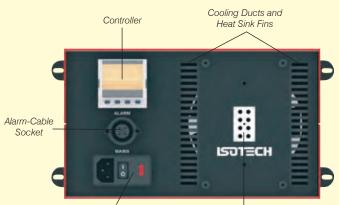

Switched Mains Cable Socket & Fuse-Holder Metal Block with Pockets for Thermocouple Junctions





Model		938
Operating 1	Temp.	0°C (or 45° to 70°C)
Ambient Ra	inge °C	2°C to 65°C
Stability		±0.03°C, Errors introduced by thermocouple loading can be removed by adjusting controller offset
Stabilising	Time	10 minutes from 44°C
Capacity		B1 6 x 6.2mm Pockets + 4.2mm pocket 130mm deep or B2 8 x 8.2mm Pockets +4.2mm pocket 76mm deep.
Alarm Facil	ities	Non-latching relay rated 5 Amps 240V
Power		100 Watts typical 100-130 or 208-240 VAC 50/60Hz
Dimensions	5	Height 228mm Width 253mm Depth 148mm
Weight		5.5kg
Accessorie	es	
935-14-54		Resistance Thermometer suits Block des UKAS Calibration at 0.01°C
935-14-55		Resistance Thermometer suits Block des UKAS Calibration at 0.01°C
935-17-32	Fan Filter environm	r - recommended for high dust nents
How to ord TRU Model	938	offed for each order

Normally uniquely specified for each order. Please discuss your exact requirements with us before ordering.





Thermocouple Reference Unit TRU Model 937

- Up to 100 Channels
- Compact, Pre Wired Thermocouples
- Operates in High Ambients

The TRU 937 (Temperature Reference Unit) supplies a stable and accurate 0°C or elevated Reference Temperatures between 45°C and 70°C.

It is a self-contained all solid state unit using peltier technology which provides maintenance free operation. The TRU 937 features rapid cool down from high ambient temperatures and is stable within 10 minutes from switch on.

An alarm will be activated should the reference temperatures deviate by more than the user definable span. Thermocouple Reference Junctions are located in a reference block and connected to their marked input and output terminals in an isothermal enclosure. The uniform temperature throughout the enclosure ensures that no thermoelectric EMFs are generated at the terminals.

One advantage of the TRU 937 is that the user need not be concerned with the supply and installation of reference junctions since it is only necessary to connect the thermocouple compensation cables to the input terminals and the measuring instrument to the output terminals of the TRU 937.

There are two models, the TRU 937/50 with up to 50 junctions with a single terminal cover door, and the TRU 397/100 which has doors on the front and rear, as shown in the photograph.



Note:

Thermocouple Referencing to ISO9000 at 0°C in Ambients up to 65°C or elevated reference temperatures between 45°C and 70°C References up to 100 Thermocouples Pre-wired Thermocouples, Compact Design



Model		937	
Operating T	emp.	0°C (or 45° to 70°C)	
Ambient Ra	nge	2°C to 65°C	
Stability		±0.03°C, Errors introd thermocouple loading by adjusting controller	can be removed
Stabilising 7	Гime	10 minutes from 44°C	
Capacity		Up to 100 Double Jun	ction Channels
Input/Outpu Connection		Klippon Terminals, typ	e 1.5 AKZ
Alarm facilit	ies	Non-latching relay rate	ed 5 Amps 240V
Power		100 Watts typical 100-130 or 208-240 VA	AC 50/60Hz
Dimensions		Height 265mm Width 253mm Depth TRU 937/100 Depth TRU 937/50	312mm 230mm
Weight		TRU 937/100 TRU 937/50	11kg 8kg
Accessorie			
935-14-55		Resistance Thermome des UKAS Calibration a	
935-17-33	Fan Filter environm	- recommended for higents	gh dust
	937 niquely spo	ecified for each order. αct requirements with ι	is before ordering.





Thermocouple Reference Unit TRUrac Model 847

С

- Rack Mounted
- Large Capacity
- Approved for Power Station Use

The TRUrac is a 0°C or elevated temperature thermocouple reference system mounted in a 19" chassis. It has been developed for situations where ambient temperature can be up to 65°C.

The reference temperature is normally set to 0°C or between 45°C and 70°C. For other temperatures please contact lsotech.

An alarm will be activated should the reference temperatures deviate by more then 0.2°C.

Inside the rack case is a high stability thermal block which has a capacity of up to 100 thermocouple channels, the probe wires being terminated at the rear of the unit on rail mounted screw terminals.

To special order a second thermal block may be fitted to allow a capacity of 200 channels in a single unit.

The customer simply connects their thermocouple wires and copper output wires to these terminals. All the thermocouple cold junctions are inserted into a metal oven block which is accurately temperature controlled.




lodel		847
Operating T	emp.	0°C (or 45° to 70°C)
Ambient Ra	nge	2°C to 65°C
Stability		±0.03°C, Errors introduced by thermocouple loading can be removed by adjusting controller offset
Stabilising 7	Гime	10 minutes from 44°C
Capacity		Up to 100 Double Junction Channels
nput/Outpu Connection		Klippon Terminals, type 1.5 AKZ
Alarm facilit	ies	Non-latching relay rated 5 Amps 240V
Power		100 Watts typical 100-130 or 208-240 VAC 50/60Hz
Dimensions 50 to 100 C		Height 400mm Width 483mm Depth 312mm
Veight		24kg
Accessorie 035-14-54	Platinum	Resistance Thermometer UKAS Calibration at 0.01°C
935-17-32	Fan Filter environm	- recommended for high dust ents
How to ord Should be s		iniquely on each order.

Rack mounted Temperature Thermocouple Referencing System Large Capacity. Approved for Power Station Use. Reference temperatures set to 0°C or between 45°C and 70°C.

Note:

Should be specified uniquely on each order. Please discuss your exact requirements with us before ordering.



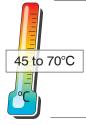

Thermocouple Reference Unit Hotbox Model 830

- Water Proof Cased
- Large Capacity up to 200 junctions
- Suitable for Power Station Use

The Hotbox is housed in a robust waterproof enclosure to IP66 incorporating bottom gland plate. Fixing lugs for wall mounting are provided. Easy access to terminal rails and oven assembly is via a hinged door.

All the thermocouple cold junctions are inserted into a metal oven block which is accurately temperature controlled.

A safety switch is fitted to cut off supplies if the temperature rises 10°C above the set point.



Model	830
Referencing Temp	45°C - 70°C
Accuracy	±0.1°C
ong Term Stability	±0.05°C per 1000 Hours
Temperature Gradient	±0.1°C between Junctions
Stabilisation Time	120 minutes
Max Ambient Temperature	50°C
Thermocouple	0 to 100 channels
Capacity	Double junction referencing
Power	250 Watts typical 100-130 or 208-240 VAC 50/60Hz
Dimensions	Height 600mm Width 600mm Depth 300mm
Veight	40kg

Note:

Numerous special versions are available and can be supplied either in their existing form or modified to customers' requirements.

Accessories


935-14-08

Platinum Resistance Thermometer. Includes UKAS calibration.

How to order

Model 830 Hotboxes are normally uniquely specified for each order. Please discuss your exact requirements with us before ordering.



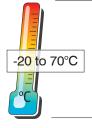

Thermocouple Reference Unit **Isobox Model 842**

- Water Proof Cased
- Large Capacity
- Approved for Power Station Use

The function of this reference unit differs from other cased systems in that the temperature of the metal block in which the thermocouple and copper leads are inserted, actually varies with ambient temperature. A separate output signal is produced which is proportional to the difference between the environmental temperature and the actual reference temperature. The output signal can be fed directly into a computer/data logger to give accurate compensation over a large ambient range. A thermal reservoir, heavily insulated, contains the reference junction probes. The reservoir temperature slowly follows the ambient temperature; an electrical compensation device is thermally integrated with the reservoir and thus senses the reservoir temperature. The device produces an output proportional to the difference between the reservoir temperature and the reference temperature (usually 0°C).

This is the signal the computer/data logger uses to compensate for the temperature of the reference probes junctions.

The output signal can be in the form of a DC mV output, 4-20 mA or from a platinum resistance thermometer. Please discuss prior to order.

Units are housed in robust weatherproof enclosures to IP66 incorporating bottom gland plate. Fixing lugs for wall mounting are provided. Easy access to terminal rails and block assembly is via a lockable hinged front door.



del	842
erencing Temp	Effectively 0°C
curacy	±0.1°C per 10° ambient span
ng Term Stability	±0.05°C per 1000 Hours
nperature Gradient	±0.1°C between Junctions
bilisation Time	10 minutes
bient Temperature	-20°C to +70°C
ermocouple	0 to 100 channels
pacity	Double junction referencing
ver	10 Watts typical 100-130 or 208-240 VAC 50/60Hz Low level D.C. Consumption 6VA typical Passive Option Available
nensions	Height 600mm Width 600mm Depth 300mm
ight	40kg

Note:

Dim

Wei

Мо

Ref

Acc

Lon

Terr

Sta

Am

The Cap Pov

Numerous special versions are available and can be supplied either in their existing form or modified to customers' requirements.


How to order

Model 842 Isoboxes are normally uniquely specified for each order. Please discuss your exact requirements with us before ordering.

Note:

Waterproof Cased Ambient Temperature Thermocouple Referencing System. Large Capacity up to 100 Thermocouples. Wide Ambient Range. Approved for Power Station Use.



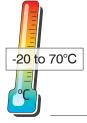

Thermocouple Reference Unit Isorac Model 844

- Rack Mounted Ambient Referencing System
- Large Capacity
- Approved for Power Station Use

The function of this reference unit differs from other rack mounted systems in that the temperature of the metal block in which the thermocouple and copper leads are inserted actually varies with ambient temperature. A separate output signal is produced which is proportional to the difference between the environmental temperature and the actual reference temperature. The output signal can be fed directly into a computer/data logger to give accurate compensation over a large ambient range.

A thermal reservoir, heavily insulated, contains the reference junction probes.

The reservoir temperature slowly follows the ambient temperature; an electrical compensation device is thermally integrated with the reservoir and thus senses the reservoir temperature.

The device produces an output proportional to the difference between the reservoir temperature and the reference temperature (usually 0°C). This is the signal the computer/data logger uses to compensate for the temperature of the reference probes junctions.

The output signal can be in the form of a DC mV output, 4-20 mA or from a platinum resistance thermometer. Please discuss prior to order.



Model No.	844
Reference Temperature	Effectively 0°C
Temperature Gradient	$\pm 0.1^{\circ}$ C between Junctions
Stability	±0.05°C per 1000 Hours
Ambient Range	-20°C to +70°C
Accuracy	$\pm 0.1^{\circ}$ C per 15°C ambient span
Stabilisation Time	10 mins
Thermocouple Capacity	0 to 100 channels double junction referencing
Power	10 Watts typical 100-130 or 208-240 VAC 50/60Hz Low level D.C. Consumption 6VA typical Passive Option Available
Dimensions	
Up to 50 channels	Height 255mm Width 483mm Depth 312mm
Weight	17.2kg
50 to 100 channels	Height 309mm Width 483mm Depth 312mm
Weight	20.4kg

Numerous special versions are available and can be supplied either in their existing form or modified to customers requirements.

How to order

Model 844 Isoracs are normally uniquely specified for each order. Please discuss your exact requirements with us before ordering.

Note:

Rack Mounted Ambient Temperature Referencing System. 19" Rack Mounting. Large capacity up to 100 Thermocouples. Approved for Power Station Use.




Reference Chart

Thermocouple Tolerances - IEC 60584 2013

B	Class 1 Class 2 Class 3	- 600 600	- +1700°C: +1700°C:	- ±0.0025 . t ±0.005 . t	or ±1.5°C or ±4.0°C
E	Class 1	- 40	+800°C:	±0.004 . t	or ±1.5°C
	Class 2	- 40	+900°C:	±0.0075 . t	or ±2.5°C
	Class 3	- 200	+40°C:	±0.015 . t	or ±2.5°C
J	Class 1 Class 2 Class 3	- 40 - 40 -	+750°C: +750°C: -	±0.004 . t ±0.0075 . t -	or ±1.5°C or ±2.5°C
K	Class 1	- 40	+1000°C:	±0.004 . t	or ±1.5°C
	Class 2	- 40	+1200°C:	±0.0075 . t	or ±2.5°C
	Class 3	- 200	+40°C:	±0.015 . t	or ±2.5°C
N	Class 1	- 40	+1000°C:	±0.004 . t	or ±1.5°C
	Class 2	- 40	+1200°C:	±0.0075 . t	or ±2.5°C
	Class 3	- 200	+40°C:	±0.015 . t	or ±2.5°C
R	Class 1 Class 2	0 0	+1600°C: +1600°C:	±1 for t<1100°C ±[1+ 0.003 x (t-1100)] for t>1100°C ±0.0025 . t	or ±1.5°C
S	Class 1 Class 2	0 0	+1600°C: +1600°C:	±1 for t<1100°C ±[1+ 0.003 x (t-1100)] for t>1100°C ±0.0025 . t	or ±1.5°C
Т	Class 1	- 40	+350°C:	±0.004 . t	or ±0.5°C
	Class 2	- 40	+350°C:	±0.0075 . t	or ±1.0°C
	Class 3	- 200	+40°C:	±0.015 . t	or ±1.0°C

t = temperature, degree Celsius, larger value applies

Thermocouple Colour Coding - IEC 60584-3

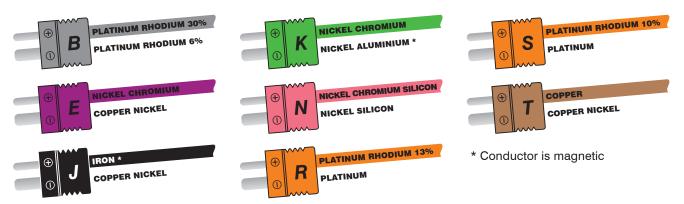



About the Laboratory

Isotech's UKAS accredited calibration laboratory, Northern Temperature Primary Laboratory (NTPL), was established in 1980 and has grown to be a full scale laboratory providing calibration to the smallest of uncertainties. Isotech was the first UKAS laboratory to be accredited to calibrate ITS-90 Fixed Point Cells. Our accredited uncertainties are now smaller than many the scope of other accredited laboratories and smaller than most National Measurement Institutes.

NTPL comprises of three physically separate laboratories, A Primary Laboratory, A Secondary Laboratory and Calibration laboratory for Industrial Products

Summary of Accreditation

Electrical

Bridges and similar instruments, resistance Calibrators, temperature simulation Resistance boxes Resistors, AC Resistors, DC Temperature indicators, electrical calibration Voltmeters, DC

Temperature

Block calibrators Fixed point cells Resistance thermometers, calibration by comparison Resistance thermometers, fixed point calibrations Temperature indicators and recorders, electrical calibration without sensor Temperature indicators and recorders, with temperature sensor(s) Thermocouples, base metal types, e.g. K, N, T Thermocouples, gold/platinum Thermocouples, platinum/rhodium types, e.g. S, R



http://www.isotechna.com

NTPL calibrate both Isotech equipment and devices from other manufacturers to the smallest of calibration uncertainties. The latest UKAS schedule can be downloaded from our website, a brief summary follows.





■ Calibration of Thermometers

Measured Quantity Instrument or Gauge	Range	Best Measurement Capability Expressed as an Expanded Uncertainty (<i>k</i> =2)	Remarks
Calibration by comparisons	-80°C to -40°C -40°C to 50°C 50°C to 156°C 156°C to 300°C 300°C to 420°C 420°C to 660°C	7.0 mk 4.0 mk 5.0 mk 6.5 mk 20 mk 35 mk	
Calibration at Fixed Points (See No	ote 1)		
BP Nitrogen TP Argon TP Mercury TP Water <i>(See Note 2)</i> MP Gallium FP Indium FP Tin FP Zinc FP Aluminium FP Silver	-195.798°C -189.3442°C -38.8344°C 0.01°C 29.7646°C 156.5985°C 231.928°C 419.527°C 660.323°C 961.78°C	5 mK 0.50 mK 0.24 mK 0.07 mK 0.15 mK 1.0 mK 1.0 mK 1.2 mK 2.0 mK 7 mK	Note 1: Suitable only for HT/SPRTs with high stability. Includes extrapolation to zero power and immersion checks. Note 2: Suitable for most SPRTs using nominal current.
■ Calibration of Thermoco	ouples		
Measured Quantity Instrument or Gauge	Range	Best Measurement Capability Expressed as an Expanded Uncertainty (<i>k</i> =2)	Remarks
Platinum thermocouples			
Calibration by comparisons	-50°C to 0°C 0°C to 50°C 50°C to 660°C 660°C to 1100°C 1100°C to 1300°C	0.5°C 0.45°C 0.4°C 0.7°C 1.7°C	Thermocouples without a cold junction will have increased uncertainty
Calibrations at fixed points			
FP Tin FP Zinc FP Aluminium FP Silver	231.928°C 419.527°C 660.323°C 961.78°C	0.4°C 0.4°C 0.4°C 0.4°C	
Other thermocouples	-196°C -80°C to 0°C 0°C to 50°C 50°C to 300°C 300°C to 420°C 420°C to 660°C 660°C to 1100°C 1100°C to 1300°C	0.3°C 0.25°C 0.1°C 0.25°C 0.30°C 0.4°C 0.8°C 2.2°C	

TP = Triple Point FP = Freezing Point MP = Melting Point BP = Boiling Point



The latest schedule can be found on the Isotech website or at www.ukas.org.

ISOTECH

Isotech Training Course

From Industry to ITS-90



Isotech Training Course, "From Industry to ITS-90"

Isotech is widely recognized as a leading provider of training for temperature calibration, Isotech's new course is now available.

- Isotech: A brief introduction and history
- Fundamental Questions: What is Temperature?
- The Two Methods of Calibration
- Equipping for Industrial Calibration
- Equipping a Secondary Laboratory
- Equipping a Primary Laboratory
- Degrees of Automation
- Understanding Sources of Error
- Traceability
- Uncertainty

Who is it aimed at?

- People who are considering starting to calibrate sensors
- Those new to Temperature Calibration
- Industrial Users seeking a fuller picture of calibration
- People who are planning to purchase new equipment

Benefits

- An understanding of temperature calibration
- The methods of calibration
- The types of equipment available
- How to select the appropriate equipment
- How to avoid common pitfalls





Where?

The course is available at Isotech and is also being offered in conjunction with Isotech Distributors across the globe. For those not able to attend the training is also offered as an internet E-learning opportunity via online modules.

The online training modules explain Isotech's work and products and provide an invaluable insight into areas such as calibrating sensors and equipping calibration laboratories.

The resources have been developed as part of the e-Learning Place. This £3million project has enabled Isotech, in partnership with the University of Liverpool and Connect Internet Solutions, to create high-quality courses incorporating multimedia features and on-line interactive assessment.

An access pack which gives 12 month access to the course along with a full set of written notes is available

The elearningPlac	e Liame
SUTECH Transport	Lane science a program a galanti rage intercentente de la substance ante de program a substance en la substance de la substance ante de program a substance a substance de la

SUITER		sector (Mark) 200 (Sector
-		
	Thereastate chain	* Finder angewend is to some a damaged in the second law values and a compared of addition water begins the compared is to some place of a statements for a some place and the source place of a statements for a some place and place and the place and the source place and the source place.
-	in the second se	A De annual esta companya de la comp
	UKAS	All radiust "the visual and in a listence relation All radiust "the visual and in a listence relation All radiust and the second relation of the radiust and the All radiust and the radiust and the radiust and the All radiust





Industrial Platinum Resistance Thermometer Tables - $R(0) = 100.00\Omega$

-200 18.52 -190 22.83 22.40 21.97 21.54 21.11 20.68 20.25 19.82 19.38 18.95 18.52 -170 31.34 30.91 30.07 28.64 29.22 28.60 28.37 27.95 27.152 27.14 20.51 53.55 53.45 25.474 25.03 52.52 52.11 110.00 60.26 65.51 65.51 65.51 65.51 65.51 65.61 66.17 66.17 66.17 66.17 66.17	°C ITS 90
-160 77.10 26.67 26.24 25.39 24.37 24.54 24.11 23.86 22.25 22.85 -170 31.34 30.01 30.42 30.07 29.64 29.27 28.80 28.57 27.95 27.52 27.10 -160 35.54 35.12 34.70 34.28 33.46 33.02 32.18 31.76 31.34 -150 39.72 39.31 38.89 38.47 38.05 37.44 37.22 38.80 36.34 33.96 35.54 -140 43.88 43.46 43.05 42.68 42.44 48.83 48.02 48.80 -110 66.19 55.79 53.85 54.97 54.35 53.75 53.45 57.41 57.01 66.80 66.26 80 68.33 67.92 67.42 66.72 66.31 65.91 65.51 65.11 64.70 64.33 -70 72.33 71.33 71.33 71.33 71.37	-200
$\begin{array}{c} -170 & 31.34 & 30.91 & 30.49 & 30.07 & 29.64 & 29.27 & 28.80 & 28.37 & 27.95 & 27.52 & 27.10 \\ -160 & 35.54 & 35.12 & 34.70 & 34.28 & 33.86 & 33.44 & 33.02 & 32.60 & 32.18 & 31.76 & 31.34 \\ -150 & 39.72 & 39.31 & 38.89 & 38.47 & 38.05 & 37.64 & 37.22 & 58.80 & 36.38 & 35.96 & 35.54 \\ -140 & 43.88 & 43.46 & 43.05 & 42.63 & 42.22 & 41.80 & 41.39 & 40.97 & 40.56 & 40.14 & 39.72 \\ -130 & 48.00 & 47.59 & 47.18 & 46.77 & 46.36 & 45.94 & 45.53 & 45.12 & 44.70 & 44.29 & 43.88 \\ -120 & 52.11 & 51.70 & 51.29 & 50.88 & 50.47 & 50.06 & 49.65 & 49.24 & 48.83 & 48.42 & 48.00 \\ -110 & 56.19 & 55.79 & 55.38 & 59.44 & 50.40 & 58.63 & 57.42 & 57.41 & 57.01 & 56.60 & 56.19 \\ -90 & 64.30 & 63.90 & 63.49 & 63.09 & 62.68 & 62.28 & 61.88 & 61.47 & 61.07 & 60.66 & 60.26 \\ -80 & 68.33 & 67.92 & 67.52 & 67.12 & 66.72 & 66.31 & 65.91 & 65.51 & 66.71 & 64.30 \\ -70 & 72.33 & 71.93 & 71.53 & 71.13 & 70.73 & 70.33 & 69.93 & 69.53 & 69.13 & 66.73 & 68.33 \\ -60 & 76.33 & 75.93 & 75.53 & 75.13 & 77.47 & 74.33 & 73.39 & 73.53 & 73.13 & 72.73 & 72.33 \\ -40 & 84.27 & 83.87 & 83.48 & 83.08 & 82.69 & 82.29 & 81.89 & 81.50 & 81.10 & 80.70 & 80.31 \\ -30 & 86.22 & 87.83 & 87.43 & 87.04 & 86.64 & 86.25 & 85.65 & 65.46 & 85.06 & 46.67 & 86.33 \\ -10 & 100.00 & 99.61 & 99.22 & 98.83 & 98.44 & 98.04 & 97.65 & 97.26 & 96.64 & 50.69 \\ 0 & 100.00 & 99.61 & 99.22 & 98.83 & 98.44 & 98.04 & 97.65 & 97.26 & 96.64 & 50.69 \\ 0 & 100.00 & 100.39 & 100.78 & 101.17 & 101.56 & 101.95 & 102.34 & 102.73 & 103.12 & 103.51 & 103.90 \\ 10 & 100.00 & 99.61 & 99.22 & 98.83 & 98.44 & 98.04 & 97.65 & 97.26 & 96.64 & 97.66 & 97.26 & 96.79 \\ 0 & 100.00 & 99.61 & 99.22 & 98.83 & 98.44 & 98.04 & 97.65 & 97.26 & 96.64 & 97.06 \\ 0 & 100.00 & 100.39 & 100.78 & 101.17 & 101.56 & 101.95 & 102.34 & 102.73 & 103.15 & 103.590 \\ 10 & 100.00 & 99.61 & 99.22 & 98.83 & 98.44 & 98.04 & 97.65 & 97.26 & 96.64 & 97.06 \\ 0 & 100.00 & 100.39 & 100.78 & 101.17 & 101.56 & 101.93 & 103.13 & 105.2 & 103.57 & 103.81 & 103.57 \\ 10 & 100.03 & 100.78 & 101.17 & 101.56 & 101.93 & 103.71 & 103.51 & 103.9$	-190
$\begin{array}{c} -160 & 35.54 & 35.12 & 34.70 & 34.28 & 33.86 & 33.44 & 33.02 & 32.60 & 32.18 & 31.76 & 31.34 \\ -150 & 39.72 & 39.31 & 38.89 & 38.47 & 38.05 & 37.64 & 37.22 & 36.80 & 36.38 & 35.96 & 35.54 \\ -140 & 43.88 & 43.64 & 43.05 & 42.63 & 42.22 & 41.80 & 41.39 & 40.97 & 40.56 & 40.14 & 39.72 \\ -130 & 48.00 & 47.59 & 47.18 & 46.77 & 46.36 & 45.94 & 45.53 & 45.12 & 44.70 & 44.29 & 48.88 \\ 120 & 52.11 & 51.70 & 51.29 & 50.88 & 50.47 & 50.06 & 49.65 & 49.24 & 48.83 & 48.42 & 48.00 \\ -110 & 56.19 & 55.79 & 55.38 & 54.97 & 54.56 & 54.15 & 53.75 & 53.34 & 52.93 & 52.52 & 52.11 \\ -100 & 60.26 & 59.85 & 59.44 & 59.04 & 58.63 & 58.23 & 57.82 & 57.41 & 57.01 & 56.60 & 66.19 \\ -30 & 68.33 & 67.92 & 67.52 & 67.12 & 66.72 & 66.31 & 65.91 & 65.51 & 65.11 & 64.70 & 64.30 \\ -70 & 72.33 & 71.93 & 77.53 & 77.13 & 77.43 & 77.433 & 77.39 & 73.53 & 73.13 & 72.73 & 72.33 \\ -50 & 80.31 & 79.91 & 79.51 & 79.11 & 78.72 & 78.32 & 77.92 & 77.52 & 77.12 & 76.73 & 76.33 \\ -40 & 84.27 & 83.87 & 83.48 & 83.08 & 82.69 & 82.29 & 81.89 & 81.50 & 81.10 & 80.70 & 80.31 \\ -30 & 88.22 & 87.83 & 87.43 & 87.43 & 87.44 & 86.44 & 65.25 & 65.54 & 65.10 & 64.46 & 8427 \\ -20 & 92.16 & 91.77 & 91.37 & 90.98 & 90.59 & 90.19 & 89.80 & 99.40 & 89.40 & 89.61 \\ 0 & 100.00 & 90.96 & 99.50 & 94.91 & 94.52 & 94.12 & 93.73 & 93.73 & 33.25 & 92.55 & 92.16 \\ 0 & 100.00 & 90.61 & 99.22 & 98.83 & 98.44 & 98.04 & 97.65 & 97.26 & 96.87 & 96.48 & 96.09 \\ \hline 0 & 100.00 & 100.39 & 100.78 & 101.17 & 101.56 & 101.95 & 102.34 & 102.73 & 103.12 & 103.51 & 103.90 \\ 10 & 103.90 & 100.78 & 101.17 & 101.56 & 101.95 & 102.34 & 102.73 & 103.312 & 103.51 & 103.90 \\ 10 & 100.00 & 100.39 & 100.78 & 101.17 & 101.56 & 101.95 & 102.34 & 102.73 & 103.52 & 92.66 \\ 70 & 170.00 & 190.42 & 106.48 & 105.07 & 104.81 & 105.71 & 105.55 & 105.64 & 126.69 & 127.08 \\ 70 & 127.08 & 127.44 & 128.47 & 128.47 & 128.47 & 128.47 & 118.54 & 118.63 & 116.77 & 117.95 \\ 110 & 142.92 & 126.67 & 127.48 & 123.49 & 124.29 & 124.54 & 1126.39 & 126.31 & 122.66 & 127.08 \\ 70 & 107.79 & 106.18 & 109.67 & 109.35 & 109.73 & 1$	-180
$\begin{array}{c} -150 & 39.72 & 39.31 & 38.89 & 38.47 & 38.05 & 37.64 & 37.22 & 38.80 & 36.38 & 35.96 & 35.54 \\ -140 & 43.88 & 43.46 & 43.05 & 42.63 & 42.22 & 41.80 & 41.39 & 40.97 & 40.56 & 40.14 & 39.72 \\ -130 & 48.00 & 47.59 & 47.18 & 46.77 & 46.36 & 45.94 & 45.53 & 45.12 & 44.70 & 44.29 & 43.88 \\ -120 & 52.11 & 51.70 & 51.29 & 50.88 & 50.47 & 50.06 & 49.65 & 49.24 & 48.83 & 48.42 & 48.00 \\ -110 & 56.19 & 55.79 & 55.38 & 54.97 & 54.56 & 54.15 & 53.75 & 53.34 & 52.52 & 52.11 \\ -100 & 60.26 & 59.85 & 59.44 & 50.90 & 58.63 & 58.23 & 57.82 & 57.41 & 57.01 & 56.60 & 56.19 \\ -90 & 64.30 & 63.90 & 63.49 & 63.09 & 62.68 & 62.28 & 61.88 & 61.47 & 61.07 & 60.66 & 60.26 \\ -80 & 68.33 & 67.92 & 67.52 & 67.12 & 66.72 & 66.31 & 65.91 & 65.51 & 65.11 & 64.70 & 64.80 \\ 710 & 72.33 & 71.93 & 71.53 & 77.13 & 77.73 & 77.43 & 77.39 & 77.53 & 77.31 & 72.73 & 72.33 \\ -50 & 76.33 & 75.93 & 75.53 & 75.13 & 77.17 & 77.82 & 77.52 & 77.12 & 77.73 & 76.33 \\ -40 & 84.27 & 83.87 & 83.48 & 83.08 & 82.69 & 82.29 & 81.89 & 81.50 & 81.10 & 80.70 & 80.31 \\ -30 & 88.22 & 87.83 & 87.43 & 87.04 & 86.64 & 86.25 & 85.56 & 85.46 & 85.06 & 84.67 & 8427 \\ -20 & 92.16 & 91.77 & 91.37 & 90.8 & 05.99 & 90.19 & 93.80 & 89.40 & 89.01 & 88.6 & 88.22 \\ -10 & 96.09 & 95.69 & 95.30 & 94.91 & 94.52 & 94.12 & 93.73 & 93.34 & 92.95 & 92.55 & 92.16 \\ 0 & 100.00 & 90.61 & 99.22 & 98.83 & 98.44 & 98.04 & 97.65 & 97.26 & 96.48 & 96.09 \\ \hline 0 & 100.00 & 100.39 & 100.78 & 101.17 & 101.56 & 110.92 & 110.23 & 110.51 & 110.30 & 111.29 \\ 111.67 & 112.06 & 112.45 & 112.83 & 113.22 & 113.61 & 114.00 & 114.38 & 114.77 & 115.15 & 115.54 \\ 0 & 110.9 & 104.29 & 104.68 & 105.07 & 105.46 & 105.85 & 106.24 & 106.63 & 107.9 & 106.89 & 107.9 \\ 20 & 107.79 & 105.18 & 106.57 & 100.54 & 102.84 & 102.73 & 103.13 & 103.52 & 130.90 \\ 10 & 103.90 & 104.29 & 104.68 & 105.07 & 105.46 & 105.85 & 106.24 & 106.63 & 107.9 & 111.29 \\ 111.67 & 112.06 & 112.45 & 112.83 & 113.22 & 113.61 & 114.00 & 114.38 & 114.77 & 115.15 & 115.54 \\ 0 & 115.94 & 115.93 & 116.31 & 116.70 & 117.49 & 117.49 & 128.29 & 128.65 & 1$	-170
$\begin{array}{c} -140 & 43.88 & 43.46 & 43.06 & 42.83 & 42.22 & 41.80 & 41.33 & 40.97 & 40.56 & 40.14 & 39.72 \\ -130 & 48.00 & 47.59 & 47.18 & 46.77 & 46.36 & 45.94 & 45.53 & 45.12 & 44.70 & 44.29 & 43.88 \\ -120 & 52.11 & 51.70 & 51.29 & 55.08 & 50.47 & 50.06 & 49.66 & 49.24 & 48.83 & 46.42 & 48.00 \\ -110 & 56.19 & 55.79 & 55.38 & 54.97 & 54.56 & 54.15 & 53.76 & 57.41 & 57.01 & 56.60 & 56.19 \\ -90 & 64.30 & 63.90 & 65.49 & 63.09 & 62.68 & 62.28 & 61.88 & 61.47 & 61.07 & 60.66 & 60.26 \\ -80 & 68.33 & 67.92 & 67.52 & 67.12 & 66.72 & 66.31 & 65.91 & 65.51 & 65.11 & 64.70 & 64.30 \\ -70 & 72.33 & 71.93 & 71.53 & 71.13 & 70.73 & 70.33 & 69.33 & 69.13 & 68.73 & 68.33 \\ -60 & 76.33 & 75.93 & 75.53 & 75.13 & 74.73 & 74.33 & 73.93 & 73.53 & 73.13 & 72.73 & 72.33 \\ -60 & 80.31 & 79.91 & 79.51 & 79.11 & 78.72 & 78.22 & 77.92 & 77.12 & 77.3 & 76.33 \\ -40 & 94.27 & 83.87 & 83.48 & 83.08 & 82.69 & 82.29 & 81.99 & 81.50 & 81.10 & 80.70 & 80.31 \\ -30 & 80.22 & 87.83 & 87.43 & 87.04 & 86.64 & 86.25 & 85.85 & 85.46 & 85.06 & 84.67 & 8427 \\ -20 & 92.16 & 91.77 & 91.37 & 90.98 & 90.59 & 90.19 & 89.80 & 99.40 & 89.40 & 88.62 & 88.22 \\ -10 & 90.00 & 90.69 & 95.30 & 94.19 & 93.73 & 93.73 & 33 & 73.33 & 92.55 & 92.16 \\ 0 & 100.00 & 100.39 & 100.78 & 101.17 & 101.56 & 101.95 & 102.34 & 102.73 & 103.12 & 103.51 & 103.90 \\ 10 & 103.90 & 104.29 & 104.68 & 105.07 & 105.46 & 105.85 & 106.24 & 106.63 & 107.02 & 107.40 & 107.79 \\ 20 & 107.79 & 103.18 & 103.57 & 103.84 & 103.61 & 114.00 & 114.38 & 114.77 & 115.15 & 115.54 \\ 40 & 115.54 & 115.93 & 116.31 & 116.70 & 117.08 & 117.47 & 117.86 & 116.24 & 118.63 & 119.01 & 119.40 \\ 50 & 119.40 & 119.78 & 112.82 & 113.61 & 114.00 & 114.38 & 114.77 & 115.15 & 115.54 \\ 10 & 132.94 & 132.84 & 133.24 & 132.84 & 133.84 & 133.71 & 33.52 & 130.90 \\ 10 & 132.90 & 134.81 & 133.57 & 133.18 & 133.57 & 133.59 & 134.33 & 134.71 \\ 90 & 134.71 & 135.09 & 135.47 & 135.85 & 136.23 & 136.61 & 136.99 & 137.37 & 138.13 & 1335.51 \\ 100 & 140.74 & 140.64 & 140.82 & 147.20 & 140.40 & 140.78 & 144.18 & 144.56 & 144.94 & 145.31 & 145.6$	-160
$\begin{array}{cccccccccccccccccccccccccccccccccccc$	-150
$\begin{array}{ c c c c c c c c c c c c c c c c c c c$	-140
-110 56.19 55.79 55.38 54.97 54.56 54.15 53.75 53.34 52.93 52.52 52.11 -100 60.26 59.85 59.44 59.04 58.63 58.23 57.82 57.41 57.01 56.60 60.26 90 64.30 63.90 62.86 62.28 61.88 61.47 61.07 64.07 64.30 -70 77.33 71.93 71.53 77.13 70.33 69.93 66.51 65.11 64.70 64.33 -60 76.33 75.53 75.13 74.73 74.33 73.93 73.53 73.13 72.73 76.33 -60 86.31 79.91 79.51 77.11 78.72 77.52 77.12 76.73 76.33 -30 88.22 87.48 80.40 80.25 85.85 85.46 85.06 84.67 84.27 -20 92.16 91.77 91.37 90.38 98.44 98.04	-130
-110 56.19 55.79 55.38 54.97 54.56 54.15 53.75 53.34 52.93 52.52 52.11 -100 60.26 59.85 59.44 59.04 58.63 58.23 57.82 57.41 57.01 56.60 60.26 90 64.30 63.90 62.86 62.28 61.88 61.47 61.07 64.07 64.30 -70 77.33 71.93 71.53 77.13 70.33 69.93 66.51 65.11 64.70 64.33 -60 76.33 75.53 75.13 74.73 74.33 73.93 73.53 73.13 72.73 76.33 -60 86.31 79.91 79.51 77.11 78.72 77.52 77.12 76.73 76.33 -30 88.22 87.48 80.40 80.25 85.85 85.46 85.06 84.67 84.27 -20 92.16 91.77 91.37 90.38 98.44 98.04	-120
$\begin{array}{ c c c c c c c c c c c c c c c c c c c$	-110
$\begin{array}{ c c c c c c c c c c c c c c c c c c c$	-100
$\begin{array}{ c c c c c c c c c c c c c c c c c c c$	-90
$\begin{array}{ c c c c c c c c c c c c c c c c c c c$	-80
$\begin{array}{ c c c c c c c c c c c c c c c c c c c$	-70
$ \begin{array}{ c c c c c c c c c c c c c c c c c c c$	-60
$\begin{array}{ c c c c c c c c c c c c c c c c c c c$	-50
$\begin{array}{ c c c c c c c c c c c c c c c c c c c$	-40
$\begin{array}{ c c c c c c c c c c c c c c c c c c c$	-30
$\begin{array}{ c c c c c c c c c c c c c c c c c c c$	-20
$ \begin{array}{ c c c c c c c c c c c c c c c c c c c$	-10
$ \begin{array}{ c c c c c c c c c c c c c c c c c c c$	0
$\begin{array}{ c c c c c c c c c c c c c c c c c c c$	0
$ \begin{array}{ c c c c c c c c c c c c c c c c c c c$	0
$\begin{array}{ c c c c c c c c c c c c c c c c c c c$	10
$\begin{array}{ c c c c c c c c c c c c c c c c c c c$	20
$\begin{array}{ c c c c c c c c c c c c c c c c c c c$	30
$ \begin{array}{ c c c c c c c c c c c c c c c c c c c$	40
$ \begin{array}{ c c c c c c c c c c c c c c c c c c c$	50
70 127.08 127.46 127.84 128.22 128.61 128.99 129.37 129.75 130.13 130.52 130.90 80 130.90 131.28 131.66 132.04 132.42 132.80 133.18 133.57 133.95 134.33 134.71 90 134.71 135.09 135.47 135.85 136.23 136.61 136.99 137.37 137.75 138.13 138.51 100 138.51 138.88 139.26 139.64 140.02 140.40 140.78 141.16 141.54 141.91 142.29 110 142.29 142.67 143.05 143.43 143.80 144.18 144.56 144.94 145.31 145.69 146.07 120 146.07 146.44 146.82 147.20 147.57 147.95 148.33 148.70 149.08 149.46 149.83 130 149.83 150.21 150.58 150.96 151.33 151.71 152.08 152.46	60
80 130.90 131.28 131.66 132.04 132.42 132.80 133.18 133.57 133.95 134.33 134.71 90 134.71 135.09 135.47 135.85 136.23 136.61 136.99 137.37 137.75 138.13 138.51 100 138.51 138.88 139.26 139.64 140.02 140.40 140.78 141.16 141.54 141.91 142.29 110 142.29 142.67 143.05 143.43 143.80 144.18 144.56 144.94 145.31 145.69 146.07 120 146.07 146.44 146.82 147.20 147.57 147.95 148.33 148.70 149.08 149.46 149.83 130 149.83 150.21 150.58 150.96 151.33 151.71 152.08 152.46 152.83 153.21 153.58 140 153.58 153.96 154.33 154.71 155.08 156.20 156.58 156.95	70
90134.71135.09135.47135.85136.23136.61136.99137.37137.75138.13138.51100138.51138.88139.26139.64140.02140.40140.78141.16141.54141.91142.29110142.29142.67143.05143.43143.80144.18144.56144.94145.31145.69146.07120146.07146.44146.82147.20147.57147.95148.33148.70149.08149.46149.83130149.83150.21150.58150.96151.33151.71152.08152.46152.83153.21153.58140153.58153.96154.33154.71155.08155.46155.83156.20156.58156.95157.33150157.33157.70158.07158.45158.82159.19159.94160.31160.68161.05160161.05161.43161.80162.17162.54162.91163.29163.66164.03164.40164.77170164.77165.14165.51165.89166.26166.63167.00167.37167.74168.11168.48180168.48168.85169.22169.96170.33170.70171.07171.43171.80172.17190172.17172.54172.91173.28173.65174.02174.38174.75175.12175.49175.86200175	80
100138.51138.88139.26139.64140.02140.40140.78141.16141.54141.91142.29110142.29142.67143.05143.43143.80144.18144.56144.94145.31145.69146.07120146.07146.44146.82147.20147.57147.95148.33148.70149.08149.46149.83130149.83150.21150.58150.96151.33151.71152.08152.46152.83153.21153.58140153.58153.96154.33154.71155.08155.46155.83156.20156.58156.95157.33150157.33157.70158.07158.45158.82159.19159.56159.94160.31160.68161.05160161.05161.43161.80162.17162.54162.91163.29163.66164.03164.40164.77170164.77165.14165.51165.89166.26166.63167.00167.37167.74168.11168.48180168.48168.85169.22169.59169.96170.33170.70171.07171.43171.80172.17190172.17172.54172.91173.28173.65174.02174.38174.75175.12175.49175.86200175.86176.22176.59176.96177.33177.69178.06178.43178.79179.16179.53 </td <td>90</td>	90
110142.29142.67143.05143.43143.80144.18144.56144.94145.31145.69146.07120146.07146.44146.82147.20147.57147.95148.33148.70149.08149.46149.83130149.83150.21150.58150.96151.33151.71152.08152.46152.83153.21153.58140153.58153.96154.33154.71155.08155.46155.83156.20156.58156.95157.33150157.33157.70158.07158.45158.82159.19159.56159.94160.31160.68161.05160161.05161.43161.80162.17162.54162.91163.29163.66164.03164.40164.77170164.77165.14165.51165.89166.26166.63167.00167.37167.74168.11168.48180168.48168.85169.22169.59169.96170.33170.70171.07171.43171.80172.17190172.17172.54172.91173.28173.65174.02174.38174.75175.12175.49175.86200175.86176.22176.59176.96177.33177.69178.06178.43178.79179.16179.53210179.53179.89180.26180.63180.99181.36181.72182.09182.46182.82183.19 </td <td>100</td>	100
120146.07146.44146.82147.20147.57147.95148.33148.70149.08149.46149.83130149.83150.21150.58150.96151.33151.71152.08152.46152.83153.21153.58140153.58153.96154.33154.71155.08155.46155.83156.20156.58156.95157.33150157.33157.70158.07158.45158.82159.19159.56159.94160.31160.68161.05160161.05161.43161.80162.17162.54162.91163.29163.66164.03164.40164.77170164.77165.14165.51165.89166.26166.63167.00167.37167.74168.11168.48180168.48168.85169.22169.59169.96170.33170.70171.07171.43171.80172.17190172.17172.54172.91173.28173.65174.02174.38174.75175.12175.49175.86200175.86176.22176.59176.96177.33177.69178.06178.43178.79179.16179.53210179.53179.89180.26180.63180.99181.36181.72182.09182.46182.82183.19220183.19183.55183.92184.28184.65185.01185.38185.74186.11186.47186.84 </td <td>110</td>	110
130149.83150.21150.58150.96151.33151.71152.08152.46152.83153.21153.58140153.58153.96154.33154.71155.08155.46155.83156.20156.58156.95157.33150157.33157.70158.07158.45158.82159.19159.56159.94160.31160.68161.05160161.05161.43161.80162.17162.54162.91163.29163.66164.03164.40164.77170164.77165.14165.51165.89166.26166.63167.00167.37167.74168.11168.48180168.48168.85169.22169.59169.96170.33170.70171.07171.43171.80172.17190172.17172.54172.91173.28173.65174.02174.38174.75175.12175.49175.86200175.86176.22176.96177.33177.69178.06178.43178.79179.16179.53210179.53179.89180.26180.63180.99181.36181.72182.09182.46182.82183.19220183.19183.55183.92184.28184.65185.01185.38185.74186.11186.47186.84230186.84187.20187.56187.93188.29188.66189.02189.38189.75190.11190.47 <td< td=""><td>120</td></td<>	120
140153.58153.96154.33154.71155.08155.46155.83156.20156.58156.95157.33150157.33157.70158.07158.45158.82159.19159.56159.94160.31160.68161.05160161.05161.43161.80162.17162.54162.91163.29163.66164.03164.40164.77170164.77165.14165.51165.89166.26166.63167.00167.37167.74168.11168.48180168.48168.85169.22169.59169.96170.33170.70171.07171.43171.80172.17190172.17172.54172.91173.28173.65174.02174.38174.75175.12175.49175.86200175.86176.22176.59176.96177.33177.69178.06178.43178.79179.16179.53210179.53179.89180.26180.63180.99181.36181.72182.09182.46182.82183.19220183.19183.55183.92184.28184.65185.01185.38189.75190.11190.47230186.84187.20187.56187.93188.29188.66189.02189.38189.75190.11190.47240190.47190.84191.20191.56191.92192.29192.65193.01193.37193.74194.10 <td< td=""><td>130</td></td<>	130
150157.33157.70158.07158.45158.82159.19159.56159.94160.31160.68161.05160161.05161.43161.80162.17162.54162.91163.29163.66164.03164.40164.77170164.77165.14165.51165.89166.26166.63167.00167.37167.74168.11168.48180168.48168.85169.22169.59169.96170.33170.70171.07171.43171.80172.17190172.17172.54172.91173.28173.65174.02174.38174.75175.12175.49175.86200175.86176.22176.59176.96177.33177.69178.06178.43178.79179.16179.53210179.53179.89180.26180.63180.99181.36181.72182.09182.46182.82183.19220183.19183.55183.92184.28184.65185.01185.38185.74186.11186.47186.84230186.84187.20187.56187.93188.29188.66189.02189.38189.75190.11190.47240190.47190.84191.20191.56191.92192.29192.65193.01193.37193.74194.10250194.10194.46194.82195.18195.55195.91196.27196.63196.99197.35197.71 </td <td>140</td>	140
160161.05161.43161.80162.17162.54162.91163.29163.66164.03164.40164.77170164.77165.14165.51165.89166.26166.63167.00167.37167.74168.11168.48180168.48168.85169.22169.59169.96170.33170.70171.07171.43171.80172.17190172.17172.54172.91173.28173.65174.02174.38174.75175.12175.49175.86200175.86176.22176.59176.96177.33177.69178.06178.43178.79179.16179.53210179.53179.89180.26180.63180.99181.36181.72182.09182.46182.82183.19220183.19183.55183.92184.28184.65185.01185.38185.74186.11186.47186.84230186.84187.20187.56187.93188.29188.66189.02189.38189.75190.11190.47240190.47190.84191.20191.56191.92192.29192.65193.01193.37193.74194.10250194.10194.46194.82195.18195.55195.91196.27196.63196.99197.35197.71	150
170164.77165.14165.51165.89166.26166.63167.00167.37167.74168.11168.48180168.48168.85169.22169.59169.96170.33170.70171.07171.43171.80172.17190172.17172.54172.91173.28173.65174.02174.38174.75175.12175.49175.86200175.86176.22176.59176.96177.33177.69178.06178.43178.79179.16179.53210179.53179.89180.26180.63180.99181.36181.72182.09182.46182.82183.19220183.19183.55183.92184.28184.65185.01185.38185.74186.11186.47186.84230186.84187.20187.56187.93188.29188.66189.02189.38189.75190.11190.47240190.47190.84191.20191.56191.92192.29192.65193.01193.37193.74194.10250194.10194.46194.82195.18195.55195.91196.27196.63196.99197.35197.71	160
180168.48168.85169.22169.59169.96170.33170.70171.07171.43171.80172.17190172.17172.54172.91173.28173.65174.02174.38174.75175.12175.49175.86200175.86176.22176.59176.96177.33177.69178.06178.43178.79179.16179.53210179.53179.89180.26180.63180.99181.36181.72182.09182.46182.82183.19220183.19183.55183.92184.28184.65185.01185.38185.74186.11186.47186.84230186.84187.20187.56187.93188.29188.66189.02189.38189.75190.11190.47240190.47190.84191.20191.56191.92192.29192.65193.01193.37193.74194.10250194.10194.46194.82195.18195.55195.91196.27196.63196.99197.35197.71	170
190172.17172.54172.91173.28173.65174.02174.38174.75175.12175.12175.49175.86200175.86176.22176.59176.96177.33177.69178.06178.43178.79179.16179.53210179.53179.89180.26180.63180.99181.36181.72182.09182.46182.82183.19220183.19183.55183.92184.28184.65185.01185.38185.74186.11186.47186.84230186.84187.20187.56187.93188.29188.66189.02189.38189.75190.11190.47240190.47190.84191.20191.56191.92192.29192.65193.01193.37193.74194.10250194.10194.46194.82195.18195.55195.91196.27196.63196.99197.35197.71	180
200175.86176.22176.59176.96177.33177.69178.06178.43178.79179.16179.53210179.53179.89180.26180.63180.99181.36181.72182.09182.46182.82183.19220183.19183.55183.92184.28184.65185.01185.38185.74186.11186.47186.84230186.84187.20187.56187.93188.29188.66189.02189.38189.75190.11190.47240190.47190.84191.20191.56191.92192.29192.65193.01193.37193.74194.10250194.10194.46194.82195.18195.55195.91196.27196.63196.99197.35197.71	190
210179.53179.89180.26180.63180.99181.36181.72182.09182.46182.82183.19220183.19183.55183.92184.28184.65185.01185.38185.74186.11186.47186.84230186.84187.20187.56187.93188.29188.66189.02189.38189.75190.11190.47240190.47190.84191.20191.56191.92192.29192.65193.01193.37193.74194.10250194.10194.46194.82195.18195.55195.91196.27196.63196.99197.35197.71	200
220183.19183.55183.92184.28184.65185.01185.38185.74186.11186.47186.84230186.84187.20187.56187.93188.29188.66189.02189.38189.75190.11190.47240190.47190.84191.20191.56191.92192.29192.65193.01193.37193.74194.10250194.10194.46194.82195.18195.55195.91196.27196.63196.99197.35197.71	210
230186.84187.20187.56187.93188.29188.66189.02189.38189.75190.11190.47240190.47190.84191.20191.56191.92192.29192.65193.01193.37193.74194.10250194.10194.46194.82195.18195.55195.91196.27196.63196.99197.35197.71	220
240190.47190.84191.20191.56191.92192.29192.65193.01193.37193.74194.10250194.10194.46194.82195.18195.55195.91196.27196.63196.99197.35197.71	230
250 194.10 194.46 194.82 195.18 195.55 195.91 196.27 196.63 196.99 197.35 197.71	240
	250
260 197.71 198.07 198.43 198.79 199.15 199.51 199.87 200.23 200.59 200.95 201.31	250
<u>270</u> 201.31 201.67 202.03 202.39 202.75 203.11 203.47 203.83 204.19 204.55 204.90	270
	280
280 204.90 205.26 205.62 205.98 206.34 206.70 207.05 207.41 207.77 208.13 208.48 290 208.48 208.84 209.20 209.56 209.91 210.27 210.63 210.98 211.34 211.70 212.05	280
	300
<u>310</u> 215.61 215.96 216.32 216.67 217.03 217.38 217.74 218.09 218.44 218.80 219.15	310
<u>320</u> 219.15 219.51 219.86 220.21 220.57 220.92 221.27 221.63 221.98 222.33 222.68 230 232.68 232.04 232.30 232.74 234.00 234.45 234.80 235.15 235.50 235.85 236.21	320
330 222.68 223.04 223.39 223.74 224.09 224.45 224.80 225.15 225.50 225.85 226.21 340 226.21 226.56 226.91 227.26 227.61 227.96 228.31 228.66 229.02 229.37 229.72	330 340
<u>340 226.21 226.56 226.91 227.26 227.61 227.96 228.31 228.66 229.02 229.37 229.72</u>	340



Industrial Platinum Resistance Thermometer Tables - $R(0) = 100.00\Omega$

°C ITS 90	0	1	2	3	4	5	6	7	8	9	10	°C ITS 90
350	229.72	230.07	230.42	230.77	231.12	231.47	231.82	232.17	232.52	232.87	233.21	350
360	233.21	233.56	233.91	234.26	234.61	234.96	235.31	235.66	236.00	236.35	236.70	360
370	236.70	237.05	237.40	237.74	238.09	238.44	238.79	239.13	239.48	239.83	240.18	370
380	240.18	240.52	240.87	241.22	241.56	241.91	242.26	242.60	242.95	243.29	243.64	380
390	243.64	243.99	244.33	244.68	245.02	245.37	245.71	246.06	246.40	246.75	247.09	390
400	247.09	247.44	247.78	248.13	248.47	248.81	249.16	249.50	249.85	250.19	250.53	400
410	250.53	250.88	251.22	251.56	251.91	252.25	252.59	252.93	253.28	253.62	253.96	410
420	253.96	254.30	254.65	254.99	255.33	255.67	256.01	256.35	256.70	257.04	257.38	420
430	257.38	257.72	258.06	258.40	258.74	259.08	259.42	259.76	260.10	260.44	260.78	430
440	260.78	261.12	261.46	261.80	262.14	262.48	262.82	263.16	263.50	263.84	264.18	440
450	264.18	264.52	264.86	265.20	265.53	265.87	266.21	266.55	266.89	267.22	267.56	450
460	267.56	267.90	268.24	268.57	268.91	269.25	269.59	269.92	270.26	270.60	270.93	460
470	270.93	271.27	271.61	271.94	272.28	272.61	272.95	273.29	273.62	273.96	274.29	470
480	274.29	274.63	274.96	275.30	275.63	275.97	276.30	276.64	276.97	277.31	277.64	480
490	277.64	277.98	278.31	278.64	278.98	279.31	279.64	279.98	280.31	280.64	280.98	490
500	280.98	281.31	281.64	281.98	282.31	282.64	282.97	283.31	283.64	283.97	284.30	500
510	284.30	284.63	284.97	285.30	285.63	285.96	286.29	286.62	286.95	287.29	287.62	510
520	287.62	287.95	288.28	288.61	288.94	289.27	289.60	289.93	290.26	290.59	290.92	520
530	290.92	291.25	291.58	291.91	292.24	292.56	292.89	293.22	293.55	293.88	294.21	530
540	294.21	294.54	294.86	295.19	295.52	295.85	296.18	296.50	296.83	297.16	297.49	540
550	297.49	297.81	298.14	298.47	298.80	299.12	299.45	299.78	300.10	300.43	300.75	550
560	300.75	301.08	301.41	301.73	302.06	302.38	302.71	303.03	303.36	303.69	304.01	560
570	304.01	304.34	304.66	304.98	305.31	305.63	305.96	306.28	306.61	306.93	307.25	570
580	307.25	307.58	307.90	308.23	308.55	308.87	309.20	309.52	309.84	310.16	310.49	580
590	310.49	310.81	311.13	311.45	311.78	312.10	312.42	312.74	313.06	313.39	313.71	590
600	313.71	314.03	314.35	314.67	314.99	315.31	315.64	315.96	316.28	316.60	316.92	600
610	316.92	317.24	317.56	317.88	318.20	318.52	318.84	319.16	319.48	319.80	320.12	610
620	320.12	320.43	320.75	321.07	321.39	321.71	322.03	322.35	322.67	322.98	323.30	620
630	323.30	323.62	323.94	324.26	324.57	324.89	325.21	325.53	325.84	326.16	326.48	630
640	326.48	326.79	327.11	327.43	327.74	328.06	328.38	328.69	329.01	329.32	329.64	640
650	329.64	329.96	330.27	330.59	330.90	331.22	331.53	331.85	332.16	332.48	332.79	650
660	332.79	333.11	333.42	333.74	334.05	334.36	334.68	334.99	335.31	335.62	335.93	660
670	335.93	336.25	336.56	336.87	337.18	337.50	337.81	338.12	338.44	338.75	339.06	670
680	339.06	339.37	339.69	340.00	340.31	340.62	340.93	341.24	341.56	341.87	342.18	680
690	342,18	342.49	342.80	343.11	343.42	343.73	344.04	344.35	344.66	344.97	345.28	690
700	345.28	345.59	345.90	346.21	346.52	346.83	347.14	347.45	347.76	348.07	348.38	700
710	348.38	348.69	348.99	349.30	349.61	349.92	350.23	350.54	350.84	351.15	351.46	710
720	351.46	351.77	352.08	352.38	352.69	353.00	353.30	353.61	353.92	354.22	354.53	720
730	354.53	354.84	355.14	355.45	355.76	356.06	356.37	356.67	356.98	357.28	357.59	730
740	357.59	357.90	358.20	358.51	358.81	359.12	359.42	359.72	360.03	360.33	360.64	740
750	360.64	360.94	361.25	361.55	361.85	362.16	362.46	362.76	363.07	363.37	363.67	750
760	363.67	363.98	364.28	364.58	364.89	365.19	365.49	365.79	366.10	366.40	366.70	760
770	366.70	367.00	367.30	367.60	367.91	368.21	368.51	368.81	369.11	369.41	369.71	770
780	369.71	370.01	370.31	370.61	370.91	371.21	371.51	371.81	372.11	372.41	372.71	780
790	372.71	373.01	373.31	373.61	373.91	374.21	374.51	374.81	375.11	375.41	375.70	790
800	375.70	376.00	376.30	376.60	376.90	377.19	377.49	377.79	378.09	378.39	378.68	800
810	378.68	378.98	379.28	379.57	379.87	380.17	380.46	380.76	381.06	381.35	381.65	810
820	381.65	381.95	382.24	382.54	382.83	383.13	383.42	383.72	384.01	384.31	384.60	820
830	384.60	384.90	385.19	385.49	385.78	386.08	386.37	386.67	386.96	387.25	387.55	830
840	387.55	387.84	388.14	388.43	388.72	389.02	389.31	389.60	389.90	390.19	390.48	840
850	390.48		000.14	000.40	000.12	000.02	000.01	000.00	000.00	000.10	000.40	850
	000.40											000

The temperature/resistance relationships used in this standard are as follows:

for the range **-200°C to 0°C**. $R_t=R_o [1 + At + Bt^2 + C (t - 100°C) t^3]$

for the range of 0°C to 850°C. $R_t = R_o (1 + At + Bt^2)$ For the quality of platinum commonly used for industrial resistance themometers the values of the constants in these equations are:



The Source for Calibration Professionals

Telephone: +1(802) 863.8050

Fax: +1(802) 863.8125 Email: sales@isotechna.com

Isothermal North America 158 Brentwood Drive, Unit #4, Colchester, VT 05446



www.isotechna.com